

Cambridge Series

- ▶ **100 - Standard Standing Height Workstation**
- ▶ **110 - Standard Sitting Height Workstation**



* Please specify cross rail location - rear or center.

The 100 and 110 Series tables are strong and versatile. Their all-welded construction uses 16 gauge formed aprons, 2" square tubing, and welded bottom cross members to create a full perimeter frame capable of supporting an evenly distributed 1,000 pound load. Micro adjusting levelers are standard and provide a stable work surface on uneven floors.

100 Series: 36" Standing Height

30" Deep

Part No.	Description
T100-4830	48" W x 30" D x 36" H
T100-6030	60" W x 30" D x 36" H
T100-7230	72" W x 30" D x 36" H
T100-9630	96" W x 30" D x 36" H

36" Deep

Part No.	Description
T100-4836	48" W x 36" D x 36" H
T100-6036	60" W x 36" D x 36" H
T100-7236	72" W x 36" D x 36" H
T100-9636	96" W x 36" D x 36" H

110 Series: 30" Sitting Height

30" Deep

Part No.	Description
T110-4830	48" W x 30" D x 30" H
T110-6030	60" W x 30" D x 30" H
T110-7230	72" W x 30" D x 30" H
T110-9630	96" W x 30" D x 30" H

36" Deep

Part No.	Description
T110-4836	48" W x 36" D x 30" H
T110-6036	60" W x 36" D x 30" H
T110-7236	72" W x 36" D x 30" H
T110-9636	96" W x 36" D x 30" H

- ▶ **200 - Cantilever Standing Height Workstation**
- ▶ **210 - Cantilever Sitting Height Workstation**



* Workstation shown with optional modesty panel.

The cantilevered design makes the 200 and 210 Series a perfect choice in applications requiring operator clearance beneath the work surface. Their all-welded construction uses 16 gauge formed aprons, 2" square tubing and welded bottom cross members to create a full perimeter frame capable of supporting an evenly distributed 500 pound load. Micro adjusting levelers are standard and provide a stable work surface on uneven floors.

200 Series: 36" Standing Height

30" Deep

Part No.	Description
T200-4830	48" W x 30" D x 36" H
T200-6030	60" W x 30" D x 36" H
T200-7230	72" W x 30" D x 36" H
T200-9630	96" W x 30" D x 36" H

36" Deep

Part No.	Description
T200-4836	48" W x 36" D x 36" H
T200-6036	60" W x 36" D x 36" H
T200-7236	72" W x 36" D x 36" H
T200-9636	96" W x 36" D x 36" H

210 Series: 30" Sitting Height

30" Deep

Part No.	Description
T210-4830	48" W x 30" D x 30" H
T210-6030	60" W x 30" D x 30" H
T210-7230	72" W x 30" D x 30" H
T210-9630	96" W x 30" D x 30" H

36" Deep

Part No.	Description
T210-4836	48" W x 36" D x 30" H
T210-6036	60" W x 36" D x 30" H
T210-7236	72" W x 36" D x 30" H
T210-9636	96" W x 36" D x 30" H

Photos may depict products containing optional features. Work surfaces are not included; see page 8-4.

www.newenglandlab.com

Cambridge Series

▶ 300 - Ergonomic Workstation



The Ergonomic 300 Series features hand-crank, hydraulic telescoping legs, and a wide range of accessories. The 300 Series is well suited for light manufacturing use and is capable of supporting an evenly distributed 500 pound load. Note: 12" height adjustment range is standard from 30" - 42" high. Other ranges are available upon request.

30" Deep		36" Deep	
Part No.	Description	Part No.	Description
T300-4830	48" W x 30" D x Adj. Ht.	T300-4836	48" W x 36" D x Adj. Ht.
T300-6030	60" W x 30" D x Adj. Ht.	T300-6036	60" W x 36" D x Adj. Ht.
T300-7230	72" W x 30" D x Adj. Ht.	T300-7236	72" W x 36" D x Adj. Ht.

▶ 400 - Hydraulic Workstation



Designed with hydraulic leveling technology, our Cambridge 400 Series enables positive height adjustment with the simple turn of a handle. Elastomer floor contacts ensure strong, steady and level operation in the most demanding applications. Total capacity is an evenly distributed 1,000 pound load. Note: 12" height adjustment range is standard from 30" - 42" high. Other ranges are available upon request.

30" Deep		36" Deep	
Part No.	Description	Part No.	Description
T400-4830	48" W x 30" D x Adj. Ht.	T400-4836	48" W x 36" D x Adj. Ht.
T400-6030	60" W x 30" D x Adj. Ht.	T400-6036	60" W x 36" D x Adj. Ht.
T400-7230	72" W x 30" D x Adj. Ht.	T400-7236	72" W x 36" D x Adj. Ht.
T400-9630	96" W x 30" D x Adj. Ht.	T400-9636	96" W x 36" D x Adj. Ht.

* Please specify cross rail location - rear or center.

▶ 500 - Adjustable Height Workstation



Pushpin construction makes the 500 Series easy to raise and lower. This system is the economic choice when height adjustment is required. Note: 12" height adjustment range is standard from 30" - 42" high. Other ranges are available upon request.

30" Deep		36" Deep	
Part No.	Description	Part No.	Description
T500-4830	48" W x 30" D x Adj. Ht.	T500-4836	48" W x 36" D x Adj. Ht.
T500-6030	60" W x 30" D x Adj. Ht.	T500-6036	60" W x 36" D x Adj. Ht.
T500-7230	72" W x 30" D x Adj. Ht.	T500-7236	72" W x 36" D x Adj. Ht.
T500-9630	96" W x 30" D x Adj. Ht.	T500-9636	96" W x 36" D x Adj. Ht.

* Workstation shown with optional casters.
* Please specify cross rail location - rear or center.

Photos may depict products containing optional features.
Work surfaces are not included; see page 8-4.

Cambridge Series

► 600 - Full Height Flexible Workstation



Designed to match our Flex Series system, the 600 Series table provides maximum flexibility. Slotted uprights start at the base and are adjustable in 1" increments. A cantilevered work surface is provided in your choice of material. For storage space, suspended cabinets or shelving can be added. Please note not all Cambridge Series table accessories work with this system due to the special post design. We can accommodate most accessory requests, so please call your New England Lab representative for more details.

78" High / 30" Deep

Part No.	Description
T600-3630	36" W x 30" D x 78" H
T600-4830	48" W x 30" D x 78" H
T600-6030	60" W x 30" D x 78" H
T600-7230	72" W x 30" D x 78" H

78" High / 24" Deep

Part No.	Description
T600-3624	36" W x 24" D x 78" H
T600-4824	48" W x 24" D x 78" H
T600-6024	60" W x 24" D x 78" H
T600-7224	72" W x 24" D x 78" H

► 650 - Low Height Flexible Workstation



Similar to the 600, the 650 Series table provides maximum flexibility where height is not required. For storage space, suspended cabinets or shelving can be added. Please note not all Cambridge Series table accessories work with this system due to the special post design. We can accommodate most accessory requests, so please call your New England Lab representative for more details.


29" High / 30" Deep

Part No.	Description
T650-3630	36" W x 30" D x 29" H
T650-4830	48" W x 30" D x 29" H
T650-6030	60" W x 30" D x 29" H
T650-7230	72" W x 30" D x 29" H

29" High / 24" Deep

Part No.	Description
T650-3624	36" W x 24" D x 29" H
T650-4824	48" W x 24" D x 29" H
T650-6024	60" W x 24" D x 29" H
T650-7224	72" W x 24" D x 29" H

► Upper and Lower Shelving

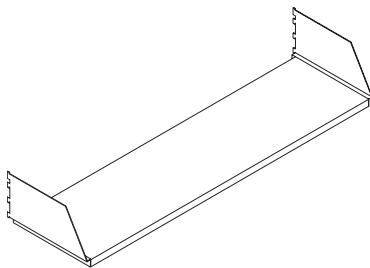

Made of 18 gauge painted steel with 5-3/4" high bookend brackets, the shelving for the 600 and 650 Series can be used both above and below the work surface. Multiple shelving units can be placed anywhere along the height of the post without interfering with the work surface support bracket.

Shelves

Part No.	Description
T605-2412	24" W x 12" D
T605-3012	30" W x 12" D
T605-3612	36" W x 12" D
T605-4812	48" W x 12" D
T605-6012	60" W x 12" D
T605-7212	72" W x 12" D
T606-2414	24" W x 14" D
T606-3014	30" W x 14" D
T606-3614	36" W x 14" D
T606-4814	48" W x 14" D
T606-6014	60" W x 14" D
T606-7214	72" W x 14" D

Part No.	Description
T607-2418	24" W x 18" D
T607-3018	30" W x 18" D
T607-3618	36" W x 18" D
T607-4818	48" W x 18" D
T607-6018	60" W x 18" D
T607-7218	72" W x 18" D
T608-2424	24" W x 24" D
T608-3024	30" W x 24" D
T608-3624	36" W x 24" D
T608-4824	48" W x 24" D
T608-6024	60" W x 24" D
T608-7224	72" W x 24" D

Please note this shelving system is different from other Cambridge Series tables, and will only fit the 600 or 650 workstations.

Photos may depict products containing optional features. Work surfaces are not included; see page 8-4.

www.newenglandlab.com

Work Surfaces

► Plastic Laminate Work Surface



Plastic Laminate Work Surfaces are the economic choice for the Cambridge table systems. It is constructed from medium density particle core, .050 high-pressure laminate, phenolic backer sheet, and a variety of edge treatments. This durable work surface is available in a wide selection of colors. Please indicate edge treatment.

*ESD Plastic Laminate with grounding strap also available.

Part No.	Description
__ 4830	48" W x 30" D
__ 6030	60" W x 30" D
__ 7230	72" W x 30" D
__ 9630	96" W x 30" D
__ 4836	48" W x 36" D
__ 6036	60" W x 36" D
__ 7236	72" W x 36" D
__ 9636	96" W x 36" D

TM = T-mold edgebanding
SE = Plastic laminate self edge
3M = 3mm PVC edgebanding
HB = 90° postformed, self edge
FB = 180° postformed, self edge

► Epoxy Resin Work Surface



Epoxy Resin Work Surfaces are the durable choice for the Cambridge table systems. Epoxy resin is resistant to heat, stains, and the corrosive effects of most laboratory chemicals. Epoxy work surface is 1" black.

Part No.	Description
ER 4830	48" W x 30" D
ER 6030	60" W x 30" D
ER 7230	72" W x 30" D
ER 9630	96" W x 30" D
ER 4836	48" W x 36" D
ER 6036	60" W x 36" D
ER 7236	72" W x 36" D
ER 9636	96" W x 36" D

► Phenolic Resin Work Surface



Phenolic Resin Work Surfaces are an outstanding choice to meet all of your laboratory needs. It is resistant to moisture, scratches and the corrosive effects of most laboratory chemicals. A variety of surface colors are available.

Part No.	Description
PH 4830	48" W x 30" D
PH 6030	60" W x 30" D
PH 7230	72" W x 30" D
PH 9630	96" W x 30" D
PH 4836	48" W x 36" D
PH 6036	60" W x 36" D
PH 7236	72" W x 36" D
PH 9636	96" W x 36" D

► Stainless Steel Work Surface



Resistant to heat, stains and corrosives, stainless steel is easy to clean and can have a visual impact on your laboratory. It is constructed from 16-gauge type 304 stainless steel with a #4 satin finish and features all metal hat channel supports. Type 316 stainless steel is also available upon request. Work surface (at left) shown with optional marine edge and backsplash.

Part No.	Description
SS 4830	48" W x 30" D
SS 6030	60" W x 30" D
SS 7230	72" W x 30" D
SS 9630	96" W x 30" D
SS 4836	48" W x 36" D
SS 6036	60" W x 36" D
SS 7236	72" W x 36" D
SS 9636	96" W x 36" D

Please contact your New England Lab representative for other sizes and custom configurations.

www.newenglandlab.com

Accessories for Cambridge Series Tables

► Casters



Casters can be added to most Cambridge Series tables to create mobile workstations. Our standard 4" diameter by 1-1/4" polyurethane wheels provide smooth rolling. All Casters are swivel type and include a dual locking mechanism for locking both swivel and wheel. Capacity is 250 pounds per caster. Priced per set of four.

Part No.	Description
T9450	Swivel Casters

► Vertical Uprights



Our Vertical Uprights are manufactured of 1-5/8" x 1-5/8" x 12 gauge channel. Vertical uprights bolt to the rear of the workstation for infinite adjustment of accessories above the work surface (12" required below the work surface for mounting). Priced per pair.

Part No.	Description
T9455	36" Uprights
T9460	48" Uprights
T9470	60" Uprights

► Flat Overhead Shelf



Our Flat Overhead Shelf is manufactured to be 1-1/4" thick with a high-pressure laminate surface, vinyl T-mold edgbanding and includes a heavy-duty welded support frame. All shelves are infinitely adjustable on our vertical uprights. Other shelving materials are available upon request.

Part No.	Description
T9200	48" W x 12" D
T9203	60" W x 12" D
T9206	72" W x 12" D
T9209	96" W x 12" D

► Sloped Overhead Shelf



The Sloped Overhead Shelf is manufactured to be 1-1/4" thick with a high-pressure plastic laminate surface, vinyl T-mold edgbanding, painted metal retainers on three sides, and includes a heavy-duty welded support frame. All shelves are infinitely adjustable on our vertical uprights. Other shelving materials are available upon request.

Part No.	Description
T9224	48" W x 12" D
T9227	60" W x 12" D
T9230	72" W x 12" D
T9233	96" W x 12" D

Photos may depict products containing optional features. Work surfaces are not included; see page 8-4.

Accessories for Cambridge Series Tables

► Lower Shelf



Our Lower Shelf is manufactured to be 1-1/4" thick with a high-pressure laminate surface and vinyl T-mold edgbanding. Other shelving materials are available upon request. Shelving depths are 14" on a 30" deep table and 17" on a 36" deep table. Center rail is required.

Part No.	Description
T9320-4814	48" W x 14" D
T9320-6014	60" W x 14" D
T9320-7214	72" W x 14" D
T9320-9614	96" W x 14" D
T9325-4817	48" W x 17" D
T9325-6017	60" W x 17" D
T9325-7217	72" W x 17" D
T9325-9617	96" W x 17" D

► Full Depth Lower Shelf



Our Full Depth Lower Shelf is manufactured to be 1-1/4" thick with a high-pressure laminate surface, and vinyl T-mold edgbanding. Other shelving materials are available upon request. Center rail is required.

Part No.	Description
T9340-4830	48" W x 30" D
T9340-6030	60" W x 30" D
T9340-7230	72" W x 30" D
T9340-9630	96" W x 30" D
T9345-4836	48" W x 36" D
T9345-6036	60" W x 36" D
T9345-7236	72" W x 36" D
T9345-9636	96" W x 36" D

► Tubular Instrument Shelf



Our heavy-duty Tubular Instrument Shelf offers ample support for heavy equipment and easily fastens to the work surface. Standard construction is 18" high x 10" deep with a 1 1/4" plastic laminate shelf and vinyl T-mold edgbanding. Other shelving materials, heights and depths are available upon request.

Part No.	Description
T9149	48" W x 12" D x 18" H
T9154	60" W x 12" D x 18" H
T9159	72" W x 12" D x 18" H
T9164	96" W x 12" D x 18" H

► Modesty Panel



Our Modesty Panel is available for added privacy at the workstation. Each panel fastens to the rear legs and is finished in decorative white melamine with black vinyl T-mold edgbanding.

Part No.	Description
T9350	48" W Modesty Panel
T9351	60" W Modesty Panel
T9352	72" W Modesty Panel
T9353	96" W Modesty Panel

Photos may depict products containing optional features. Work surfaces are not included; see page 8-4.

Accessories for Cambridge Series Tables

► Light Frame Assembly



Our Overhead Light Frame assembly includes a dual bulb light fixture (36" or 48") with an on/off switch and a 9 foot power cord. The frame is welded construction and is infinitely adjustable on our vertical uprights. The 48" frame assembly contains a 36" light fixture; all other sizes contain a 48" light fixture.

Part No.	Description
T9500	48" Light Frame
T9505	60" Light Frame
T9510	72" Light Frame
T9515	96" Light Frame

► Task Light



Lights can be mounted to the underside of any of our shelves. An on/off switch and a 9 foot power cord are included. (Shelf not included)

Part No.	Description
T9461	36" Task Light
T9471	48" Task Light

► Electrical Strip Assembly



Our standard Electrical Strip is 15 amperes, 125 volts with an on/off switch, circuit breaker, power indicator light, and a 6 foot power cord. Color is grey. Mounting rail is included and is infinitely adjustable on our vertical uprights.

Part No.	Description
T9648	48" L (8) Outlets
T9660	60" L (8) Outlets
T9672	72" L (8) Outlets
T9696	96" L (8) Outlets

► Book Stop Rail



Our Book Stop Rail is manufactured of heavy gauge, formed steel and mounts to uprights at the desired height above the shelf.

Part No.	Description
T9545	48" Book Stop
T9550	60" Book Stop
T9555	72" Book Stop
T9560	96" Book Stop

Photos may depict products containing optional features. Work surfaces are not included; see page 8-4.

Accessories for Cambridge Series Tables

► Bookends



Bookends are fabricated of 13-gauge steel with painted finish and are mechanically fastened to the underside of the shelf. Priced per pair.

Part No.	Description
T9362	5" H x 11-3/4" D

► Monitor Arm



Our standard Monitor Arm comes with a 13" x 13" tray with 360° swivel rotation and a dual pivot arm assembly. A clip-on keyboard tray (not shown) is optional.

Part No.	Description
T9621	Monitor Arm Assembly
T9622	Clip-on Keyboard Tray

► Keyboard Tray



Our pullout Keyboard Tray provides 7" height adjustment, +15°/ -15° tilt adjustment, 360° swivel, and a 23" x 11" keyboard platform.

Part No.	Description
T7400	Keyboard Tray
T7401	Keyboard Tray with Mouse Extension

Accessories for Cambridge Series Tables

► CPU Holder



The CPU holds most computer towers and pulls out for easy access. A 1/2" lip at the rear prevents the CPU from sliding backwards when the holder is pulled out. 20" H x 10 1/2" W x 22" D.

Part No.	Description
T7410	CPU Holder

► Pencil Drawer



The pencil drawer is 2 1/2" high and can be positioned anywhere between the legs of the table. The drawer comes equipped with full extension ball bearing slides.

Part No.	Description
T7415	Pencil Drawer

► Pull Out Writing Tablet



The writing tablet is 2 1/2" high and can be positioned anywhere between the legs of the table. The unit has a steel writing surface on full extension ball bearing slides.

Part No.	Description
T7420	Pull Out Writing Tablet


Photos may depict products containing optional features. Work surfaces are not included; see page 8-4.

Suspended Cabinets for Cambridge Series Tables


► Suspended Metal Cabinets



The Cambridge Series offers a wide variety of suspended storage solutions. The 22" deep painted metal cabinets are available in two widths, two heights, and a range of drawer and door/drawer configurations. All suspended metal cabinets hang on pre-installed metal cross rails, allowing "no tool" installation and easy side-to-side relocation. All cabinets are supplied with finished backs. A wide range of colors are available from New England Lab's color chart.

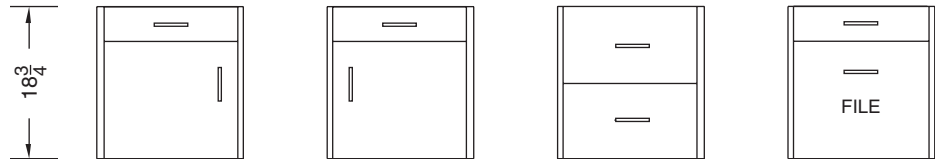

Part No.	Width	Part No.	Width	Part No.	Width	Part No.	Width
8511-15	15" W	8611-15	15" W	8099-15	15" W	8911-15	15" W
8511-18	18" W	8611-18	18" W	8099-18	18" W	8911-18	18" W
8511-24	24" W	8611-24	24" W	8099-24	24" W	8911-24	24" W

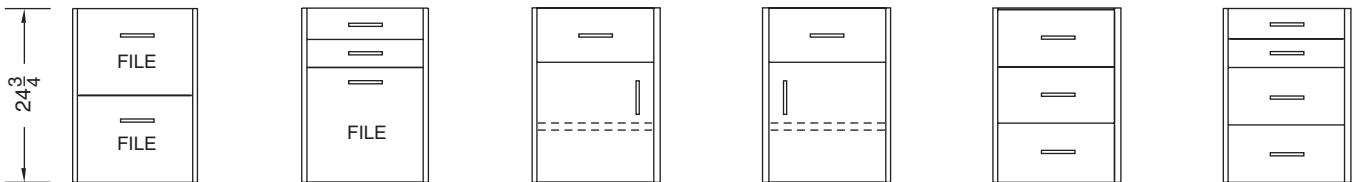

Part No.	Width	Part No.	Width	Part No.	Width	Part No.	Width	Part No.	Width	Part No.	Width
9901-15	15" W	9921-15	15" W	9511-15	15" W	9611-15	15" W	9033-15	15" W	9044-15	15" W
9901-18	18" W	9921-18	18" W	9511-18	18" W	9611-18	18" W	9033-18	18" W	9044-18	18" W
9901-24	24" W	9921-24	24" W	9511-24	24" W	9611-24	24" W	9033-24	24" W	9044-24	24" W

► Suspended Metal Cabinets



For maximum storage, wider models are available in four different widths, two heights, and a variety of drawer and door/drawer configurations. A wide range of colors are available from New England Lab's color chart.



Part No.	Width	Part No.	Width	Part No.	Width	Part No.	Width
8411-30	30" W	9433-30	30" W	9444-30	30" W	9411-30	30" W
8411-36	36" W	9433-36	36" W	9444-36	36" W	9411-36	36" W
8411-42	42" W	9433-42	42" W	9444-42	42" W	9411-42	42" W
8411-48	48" W	9433-48	48" W	9444-48	48" W	9411-48	48" W


Cambridge Series

► Suspended Cabinets - Plastic Laminate




For economical storage solutions, the Cambridge Series offers a variety of suspended plastic laminate cabinets. All cabinets are constructed of white thermoset melamine board with black vinyl T-mold edging and measure 18" wide by 22" deep.


► 6" Drawer Units

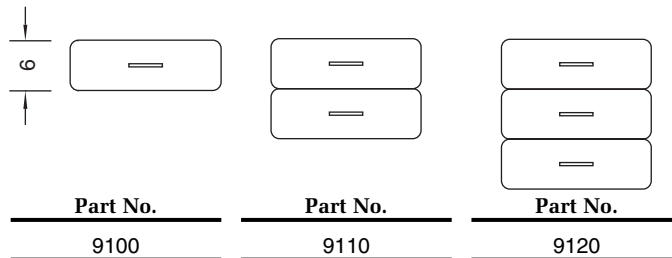

► 12" Drawer Units

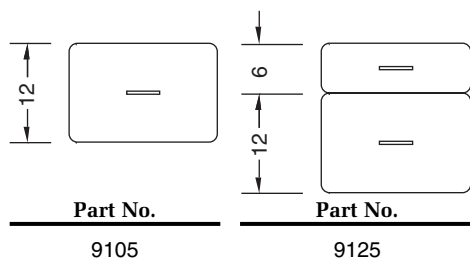

► Cupboard Units

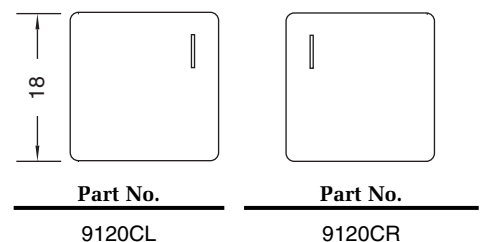

Specialty Tables

► PC Cart



The PC Cart is a great way to make your computer mobile. The PC Cart features durable locking casters, a push handle, an adjustable keyboard tray with a swing out mouse pad, and a bottom shelf with room for a computer tower and most printers. Cart measures 31" wide, 33" high and 24" deep.

Part No.	Description
T850	PC Cart

► Balance Table



New England Lab offers a quality balance table utilizing the superior chemical resistance of epoxy resin. This heavy table provides vibration free stability and maximum damping effect for analytical instrument scales. Top and legs are 2 1/2" thick, made of laminated epoxy. 1 1/2" black iron pipe provides additional stability. Standard tables are 30" high. Custom sizes are available. Tables are available in black.

Part No.	Description
T820-2424	24" W x 24" D x 30" H
T820-3024	30" W x 24" D x 30" H
T820-3030	30" W x 30" D x 30" H
T820-3630	36" W x 30" D x 30" H

Clean Tables

► Standard Make-Up Tables



Clean Table with H-Frame

Built for cleanliness and stability, these tables feature a seamless 14-gauge stainless steel work surface, with support structures, posts and sturdy construction throughout. Units are shipped knocked down and can be assembled in minutes without tools.

All Clean Tables are 34" high with stationary posts and leveling feet.

30" Deep Tables

Part No.	Length (Inches)	Post Material	Bottom Shelf Material	Bottom 3-Sided Frame Material	Bottom H-Frame Material
T720-6030	60	SS	SS	-	-
T720-7230	72	SS	SS	-	-
T720-9630	96	SS	SS	-	-
T721-6030	60	Chrome	Galvanized	-	-
T721-7230	72	Chrome	Galvanized	-	-
T721-9630	96	Chrome	Galvanized	-	-
T722-6030	60	SS	-	SS	-
T722-7230	72	SS	-	SS	-
T722-9630	96	SS	-	SS	-
T723-6030	60	SS	-	-	SS
T723-7230	72	SS	-	-	SS
T723-9630	96	SS	-	-	SS



Clean Table with 3-Sided Frame

36" Deep Tables

Part No.	Length (Inches)	Post Material	Bottom Shelf Material	Bottom 3-Sided Frame Material	Bottom H-Frame Material
T720-6036	60	SS	SS	-	-
T720-7236	72	SS	SS	-	-
T720-9636	96	SS	SS	-	-
T721-6036	60	Chrome	Galvanized	-	-
T721-7236	72	Chrome	Galvanized	-	-
T721-9636	96	Chrome	Galvanized	-	-
T722-6036	60	SS	-	SS	-
T722-7236	72	SS	-	SS	-
T722-9636	96	SS	-	SS	-
T723-6036	60	SS	-	-	SS
T723-7236	72	SS	-	-	SS
T723-9636	96	SS	-	-	SS

Clean Tables

► Mobile Clean Tables



All Mobile Clean Tables are 34" high. Mobile units have special posts and four 5" poly casters - 2 swivel/2 brake.

30" Deep Mobile Tables

Part No.	Length (Inches)	Post Material	Bottom Shelf Material	Bottom 3-Sided Frame Material	Bottom H-Frame Material
T730-6030	60	SS	SS	-	-
T730-7230	72	SS	SS	-	-
T730-9630	96	SS	SS	-	-
T731-6030	60	Chrome	Galvanized	-	-
T731-7230	72	Chrome	Galvanized	-	-
T731-9630	96	Chrome	Galvanized	-	-
T732-6030	60	SS	-	SS	-
T732-7230	72	SS	-	SS	-
T732-9630	96	SS	-	SS	-
T733-6030	60	SS	-	-	SS
T733-7230	72	SS	-	-	SS
T733-9630	96	SS	-	-	SS

► Clean Tables With Overhead



All Clean Tables with Overhead have rear stationary posts approximately 88" high. Posts are fitted with adjustable leveling feet. Cantilever shelves for overhead models are sold separately.

30" Deep Tables

Part No.	Length (Inches)	Post Material	Bottom Shelf Material	Bottom 3-Sided Frame Material	Bottom H-Frame Material
T740-6030	60	SS	SS	-	-
T740-7230	72	SS	SS	-	-
T740-9630	96	SS	SS	-	-
T741-6030	60	Chrome	Galvanized	-	-
T741-7230	72	Chrome	Galvanized	-	-
T741-9630	96	Chrome	Galvanized	-	-
T742-6030	60	SS	-	SS	-
T742-7230	72	SS	-	SS	-
T742-9630	96	SS	-	SS	-