



CASEWORK CONSTRUCTION DETAILS

Overview

Base Cabinets

Each cabinet is fabricated from prime quality furniture grade cold rolled or type 304-4 or 316-4 stainless steel. All components are die formed, notched, keyed and tightly fitted to form a solidly constructed electro spot welded rigid unit.

Each base cabinet incorporates a fully enclosed 4" high by 3" deep toe space at the front of the cabinet base. All base cabinets are provided with four 3/8" diameter leveler bolts accessible through 3/4" diameter pre-punched holes in the cabinet base. Each cabinet is provided with black press plugs to cover access holes after installation. All sink cabinets have satin coated (galvanized) steel bottoms, powder coated to match cabinets.

All base cabinets have full height removable back panels, except for sink cabinets which have 8" height removable back panels to facilitate plumbing. Removable backs are optional on drawer units. All cabinets have been designed to be vermin resistant.

Base cabinet doors and drawers are interchangeable in the field without special tools. Doors open full width without center posts on 30", 36", 42" and 48" wide base cabinets.

Drawers

Full access drawer bodies are formed with a one piece, coved bottom for effective cleaning, and a reinforcing bend on top edges. Depth of drawers are equal to the body depth of a standard cabinet less two inches. Drawer fronts are 3/4" thick pan construction. The space between the inner and outer drawer front is filled with a sound deadening core. Full extension ball bearing drawer runners have integral stops to prevent accidental withdrawal. Other drawer runner options such as soft closing, and interlocking are available.

Doors

Hinged metal doors feature a sturdy 3/4" thick double-pan construction, with a sound-dampening core housed in the space between the panels. Doors are mounted on high performance 2-1/2" five knuckle 14 gauge hinges, attached with four screws, and are held firmly closed with adjustable nylon roller catches. Framed glass doors are supplied with 1/8" (3mm) clear float glass.

Sliding doors are similar to hinged doors. Each sliding door is top hung with two 1" diameter nylon ball bearing rollers in a steel galvanized double track. The bottom of each door runs in a double PVC track. Rubber bumpers and stops are supplied with all doors.

Unframed glass doors are 1/4" (6mm) clear float glass with ground edges and integral pulls, set in an extruded aluminum shoe with nylon wheels. Doors slide in a top and bottom extruded aluminum track.

All glass supplied in floor cabinets (sliding frameless, hinged framed, sliding framed glass door units) is 1/4" (6mm) clear float glass.

Shelves

Shelves are fully adjustable on 1" increments with all edges formed down and back 1". Pull out shelves are suspended the same as standard drawers.

Hardware

Hardware is not included with products unless specified.

Coating Finish

All parts are pre-treated with an iron phosphate multi-stage process that includes a final sealer and reverse osmosis water rinse. The result is excellent adhesion between the metal and the coating. All exposed surfaces, including door and drawer head interiors, are electro-statically coated with a thermosetting laboratory grade powder coating. A broad palette of colors are available including multiple color schemes. Our coating is quality controlled to ensure consistent color match.

Quality

Our commitment to quality through our ISO Quality Management System assures the reliable, long life of our laboratory furniture.





STEEL LABORATORY FURNITURE CONSTRUCTION

Built to withstand tough laboratory environments, Mott's steel cabinets set the standard for durability and strength. Our casework can be part of fixed installations or components of mobile, suspended or flexible systems. We offer extensive door and drawer configurations, a wide variety of options, colors and finishes - providing you with a variety of layout options to maximize storage in your lab.

Mott's steel casework can meet the most challenging specification. All components have been tested to SEFA 8 Recommended Practices to ensure many years of intensive use in demanding environments. Our commitment to quality through the ISO Quality Management System assures the reliable, long life of your laboratory furniture.

We manufacture our basic cabinet as outlined below. However, today's laboratories often require customization, and to meet your unique needs, we can custom design and manufacture to meet your exact requirements. Simply contact us with your needs.

Drawers - Drawer bodies are one piece, reinforced and available in painted cold rolled steel or stainless steel. Drawer fronts are double pan construction and include sound deadening material.

Field Adjustable Design - Standard base cabinet doors and drawers are interchangeable in the field without the use of special tools. This provides unique future adaptability for changing laboratory needs.

Adjustable Shelves - Cabinet shelves are conveniently adjustable on 1" increments for maximum flexibility.

Protective Barrier - Casework can also be constructed from Sigma Barrier™, satin galvanized steel that provides an additional level of protection against rust and corrosion in high moisture areas.

Sink Cabinets - Supplied with Sigma Barrier™ satin galvanized toe kick base and cabinet bottom to further prevent rust developing from commonly wet operating conditions. Cabinet Construction - Each cabinet is fabricated with quality workmanship using furniture grade cold rolled steel with high pre and post-consumer recycled content. All components are die formed, notched, welded and tightly fitted to provide maximum strength, rigidity, and consistent appearance.



Coating Finish and Color - A thermosetting laboratory grade powder coating delivers exceptional chemical and scratch resistance, excellent hardness, adhesion, and a quality luster. The long lasting finish is available with 15 standard color options, including multiple color schemes, and color matching capabilities.

Full Height Removable Back Panels

- Standard door base cabinet backs are easily removable without the use of tools. This provides effective vermin control, and very importantly, increased access to services.

Drawer Suspension - An array of drawer suspension systems are available: 100lbs to 200lbs load capacity and features such as full extension, soft-close and selfclosing.

Handles and Hinges - A wide selection of handles and hinges are available to meet specific aesthetic and functional requirements.

Positive Door Catch - Optional self closing catch on inset cabinets that prevents rebounding. The body of the catch is built into the cabinet and the latching pin is mounted to the door for secure, solid closure. Standard on overlay cabinet design.

Doors - Double pan construction with sound deadening material and durable five knuckle hinges are bolted through the cabinet frame and are standard on all doors.





SEFA 8 LABORATORY FURNITURE

Recommended Practices

The Scientific Equipment and Furniture Association (SEFA) is a voluntary international trade association representing members of the laboratory furniture, casework, fume hood, and related laboratory equipment industry. Mott Manufacturing is a SEFA member.

The SEFA mission is to provide a forum for members to identify strategies to improve the quality, safety and timely completion of laboratory facilities to meet customer requirements.

SEFA has developed, through its membership, several recommended practices that are designed to promote better understanding between manufacturers and purchasers and to assist purchasers in selecting and specifying the proper product to meet the user's needs.

The SEFA 8 Recommended Practices provide manufacturers, specifiers, and users tools for evaluating the safety, durability, and structural integrity of laboratory casework.

Mott's line of laboratory furniture meets the performance requirements of the SEFA 8 Recommended Practices. A copy of the SEFA 8 Recommended Practices can be obtained from Mott Manufacturing complete with a certificate confirming the Sigma Systems line of laboratory furniture compliance to this performance standard.

The following is a summary of some of the SEFA 8 performance standards required of a standard 48" wide standing height base cabinet.

Cabinet Physical Strength Performance Tests

Cabinet Load Test

A 2,000 pound test load is mounted on top of a cabinet to evaluate its structural integrity and load bearing capability.

Cabinet Concentrated Load Test

A 200 pound center test load in mounted on top of a cabinet to evaluate its functional characteristics.

Cabinet Torsional Strength Test

A 200 pound test load is mounted on top of a cabinet which has been raised four inches off the floor with supports under three of the four corners to evaluate the cabinet's structural integrity.

Door Hinge Test

A 200 pound test load is hung from a door as it is cycled through the hinge full arc to demonstrate the door and its

hardware durability.

Door Impact Test

A 240 inch-pound impact test force is applied to the door face to evaluate its impact resistance.

Door Cycle Test

A door is swing tested through 90 degrees for 100,000 cycles to evaluate the hinge durability.

Drawer Load Test

A 150 pound test load is hung from an opened drawer to demonstrate its ability to support a point load at the drawer front and confirm the attachment of the drawer head to the body.

Drawer Cycle Test

A 100 pound static test load is distributed evenly in the drawer and the drawer is cycled 50,000 times from the fully open to fully closed position to confirm a reasonable life of the drawer slide mechanism under full load.

Drawer and Door Pull Test

A 50 pound test load is applied to the pull to evaluate the strength of the pull and its mounting hardware.

Drawer Impact Test

A 10 pound test load is dropped into an open drawer to evaluate the drawer and slide mechanism impact resistance.

Cabinet Surface Coating Performance Tests

Chemical Spot Test

The coated surface is exposed to 49 specified chemical reagents commonly found in laboratory applications to evaluate its resistance to chemical spills.

Hot Water Test

The coated surface is tested to evaluate its resistance to hot water.

Impact Test

A one pound ball is dropped on the surface to evaluate the coating's ductility.

Adhesion Test

A scratch test conducted according to ASTM D2197-86 is conducted to confirm adequate coating adhesion.

Hardness Test

A 4H pencil lead is used to determine the hardness of the coating.







TYPICAL WORK AREA

Details


The following details are provided to illustrate how base cabinets, filler panels and countertops combine to construct a laboratory work area.

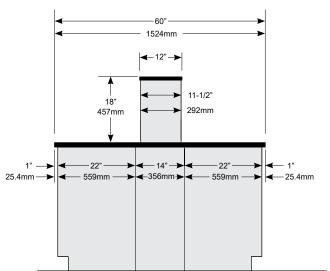

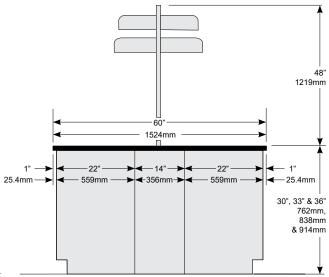


Island type configuration with 60" wide counter top, base cabinets, 14" filler panel and reagent shelf. A large space for services is concealed behind the 14" filler panel.

Island type configuration with 60" wide counter top, base cabinets and 14" filler panel. Optional center post assembly with shelving. A large space for services is concealed behind the 14" filler panel.

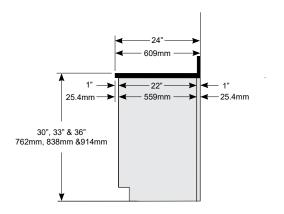

Standard base cabinet configuration with 24" wide counter top, base cabinet and filler 1" panel.

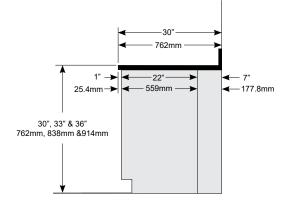

Base cabinet configuration with 30" wide counter top, base cabinet and filler 7" panel. The area covered by the 7" filler panel provides space for required services.





CABINET STYLES

The full range of Mott cabinetry is available in the following styles:

Our steel cabinets can be fitted with your choice of door and drawer fronts made from painted steel, stainless steel, wood, plastic laminate or phenolic resin. Typically designed with overlay construction, inset fronts can also be provided. Our flexibility in design allows you to create a unique lab that appeals to aesthetics, while also meeting the varied preferences and complex needs of the laboratory user.

Inset Cabinets

Powder Coated Steel Inset Cabinet



The long-lasting powder coated finish, with over 15 standard color options, including multiple color schemes can also be matched to custom colors. Standard construction.

Stainless Steel Inset Cabinet



Type 304 or 316 stainless steel are manufactured to a No.4 finish. Use option code 51 for type 304 stainless steel. Use option code 52 for type 316 stainless steel.

Architectural Overlay Cabinets

Powder Coated Steel Architectural Overlay Cabinet



To specify Architectural Overlay design use option code OL. The long-lasting powder coated finish, with over 15 standard color options, including multiple color schemes can also be matched to custom colors.

Stainless Steel Architectural Overlay Cabinet



To specify Architectural Overlay design use option code OL. Type 304 or 316 stainless steel are manufactured to a No.4 finish. Use option code 51 for type 304 stainless steel. Use option code 52 for type 316 stainless steel.

Architectural Overlay Cabinets With Wood Veneer, Plastic Laminate, Phenolic Resin Fronts



Wood Veneer Fronts: Combine the strength of our steel casework with the aesthetic warmth of wood fronts. We can provide a wide range of wood veneers in various cuts, matching and sequencing options, for drawer and door fronts. We offer a variety of stains and custom staining to match existing casework. Typically designed with overlay construction, inset fronts can be provided.

Plastic Laminate Fronts: Integrate the durability of our steel cabinet construction with the endless combinations of high-pressure laminate fronts. Typically designed with overlay construction, inset fronts can also be provided.

Phenolic Resin Fronts: Pair our rugged steel cabinet construction with our wide selection of phenolic resin fronts to create a resilient cabinet for harsh environments. Typically designed with overlay construction, inset fronts can also be provided.

SIGMA BARRIER™ COATING



Mott's Sigma Barrier™ coating provides excellent moisture and corrosion resistance in applications that involve high levels of wetness and humidity. Our full range of powder coated products are available with Sigma Barrier™ coating. It is standard on the bottom of all sink cabinets to protect against rusting in moist, high abuse areas. To specify Sigma Barrier™ use option code 53. Contact Mott for detailed Sigma Barrier™ test results.





OPTIONS

Door/Drawer Handles

Several different styles of door/drawer handles that are offered to suit specific customer requirements. Please use the option numbers below to specify the handle of your choice.

* Provide complete specification and source for all special options ordered.

Option Code	Description
16	Ship handles not installed
18	Horizontally mounted handles
29	Other customer specified hinge*
69	Other customer specified handle*
6A	Recessed stainless steel tapered pull



OPTION 62 (Standard) 4" Brushed Aluminum Pull



OPTION 60 4" Black Nylon Pull



OPTION 61
4" Stainless Steel Wire Pull
OPTION 68
4" Aluminum Wire Pull



OPTION 63 4" Chrome Pull



OPTION 64 5" Black Flush ABS Pull 64B Beige, 64G Gray, 64W White



4" Flush Aluminum Pull For Sliding Doors



OPTION 65
4" Flush Aluminum Pull
For Hinged Doors and
Drawers



OPTION 664" Flush Aluminum Painted Black Pull For Sliding Doors



OPTION 66 4" Flush Aluminum Painted Black Pull For Hinged Doors and Drawers



OPTION 67Flush Full Width Pull

Door Hinges

Mott cabinets come with durable 14 gauge, five knuckle hinges. Hinges are available in black powder coat, chrome plated or stainless steel finishes. Please use the option numbers below to specify the hinge of your choice.



OPTION 21 Chrome Hinge (Standard)

OPTION 22 Black Hinge

OPTION 23 Stainless Steel Hinge



Overlay Hinge for cabinets with steel fronts



Overlay Hinge for cabinets with wood fronts

OPTION 20 Dull Chrome Overlay Hinge OPTION 24 Stainless Steel Overlay Hinge

Drawer Body & Suspension

Drawers are available in powder coated cold rolled steel or stainless steel with the available options.

Option	Description
00	Modified door/drawer configuration
50	304 stainless steel drawer bodies on powder coated steel cabinet
26	200lb full extension drawer runners
F26	File drawer with 200lb full extension drawer runners
I01	100lb full extension drawer runners with interlock
126	200lb full extension drawer runners with interlock
S01	100lb full extension soft close drawer runners
S26	200lb full extension soft close drawer runners
59	Customer specified drawer runners (customer details required)





OPTIONS

Lock and Key Options

Option	Description
30	Standard locks on doors/drawers supplied and factory installed*
31	Provisions for standard locks on doors/drawers (holes & tabs)*
37	Hasp Lock (padlocks not included)
37A	Hasp Lock with stainless steel plate (pad lock not included)
38	Customer specified locks on doors/drawers supplied and factory installed**
39	Customer specified locks on doors/drawers supplied and not factory installed**
BS	Best locks (housing only), core not provided
LK5	5 pin tumbler with non-removable core, three level keying
LK6	5 pin deadbolt with non-removable core, three level keying

Notes: *Standard lock is a disc tumbler lock. Each cabinet is keyed differently per unit unless otherwise specified.


All Locks Keyed Different Per Unit



Locks are keyed different per unit from other units.


Note: This is the standard key schedule with disc tumbler lock supplied for option code 30 unless otherwise specified.

All Locks Keyed Alike

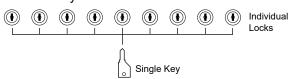


All locks are keyed alike, each lock can be opened with a single key.

Locks Keyed Alike in Sets

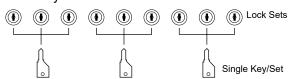


Locks are keyed alike in sets. Only one key will open one set of locks within the set. Up to 229 key changes for disc are available.

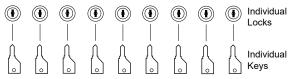


Locks are keyed different from other locks. Up to 229 key changes for disc are available.

Master Keyed Locks

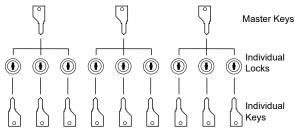

All locks are keyed different from all other locks. Each set within a group can be opened with a master key. Up to 229 key changes for disc are available.

Pin Locks Keyed Different, Grand Master & Master Keyed

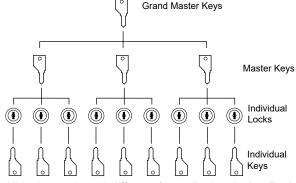

All locks are keyed different from all other locks. Each set within a group can be opened with a master key. All groups can be opened with a grand master key. 400 key changes and up to 10 master keys.

^{**}Provide complete specification and source for all special options ordered.





BASE CABINETS

34-3/4" HIGH X 22" DEEP

Open Units



Full Door Units







30" and wider

Door/Drawer Units



up to 24" wide



30" and wider



up to 24" wide

1 1 1		-	-	-
		-	-	-
	1			

30" and wider



Width	Item Number
30"	1310093
36"	1410093
42"	1810093
48"	1510093
1000mm	1710093

Width	Item Number
12"	1010000
15"	1910000
18"	1110000
24"	1210000
30"	1310000
36"	1410000
42"	1810000
48"	1510000
513mm	1610000
1000mm	1710000

Width	Right Hinged	Left Hinged	# of Doors
12"	1010011	1010211	1
15"	1910011	1910211	1
18"	1110011	1110211	1
24"	1210011	1210211	1
30"	1310011	-	2
36"	1410011	-	2
42"	1810011	-	2
48"	1510011	-	2
513mm	1610011	1610211	1
1000mm	1710011	-	2

Width	Right Hinged	Left Hinged	# of Doors
12"	1010022	1010222	1
15"	1910022	1910222	1
18"	1110022	1110222	1
24"	1210022	1210222	1
30"	1310022	-	2
36"	1410022	-	2
42"	1810022	-	2
48"	1510022	-	2
513mm	1610022	1610222	1
1000mm	1710022	-	2

Width	Right Hinged	Left Hinged	# of Doors
12"	1010033	1010233	1
15"	1910033	1910233	1
18"	1110033	1110233	1
24"	1210033	1210233	1
30"	1310033	-	2
36"	1410033	-	2
42"	1810033	-	2
48"	1510033	-	2
513mm	1610033	1610233	1
1000mm	1710033	-	2



Width	Item Number
30"	1310096
36"	1410096
42"	1810096
48"	1510096
1000mm	1710096





BASE CABINETS

34-3/4" HIGH X 22" DEEP

Door/Drawer Units



Width	Right Hand	Left Hand
30"	1310032	1310023
36"	1410032	1410023
42"	1810032	1810023
48"	1510032	1510023
1000mm	1710032	1710023

- · Left hand door/drawer unit shown
- Right hand reversed



Width	Right Hand	Left Hand
30"	1310021	1310012
36"	1410021	1410012
42"	1810021	1810012
48"	1510021	1510012
1000mm	1710021	1710012

- · Left hand door/drawer unit shown
- · Right hand reversed



Width	Right Hand	Left Hand
30"	1310031	1310013
36"	1410031	1410013
42"	1810031	1810013
48"	1510031	1510013
1000mm	1710031	1710013

- Left hand door/drawer unit shown
- Right hand reversed



Width	Right Hand	Left Hand
30"	1310041	1310014
36"	1410041	1410014
42"	1810041	1810014
48"	1510041	1510014
1000mm	1710041	1710014

- · Left hand door/drawer unit shown
- · Right hand reversed



Width	Right Hand	Left Hand
30"	1310051	1310015
36"	1410051	1410015
42"	1810051	1810015
48"	1510051	1510015
1000mm	1710051	1710015

- · Left hand door/drawer unit shown
- · Right hand reversed



Width	Right Hand	Left Hand
30"	1310061	1310016
36"	1410061	1410016
42"	1810061	1810016
48"	1510061	1510016
1000mm	1710061	1710016

- · Left hand door/drawer unit shown
- · Right hand reversed



Width	Right Hand	Left Hand
30"	1310071	1310017
36"	1410071	1410017
42"	1810071	1810017
48"	1510071	1510017
1000mm	1710071	1710017

- · Left hand door/drawer unit shown
- Right hand reversed



Width	Right Hand	Left Hand
30"	1310081	1310018
36"	1410081	1410018
42"	1810081	1810018
48"	1510081	1510018
1000mm	1710081	1710018

- · Left hand door/drawer unit shown
- · Right hand reversed
- · If locks are required, the lowest drawer of the cabinet will require a slot in the bottom of the drawer body to accommodate the lock cam



Width	Right Hand	Left Hand
30"	13100B1	131001B
36"	14100B1	141001B
42"	18100B1	181001B
48"	15100B1	151001B
1000mm	17100B1	171001B





Width	Right Hand	Left Hand
30"	1310942	1310924
36"	1410942	1410924
42"	1810942	1810924
48"	1510942	1510924
1000mm	1710942	1710924

- · Left hand door/drawer unit shown
- · Right hand reversed



1710042 · Left hand door/drawer unit shown

1810042

1510042

· Right hand reversed

42"

48"

1000mm

1810024

1510024

1710024





BASE CABINETS

34-3/4" HIGH X 22" DEEP

Door/Drawer Units



Width	Right Hand	Left Hand
30"	1310052	1310025
36"	1410052	1410025
42"	1810052	1810025
48"	1510052	1510025
1000mm	1710052	1710025

- · Left hand door/drawer unit shown
- Right hand reversed



Width	Right Hand	Left Hand
30"	1310072	1310027
36"	1410072	1410027
42"	1810072	1810027
48"	1510072	1510027
1000mm	1710072	1710027

- · Left hand door/drawer unit shown
- Right hand reversed



Width	Right Hand	Left Hand
30"	1310043	1310034
36"	1410043	1410034
42"	1810043	1810034
48"	1510043	1510034
1000mm	1710043	1710034

- · Left hand door/drawer unit shown
- · Right hand reversed



Width	Right Hand	Left Hand
30"	1310053	1310035
36"	1410053	1410035
42"	1810053	1810035
48"	1510053	1510035
1000mm	1710053	1710035

- · Left hand door/drawer unit shown
- · Right hand reversed



Width	Right Hand	Left Hand
30"	1310073	1310037
36"	1410073	1410037
42"	1810073	1810037
48"	1510073	1510037
1000mm	1710073	1710037

- · Left hand door/drawer unit shown
- · Right hand reversed



Width	Right Hand	Left Hand
30"	1310062	1310026
36"	1410062	1410026
42"	1810062	1810026
48"	1510062	1510026
1000mm	1710062	1710026

- · Left hand door/drawer unit shown
- · Right hand reversed



Width	Right Hand	Left Hand
30"	1310082	1310028
36"	1410082	1410028
42"	1810082	1810028
48"	1510082	1510028
1000mm	1710082	1710028

- · Left hand door/drawer unit shown
- · Right hand reversed
- If locks are required, the lowest drawer of the cabinet will require a slot in the bottom of the drawer body to accommodate the lock cam



Width	Right Hand	Left Hand
30"	1310063	1310036
36"	1410063	1410036
42"	1810063	1810036
48"	1510063	1510036
1000mm	1710063	1710036

- Left hand door/drawer unit shown
- · Right hand reversed



Width	Right Hand	Left Hand
30"	1310083	1310038
36"	1410083	1410038
42"	1810083	1810038
48"	1510083	1510038
1000mm	1710083	1710038

- · Left hand door/drawer unit shown
- · Right hand reversed
- If locks are required, the lowest drawer of the cabinet will require a slot in the bottom of the drawer body to accommodate the lock cam





BASE CABINETS

34-3/4" HIGH X 22" DEEP

Drawer Units



up to 24" wide



30" and wider



up to 24" wide



30" and wider



up to 24" wide

-	-
-	
-	-
-	-
-	-

30" and wider



up to 24" wide

	-
-	
-	-
-	-
-	-

30" and wider



up to 24" wide

-	
-	-
-	-
	-
-	-
	-

30" and wider

Width	Item Number
12"	1010090
15"	1910090
18"	1110090
24"	1210090
30"	13100NN
36"	14100NN
42"	18100NN
48"	15100NN
513mm	1610090
1000Mm	17100NN

Width	Item Number
12"	1010044
15"	1910044
18"	1110044
24"	1210044
30"	1310044
36"	1410044
42"	1810044
48"	1510044
513mm	1610044
1000mm	1710044

Width	Item Number	
15"	19100BB	
18"	11100BB	
24"	12100BB	
30"	13100BB	
36"	14100BB	
42"	18100BB	
48"	15100BB	
513mm	16100BB	
1000mm	17100BB	

Width	Item Number
15"	1910055
18"	1110055
24"	1210055
30"	1310055
36"	1410055
42"	1810055
48"	1510055
513mm	1610055
1000mm	1710055

Width	Item Number
15"	1910066
18"	1110066
24"	1210066
30"	1310066
36"	1410066
42"	1810066
48"	1510066
513mm	1610066
1000mm	1710066





BASE CABINETS

34-3/4" HIGH X 22" DEEP

Drawer Units



up to 24" wide



30" and wider



up to 24" wide



30" and wider

Width	Item Number
15"	1910077
18"	1110077
24"	1210077
30"	1310077
36"	1410077
42"	1810077
48"	1510077
513mm	1610077
1000mm	1710077

Width	Item Number
15"	1910088
18"	1110088
24"	1210088
30"	1310088
36"	1410088
42"	1810088
48"	1510088
513mm	1610088
1000mm	1710088

 If locks are required, the lowest drawer of the cabinet will require a slot in the bottom of the drawer body to accommodate the lock cam



Width	Item Number
30"	1310A91
36"	1410A91
42"	1810A91
48"	1510A91
1000mm	1710A91



Width	Item Number
30"	1310099
36"	1410099
42"	1810099
48"	1510099
1000mm	1710099



Width	Item Number
30"	1310090
36"	1410090
42"	1810090
48"	1510090
1000mm	1710090



Width	Item Number
30"	1310B44
36"	1410B44
42"	1810B44
48"	1510B44
1000mm	1710B44



Item Number
1310097
1410097
1810097
1510097
1710097



Width	Right Hand	Left Hand
30"	1310094	1310092
36"	1410094	1410092
42"	1810094	1810092
48"	1510094	1510092
1000mm	1710094	1710092

- · Left hand drawer unit shown
- · Right hand reversed





BASE CABINETS

34-3/4" HIGH X 22" DEEP

Drawer Units



Width	Item Number
30"	1310098
36"	1410098
42"	1810098
48"	1510098
1000mm	1710098



Width	Item Number
30"	1310091
36"	1410091
42"	1810091
48"	1510091
1000mm	1710091



Width	Item Number
30"	1310039
36"	1410039
42"	1810039
48"	1510039
1000mm	1710039



Width	Item Number
30"	1310095
36"	1410095
42"	1810095
48"	1510095
1000mm	1710095



Width	Item Number
30"	13100CC
36"	14100CC
42"	18100CC
48"	15100CC
1000mm	17100CC



Width	Right Hand	Left Hand
30"	1310019	1310029
36"	1410019	1410029
42"	1810019	1810029
48"	1510019	1510029
1000mm	1710019	1710029

- · Left hand drawer unit shown
- · Right hand reversed

File Drawer Units







30" and wider



up to 24" wide



30" and wider

Width	Item Number
12"	1010133
15"	1910133
18"	1110133
24"	1210133
30"	1310133
36"	1410133
42"	1810133
48"	1510133
513mm	1610133
1000mm	1710133

Width	Item Number
12"	1010100
15"	1910100
18"	1110100
24"	1210100
30"	1310100
36"	1410100
42"	1810100
48"	1510100
513mm	1610100
1000mm	1710100





BASE CABINETS

34-3/4" HIGH X 22" DEEP

File Drawer Units



up to 24" wide



30" and wider

Width	Item Number
12"	1010111
15"	1910111
18"	1110111
24"	1210111
30"	1310111
36"	1410111
42"	1810111
48"	1510111
513mm	1610111
1000mm	1710111



up to 24" wide



30" and wider

Width	Item Number
12"	1010122
15"	1910122
18"	1110122
24"	1210122
30"	1310122
36"	1410122
42"	1810122
48"	1510122
513mm	1610122
1000mm	1710122



up to 24" wide



30" and wider

Width	Item Number
12"	10101AA
15"	19101AA
18"	11101AA
24"	12101AA
30"	13101AA
36"	14101AA
42"	18101AA
48"	15101AA
513mm	16101AA
1000mm	17101AA





BASE CABINETS

34-3/4" HIGH X 22" DEEP

Waste/Storage Bin Units





up to 24" wide

30" and wider





Tilt Out Unit (shown with optional bin)

Pull Out Unit

Width	Tilt Out Units	Pull Out Units	Bin Capacity
15"	1910500	1910511	1
18"	1110500	1110511	1
24"	1210500	1210511	1
30"	1310500	1310511	2
36"	1410500	1410511	2
42"	1810500	1810511	2
48"	1510500	1510511	2
513mm	1610500	1610511	1

· Add option code 84 to order stainless steel waste/storage bin, or add option code 85 for PVC waste/storage bin for Tilt Out units

1710511

Tilt out units are hinged at the bottom


1710500

1000mm

Pull out units include full extension runners

Diagonal Corner Units



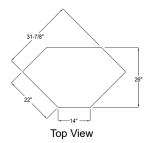

Description	Right Hinged	Left Hinged	
Full Door Unit	1017011	1017111	
False Panel Unit	1017422	1017922	

- · Left hinged, full door corner unit shown
- Add option code 89 to order a Lazy Susan


Sliding Door Units



Width	Item Number
30"	1310711
36"	1410711
42"	1810711
48"	1510711
1000mm	1710711

Blind Corner Units



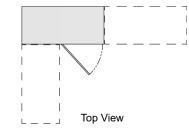

Width	Door on Right	Door on Left
30"	1310311	1310411
36"	1410311	1410411
42"	1810311	1810411
48"	1510311	1510411
1000mm	1710311	1710411

- · Right door blind corner unit shown, left door reversed
- · Center stile located at the middle point





BASE CABINETS

34-3/4" HIGH X 22" DEEP

Sink Units





up to 24" wide

30" and wider





up to 24" wide

30" and wider

Width **Right Hinged** Left Hinged # of Doors 12" 1013011 1013211 1 15" 1913011 1913211 18" 1113011 1113211 1 24" 1213011 1213211 2 30" 1313011 36" 2 1413011 42" 1813011 2 48" 1513011 2 1613011 1613211 513mm 1 1000mm 1713011 2

Width	Right Hinged	Left Hinged	# of Doors
12"	1018022	1018222	1
15"	1918022	1918222	1
18"	1118022	1118222	1
24"	1218022	1218222	1
30"	1318022	-	2
36"	1418022	-	2
42"	1818022	-	2
48"	1518022	-	2
513mm	1618022	1618222	1
1000mm	1718022	-	2

- · Add option code 90 to order Sink Hanger Supports
- Sink units are constructed with satin coated (galvanized) bottom, powder coated to match cabinet

End Sink Units




Description	54" Wide	58" Wide
Full Door Unit	9900017	9900018
False Panel Unit	9900021	9900022

- Full door end sink unit shown
- Add option code 90 to order Sink Hanger Supports
- Sink units are constructed with satin coated (galvanized) bottom, powder coated to match cabinet

Diagonal Corner Sink Units



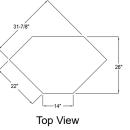

Description	Right Hinged	Left Hinged	
Full Door Unit	1017211	1017511	
False Panel Unit	1017722	1017622	

- · Left hinged, full door corner sink unit shown
- Add option code 90 to order Sink Hanger Supports
- Sink units are constructed with satin coated (galvanized) bottom, powder coated to match cabinet





BASE CABINETS

28-3/4" HIGH X 22" DEEP

Open Units



Full Door Units







30" and wider

Door/Drawer Units



up to 24" wide



30" and wider



up to 24" wide



30" and wider



Width	Item Number
30"	1320093
36"	1420093
42"	1820093
48"	1520093
1000mm	1720093

Width	Item Number
12"	1020000
15"	1920000
18"	1120000
24"	1220000
30"	1320000
36"	1420000
42"	1820000
48"	1520000
513mm	1620000
1000mm	1720000

Width	Right Hinged	Left Hinged	# of Doors
12"	1020011	1020211	1
15"	1920011	1920211	1
18"	1120011	1120211	1
24"	1220011	1220211	1
30"	1320011	-	2
36"	1420011	-	2
42"	1820011	-	2
48"	1520011	-	2
513mm	1620011	1620211	1
1000mm	1720011	-	2

Width	Right Hinged	Left Hinged	# of Doors
12"	1020022	1020222	1
15"	1920022	1920222	1
18"	1120022	1120222	1
24"	1220022	1220222	1
30"	1320022	-	2
36"	1420022	-	2
42"	1820022	-	2
48"	1520022	-	2
513mm	1620022	1620222	1
1000mm	1720022	-	2

Width	Item Number	Left Hinged	# of Doors
12"	1020033	1020233	1
15"	1920033	1920233	1
18"	1120033	1120233	1
24"	1220033	1220233	1
30"	1320033	-	2
36"	1420033	-	2
42"	1820033	-	2
48"	1520033	-	2
513mm	1620033	1620233	1
1000mm	1720033	-	2



Width	Item Number
30"	1320096
36"	1420096
42"	1820096
48"	1520096
1000mm	1720096





BASE CABINETS

28-3/4" HIGH X 22" DEEP

Door/Drawer Units



Width	Right Hand	Left Hand
30"	1320032	1320023
36"	1420032	1420023
42"	1820032	1820023
48"	1520032	1520023
1000mm	1720032	1720023

- Left hand door/drawer unit shown
- Right hand reversed



Width	Right Hand	Left Hand
30"	1320942	1320924
36"	1420942	1420924
42"	1820942	1820924
48"	1520942	1520924
1000mm	1780942	1720924

- · Left hand door/drawer unit shown
- Right hand reversed



Width	Right Hand	Left Hand
30"	1320021	1320012
36"	1420021	1420012
42"	1820021	1820012
48"	1520021	1520012
1000mm	1720021	1720012

- Left hand door/drawer unit shown
- Right hand reversed



Width	Right Hand	Left Hand
30"	1320042	1320024
36"	1420042	1420024
42"	1820042	1820024
48"	1520042	1520024
1000mm	1720042	1720024

- · Left hand door/drawer unit shown
- · Right hand reversed



Width	Right Hand	Left Hand
30"	1320031	1320013
36"	1420031	1420013
42"	1820031	1820013
48"	1520031	1520013
1000mm	1720031	1720013

- Left hand door/drawer unit shown
- · Right hand reversed



Width	Right Hand	Left Hand
30"	1320052	1320025
36"	1420052	1420025
42"	1820052	1820025
48"	1520052	1520025
1000mm	1720052	1720025

- · Left hand door/drawer unit shown
- · Right hand reversed



Width	Right Hand	Left Hand
30"	1320041	1320014
36"	1420041	1420014
42"	1820041	1820014
48"	1520041	1520014
1000mm	1720041	1720014

- Left hand door/drawer unit shown
- · Right hand reversed



Width	Right Hand	Left Hand
30"	1320062	1320026
36"	1420062	1420026
42"	1820062	1820026
48"	1520062	1520026
1000mm	1720062	1720026

- · Left hand door/drawer unit shown
- Right hand reversed



Width	Right Hand	Left Hand
30"	1320051	1320015
36"	1420051	1420015
42"	1820051	1820015
48"	1520051	1520015
1000mm	1720051	1720015

- · Left hand door/drawer unit shown
- · Right hand reversed



Width	Right Hand	Left Hand
30"	1320072	1320027
36"	1420072	1420027
42"	1820072	1820027
48"	1520072	1520027
1000mm	1720072	1720027

- · Left hand door/drawer unit shown
- Right hand reversed



Width	Right Hand	Left Hand
30"	1320061	1320016
36"	1420061	1420016
42"	1820061	1820016
48"	1520061	1520016
1000mm	1720061	1720016

- · Left hand door/drawer unit shown
- · Right hand reversed



Width	Dight Hand	Left Hand
wiatii	Right Hand	Leit Hailu
30"	1320043	1320034
36"	1420043	1420034
42"	1820043	1820034
48"	1520043	1520034
1000mm	1720043	1720034

- · Left hand door/drawer unit shown
- · Right hand reversed





BASE CABINETS

28-3/4" HIGH X 22" DEEP

Door/Drawer Units



Width	Right Hand	Left Hand
30"	1320053	1320035
36"	1420053	1420035
42"	1820053	1820035
48"	1520053	1520035
1000mm	1720053	1720035

- · Left hand door/drawer unit shown
- · Right hand reversed



Width	Right Hand	Left Hand
30"	1320063	1320036
36"	1420063	1420036
42"	1820063	1820036
48"	1520063	1520036
1000mm	1720063	1720036

- · Left hand door/drawer unit shown
- · Right hand reversed



Width	Right Hand	Left Hand
30"	1320073	1320037
36"	1420073	1420037
42"	1820073	1820037
48"	1520073	1520037
1000mm	1720073	1720037

- · Left hand door/drawer unit shown
- · Right hand reversed

Drawer Units



up to 24" wide



30" and wider



up to 24" wide



30" and wider



up to 24" wide



30" and wider

Width	Item Number
12"	1020044
15"	1920044
18"	1120044
24"	1220044
30"	1320044
36"	1420044
42"	1820044
48"	1520044
513mm	1620044
1000mm	1720044

Width	Item Number	
12"	10200MM	
15"	19200MM	
18"	11200MM	
24"	12200MM	
30"	13200MM	
36"	14200MM	
42"	18200MM	
48"	15200MM	
513mm	16200MM	
1000mm	17200MM	

Width	Item Number
12"	1020055
15"	1920055
18"	1120055
24"	1220055
30"	1320055
36"	1420055
42"	1820055
48"	1520055
513mm	1620055
1000mm	1720055





BASE CABINETS

28-3/4" HIGH X 22" DEEP

Drawer Units



up to 24" wide



30" and wider



up to 24" wide



30" and wider

Width	Item Number
15"	1920066
18"	1120066
24"	1220066
30"	1320066
36"	1420066
42"	1820066
48"	1520066
513mm	1620066
1000mm	1720066

Width	Item Number	
15"	1920077	
18"	1120077	
24"	1220077	
30"	1320077	
36"	1420077	
42"	1820077	
48"	1520077	
513mm	1620077	
1000mm	1720077	



Width	Item Number
30"	1320097
36"	1420097
42"	1820097
48"	1520097
1000mm	1720097



Width	Item Number
30"	13200LL
36"	14200LL
42"	18200LL
48"	15200LL
1000mm	17200LL



Width	Item Number
30"	1320098
36"	1420098
42"	1820098
48"	1520098
1000mm	1720098



Width	Item Number
30"	1320091
36"	1420091
42"	1820091
48"	1520091
1000mm	1720091



Width	Item Number
30"	1320099
36"	1420099
42"	1820099
48"	1520099
1000mm	1720099



Width	Right Hand	Left Hand
30"	1320094	1320092
36"	1420094	1420092
42"	1820094	1820092
48"	1520094	1520092
1000mm	1720094	1720092

- Left hand door/drawer unit shown
- Right hand reversed





BASE CABINETS

28-3/4" HIGH X 22" DEEP

Drawer Units



Width	Item Number
30"	1320095
36"	1420095
42"	1820095
48"	1520095
1000mm	1720095



Width	Item Number
30"	1320B44
36"	1420B44
42"	1820B44
48"	1520B44
1000mm	1720B44

File Drawer Units







30" and wider

Width	Item Number	
12"	1020122	
15"	1920122	
18"	1120122	
24"	1220122	
30"	1320122	
36"	1420122	
42"	1820122	
48"	1520122	
513mm	1620122	
1000mm	1720122	



up to 24" wide



30" and wider

Item Number
1020177
1920177
1120177
1220177
1320177
1420177
1820177
1520177
1620177
1720177



up to 24" wide



30" and wider

Width	Item Number
12"	1020188
15"	1920188
18"	1120188
24"	1220188
30"	1320188
36"	1420188
42"	1820188
48"	1520188
513mm	1620188
1000mm	1720188





BASE CABINETS

28-3/4" HIGH X 22" DEEP

Waste/Storage Bin Units





up to 24" wide

30" and wider







Width	Tilt Out Units	Pull Out Units	Bin Capacity
15"	1920500	1920511	1
18"	1120500	1120511	1
24"	1220500	1220511	1
30"	1320500	1320511	2
36"	1420500	1420511	2
42"	1820500	1820511	2
48"	1520500	1520511	2
513mm	1620500	1620511	1
1000mm	1720500	1720511	2

- · Add option code 84 to order stainless steel waste/storage bin, or add option code 85 for PVC waste/storage bin for Tilt Out units
- Tilt out units are hinged at the bottom
- Pull out units include full extension runners


Sliding Door Units



Width	Item Number
30"	1320711
36"	1420711
42"	1820711
48"	1520711
1000mm	1720711

Blind Corner Units



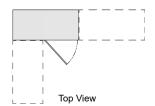


Width	Door on Right	Door on Left
30"	1320311	1320411
36"	1420311	1420411
42"	1820311	1820411
48"	1520311	1520411
1000mm	1720311	1720411

- · Right door blind corner unit shown, left door reversed
- · Center stile located at the middle point

Diagonal Corner Units



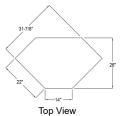

Description	Right Hinged	Left Hinged	
Full Door Unit	1027011	1027111	
False Panel Unit	1027422	1027922	

- · Left hinged, full door corner unit shown
- · Add option code 89 to order a Lazy Susan





BASE CABINETS

28-3/4" HIGH X 22" DEEP

Sink Units





up to 24" wide



up to 24" wide

30" and wider

Width **Right Hinged** Left Hinged # of Doors 12" 1023011 1023211 1 15" 1923011 1923211 18" 1123011 1123211 1 24" 1223011 1223211 30" 1323011 2 2 36" 1423011 42" 1823011 2 48" 1523011 2 513mm 1623011 1623211 1 1000mm 2 1723011

Width	Right Hinged	Left Hinged	# of Doors
12"	1028022	1028222	1
15"	1928022	1928222	1
18"	1128022	1128222	1
24"	1228022	1228222	1
30"	1328022	-	2
36"	1428022	-	2
42"	1828022	-	2
48"	1528022	-	2
513mm	1628022	1628222	1
1000mm	1728022	-	2

- · Add option code 90 to order Sink Hanger Supports
- Sink units are constructed with satin coated (galvanized) bottom, powder coated to match cabinet

End Sink Units



Sink Support Kit



Description	54" Wide	58" Wide
Full Door Unit	9900019	9900020
False Panel I Init	9900023	9900024

- · Full door end sink unit shown
- Add option code 90 to order Sink Hanger Supports
- · Sink units are constructed with satin coated (galvanized) bottom, powder coated to match cabinet

Width	Item Number
12"	SSK1012
15"	SSK1015
18"	SSK1018
24"	SSK1024
30"	SSK1030
36"	SSK1036
42"	SSK1042
48"	SSK1048
54"	SSK1054
58"	SSK1058
Corner Unit	SSK0000

- Metal frame to support and position an undermount sink
- Sink Support Kit for Corner Unit is 26" deep and not exactly as shown





BASE CABINETS

31-3/4" HIGH X 22" DEEP

Open Units



Full Door Units







30" and wider

Door/Drawer Units



up to 24" wide



30" and wider



-		-
]	-
		-
		-

Width	Right Hand	Left Hand
30"	1380942	1380924
36"	1480942	1480924
42"	1880942	1880924
48"	1580942	1580924
1000mm	1780942	1780924

- Left hand door/drawer unit shown
- · Right hand reversed

Width	Item Number
12"	1080000
15"	1980000
18"	1180000
24"	1280000
30"	1380000
36"	1480000
42"	1880000
48"	1580000
513mm	1680000
1000mm	1780000

Width	Right Hinged	Left Hinged	# of Doors
12"	1080011	1080211	1
15"	1980011	1980211	1
18"	1180011	1180211	1
24"	1280011	1280211	1
30"	1380011	-	2
36"	1480011	-	2
42"	1880011	-	2
48"	1580011	-	2
513mm	1680011	1680211	1
1000mm	1780011	-	2

Width	Right Hinged	Left Hinged	# of Doors
12"	1080022	1080222	1
15"	1980022	1980222	1
18"	1180022	1180222	1
24"	1280022	1280222	1
30"	1380022	-	2
36"	1480022	-	2
42"	1880022	-	2
48"	1580022	-	2
513mm	1680022	1680222	1
1000mm	1780022	-	2

Width	Item Number	
30"	1380093	
36"	1480093	
42"	1880093	
48"	1580093	
1000mm	1780093	



Width	Right Hand	Left Hand
30"	1380042	1380024
36"	1480042	1480024
42"	1880042	1880024
48"	1580042	1580024
1000mm	1780042	1780024

- Left hand door/drawer unit shown
- · Right hand reversed





BASE CABINETS

31-3/4" HIGH X 22" DEEP

Drawer Units



up to 24" wide



30" and wider



up to 24" wide



30" and wider



up to 24" wide



30" and wider



Width	Item Number
30"	1380090
36"	1480090
42"	1880090
48"	1580090
1000mm	1780090



Width	Item Number
30"	1380099
36"	1480099
42"	1880099
48"	1580099
1000mm	1780099

Width	Item Number
15"	1980090
18"	1180090
24"	1280090
30"	13800NN
36"	14800NN
42"	18800NN
48"	15800NN
513mm	1680090
1000mm	17800NN

Width	Item Number
12"	1080044
15"	1980044
18"	1180044
24"	1280044
30"	1380044
36"	1480044
42"	1880044
48"	1580044
513mm	1680044
1000mm	1780044

Width	Item Number
12"	1080055
15"	1980055
18"	1180055
24"	1280055
30"	1380055
36"	1480055
42"	1880055
48"	1580055
513mm	1680055
1000mm	1780055



Width	Item Number
30"	1380097
36"	1480097
42"	1880097
48"	1580097
1000mm	1780097





Width

15"

18"

24"

30"

36"

42"

48"

513mm

Width

15"

18"

24"

30"

36"

42"

48"

513mm

1000mm

1000mm

Item Number

1980100

1180100

1280100

1380100

1480100

1880100

1580100

1680100

1780100

1980111

1180111

1280111

1380111

1480111

1880111

1580111

1680111

1780111

Item Number

QUALITY BY DESIGN

BASE CABINETS

31-3/4" HIGH X 22" DEEP

File Drawer Units



up to 24" wide

30" and wider





up to 24" wide

30" and wider

Sliding Door Units



Width	Item Number
30"	1380711
36"	1480711
42"	1880711
48"	1580711
1000mm	1780711

Waste/Storage Bin Units





up to 24" wide

30" and wider



optional bin)





Pull Out Unit

Width	Tilt Out Units	Pull Out Units	Bin Capacity
15"	1980500	1980511	1
18"	1180500	1180511	1
24"	1280500	1280511	1
30"	1380500	1380511	2
36"	1480500	1480511	2
42"	1880500	1880511	2
48"	1580500	1580511	2
513mm	1680500	1680511	1
1000mm	1780500	1780511	2

- · Add option code 84 to order stainless steel waste/storage bin, or add option code 85 for PVC waste/storage bin for Tilt Out units
- Tilt out units are hinged at the bottom
- Pull out units include full extension runners





48"

513mm

1000mm

1000mm

1583011

1683011

1783011

QUALITY BY DESIGN

BASE CABINETS

31-3/4" HIGH X 22" DEEP

2

1

2

Sink Units





up to 24" wide

30" and wider






up to 24" wide

30" and wider

End Sink Units



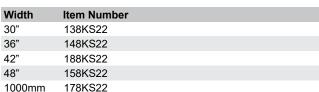

Pocket Door Sink Units







Base bottom in down position

Base bottom in knee space position

Width	Right Hinged	Left Hinged	# of Doors
12"	1083011	1083211	1
15"	1983011	1983211	1
18"	1183011	1183211	1
24"	1283011	1283211	1
30"	1383011	-	2
36"	1483011	-	2
42"	1883011	-	2

1683211

Width	Right Hinged	Left Hinged	# of Doors	
12"	1088022	1088222	1	
15"	1988022	1988222	1	
18"	1188022	1188222	1	
24"	1288022	1288222	1	
30"	1388022	-	2	
36"	1488022	-	2	
42"	1888022	-	2	
48"	1588022	-	2	
513mm	1688022	1688222	1	

- · Add option code 90 to order Sink Hanger Supports
- Sink units are constructed with satin coat (galvanized) inner horizontal bottom base plate as standard
- Description 54" Wide 58" Wide 9980018 Full Door Unit 9980017 9980021 9980022 False Panel Unit
- · Full door end sink unit shown

1788022

- · Add option code 90 to order Sink Hanger Supports
- · Sink units are constructed with satin coated (galvanized) bottom, powder coated to match cabinet
- · Incorporates pocket doors and fold-up base to convert from sink cabinet design to knee space design
- Sink units are constructed with satin coated (galvanized) bottom, powder coated to match cabinet
- · 42" wide cabinets and larger are ADA compliant






BASE CABINETS

31-3/4" HIGH X 22" DEEP

Diagonal Corner Sink Units




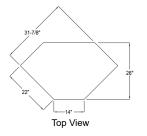

Description	Right Hinged	Left Hinged
Full Door Unit	1087211	1087511
False Panel Unit	1087722	1087622

- · Left hinged, full door corner sink unit shown
- Add option code 90 to order Sink Hanger Supports
- Sink units are constructed with satin coated (galvanized) bottom, powder coated to match cabinet

Knee Space Sink Unit Dimensions

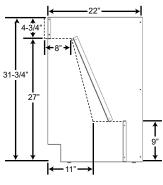

- Gable Legs are 1-1/2" wide
- · Knee Space ships unassembled
- · Accommodates a 6" deep sink
- 30" wide units with gable(s) are not ADA compliant
- · See dotted line in figure to the right for ADA requirements

Knee Space Sink Unit



Width	Item Number
30"	13K3011
33"	10K3011
36"	14K3011
42"	18K3011
48"	15K3011
1000mm	17K3011

Center Panel Only



Width	Item Number
30"	13KP011
33"	10KP011
36"	14KP011
42"	18KP011
48"	15KP011
1000mm	17KP011

Left Leg With Center Panel



Width	Item Number
30"	13L3011
33"	10L3011
36"	14L3011
42"	18L3011
48"	15L3011
1000mm	17L3011

Right Leg With Center Panel



Width	Item Number
30"	13KR011
33"	10KR011
36"	14KR011
42"	18KR011
48"	15KR011
1000mm	17KR011





MOBILE BASE CABINETS

Mobile Cabinets



Our steel cabinets can be fitted with your choice of door and drawer fronts made from painted steel, stainless steel, wood, plastic laminate or phenolic resin in inset or architectural overlay design.

Heavy duty locking 3-1/2" twin wheel casters provide a total load rating (including mobile cabinet and work top) of 500lbs. Tip resisting counter weights are provided to aid in safe operation. Resulting mobility allows the user to quickly and easily move to another workstation, as tasks change. Use option code 77 for casters on cabinets.

OPTIONS

See Steel Specification & Options section for information on the following options:

- · Door and drawer handles
- Door hinges
- · Drawer body and suspension
- · Lock and key options





BASE CABINET OPTIONS

Option	Description
00	Modified door/drawer configuration
02	Label holders on all doors/drawers
03	Numeric or alpha identification disks on all doors/drawers (customer details required)
04	Standard louver configuration on doors (6 louvers at 3" wide x 5-1/2" high)
05	Non-standard louver configuration on doors (customer drawing required)
07	Cut-outs for electrical and/or plumbing services (customer drawing required)
80	Security proof panel(s)
09	Stainless steel waste flap
10	One fixed drawer divider per drawer (centered in drawer)
11	Multiple drawer dividers per drawer (customer drawing required)
12	Pull out writing board * - cold rolled steel
13	Pull out writing board * - cold rolled steel & high pressure laminate writing surface
14	Ventilation grill in toe kick base
15	Finished back on cabinet
16	Ship handles not installed
17	Nylon levelers
17C	Nylon levelers with counter weight
18	Horizontally mounted handles on doors
19	Full width drawer on 5" high apron
26	200lb full extension drawer runners
40	One extra cabinet shelf
41	Two extra cabinet shelves
42	One interior cabinet shelf on full extension runners (replaces fixed shelf)
43	Two interior cabinet shelves on full extension runners (replaces fixed shelves)
44	Three interior cabinet shelves on full extension runners (replaces fixed shelves)
48	Keyboard tray (replaces drawer)
50	304 stainless steel drawer bodies in power coated cold rolled steel cabinet
51	Component made from 304 stainless steel
52	Component made from 316 stainless steel
53	Galvenneal steel cabinet construction (Sigma Barrier™)
54	304 stainless steel drawer bodies in galvenneal steel cabinet (Sigma Barrier™)
56	Other customer specified material
59	Customer specified runners (customer details required)
70	Front toe kick
71	Toe kick - front & left side
72	Toe kick - front & right side
70	Toe kick - front & 2 sides
73	

Note: * 30" wide cabinets and larger receive two pull out writing boards, one pull out board is provided on full width drawer cabinets.

Option	Description
74A	Toe kick - 1-1/2" on back
75	Toe kick - front, back & two sides
76	Replace toe kick with four square legs
77	Replace toe kick with four casters
78	Replace toe kick with customer specified casters
79	Other customer specified base or toe kick
84	Stainless steel bin insert
85	PVC bin insert
89	Lazy Susan storage unit for corner cabinets
90	Sink supports - for under counter mounted sinks
AD	Add - a - drawer design (6" high construction, applies to top drawer)
BF	Replace toe kick with 4" bullet feet
C25	Replace toe kick with 250lb casters
C35	Replace toe kick with 350lb casters
DC2	Magnetic door catch
F1	Welded finished top
F2	Welded finished top & back
F26	200lb full extension runners on file drawers
FCW	Front corners of doors and drawers welded
I01	100lb full extension runners with interlock
126	200lb full extension runners with interlock
LA	Bilingual label (French/English)
LF	Longitudinal file hangers on file drawers
NB	No base, cabinet to sit on floor, levelers included
NC	No base, cabinet to sit on 4" high curb, levelers not included
NF	No cabinet floor and toe kick is integral with door
OL	Overlay design
PC	Foot pedal cutouts centered
PDC	Positive door catch (standard on overlay cabinets)
PFE	Pull out writing board * on full extension runners
PL	Foot pedal cutouts left hand side
PR	Foot pedal cutouts right hand side
PW	Pre-wired (customer details required)
RB	Removable back on drawer cabinet
SK	Skid item
S01	100lb full extension soft close drawer runners
S26	200lb full extension soft close drawer runners
T1	Replace one shelf with one pull out tray on full extension runners with 2" high lip
T2	Replace two shelves with two pull out trays on full extension runners with 2" high lip
Т3	Replace three shelves with three pull out trays on full extension runners with 2" high lip
WF	Wood fronts - see specifications
wo	Wood fronts provided and installed by others





24" HIGH X 14" DEEP

Open Units



Width	Item Number
15"	5770010
18"	5170010
24"	5270010
30"	5370010
36"	5470010
42"	5670010
48"	5570010
513mm	5970010
1000mm	5870010

Hinged Framed Glass Door Units





30" and wider

Width	Right Hinged	Left Hinged	# of Doors
15"	5770030	5771030	1
18"	5170030	5171030	1
24"	5270030	5271030	1
30"	5370030	-	2
36"	5470030	-	2
42"	5670030	-	2
48"	5570030	-	2
513mm	5970030	5971030	1
1000mm	5870030	-	2

Hinged Door Units





Width	Right Hinged	Left Hinged	# of Doors
15"	5770060	5771060	1
18"	5170060	5171060	1
24"	5270060	5271060	1
30"	5370060	-	2
36"	5470060	-	2
42"	5670060	-	2
48"	5570060	-	2
513mm	5970060	5971060	1
1000mm	5870060	-	2

Sliding Glass Door Units



Width	Item Number
30"	5370020
36"	5470020
42"	5670020
48"	5570020
1000mm	5870020

Sliding Door Units



Width	Item Number
30"	5370050
36"	5470050
42"	5670050
48"	5570050
1000mm	5870050

Sliding Framed Glass Door Units



Width	Item Number
30"	5370040
36"	5470040
42"	5670040
48"	5570040
1000mm	5870040





24" HIGH X 14" DEEP

Hinged Framed Glass Door

Door on Left

5370430

5470430

5670430

5570430

5870430

Blind Corner Units

Door on Right

5370330

5470330


5670330

5570330

5870330

Blind Corner Units



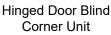




Hinged Framed Glass Door Blind Corner Unit



- · Right door corner units shown, left door reversed
- · Blind portion of cabinet is 14" wide

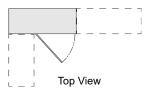

WALL CABINETS

30" HIGH X 14" DEEP

Open Units



Hinged Framed Glass Door Units



up to 24" wide



30" and wider

Hinged Door Units



up to 24" wide



Width	Item Number
15"	5710010
18"	5110010
24"	5210010
30"	5310010
36"	5410010
42"	5610010
48"	5510010
513mm	5910010
1000mm	5810010

Width	Right Hinged	Left Hinged	# of Doors
15"	5710030	5711030	1
18"	5110030	5111030	1
24"	5210030	5211030	1
30"	5310030	-	2
36"	5410030	-	2
42"	5610030	-	2
48"	5510030	-	2
513mm	5910030	5911030	1
1000mm	5810030	-	2

Width	Right Hinged	Left Hinged	# of Doors
15"	5710060	5711060	1
18"	5110060	5111060	1
24"	5210060	5211060	1
30"	5310060	-	2
36"	5410060	-	2
42"	5610060	-	2
48"	5510060	-	2
513mm	5910060	5911060	1
1000mm	5810060	-	2





30" HIGH X 14" DEEP

Sliding Glass Door Units



Width	Item Number
30"	5310020
36"	5410020
42"	5610020
48"	5510020
1000mm	5810020

Sliding Framed Glass Door Units



Width	Item Number
30"	5310040
36"	5410040
42"	5610040
48"	5510040
1000mm	5810040


Sliding Door Units



Width	Item Number
30"	5310050
36"	5410050
42"	5610050
48"	5510050
1000mm	5810050

Blind Corner Units



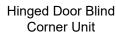




Hinged Framed Glass Door Blind Corner Unit

Width	Hinged Door Blind Corner Units		Hinged Framed Glass Door Blind Corner Units	
	Door on Right	Door on Left	Door on Right	Door on Left
30"	5310360	5310460	5310330	5310430
36"	5410360	5410460	5410330	5410430
42"	5610360	5610460	5610330	5610430
48"	5510360	5510460	5510330	5510430
1000mm	5810360	5810460	5810330	5810430

- · Right door corner units shown, left door reversed
- Blind portion of cabinet is 14" wide

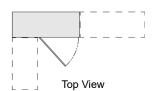





36" HIGH X 14" DEEP

Open Units



Hinged Framed Glass Door Units





up to 24" wide

30" and wider

Hinged Door Units





up to 24" wide

30" and wider

Sliding Glass Door Units



Width	Item Number
30"	5320020
36"	5420020
42"	5620020
48"	5520020
1000mm	5820020

Sliding Door Units



Width	Item Number
30"	5320050
36"	5420050
42"	5620050
48"	5520050
1000mm	5820050

Width **Item Number** 15" 5720010 18" 5120010 24" 5220010 30" 5320010 36" 5420010 42" 5620010 48" 5520010 513mm 5920010 1000mm 5820010

Width	Right Hinged	Left Hinged	# of Doors
15"	5720030	5721030	1
18"	5120030	5121030	1
24"	5220030	5221030	1
30"	5320030	-	2
36"	5420030	-	2
42"	5620030	-	2
48"	5520030	-	2
513mm	5920030	5921030	1
1000mm	5820030	-	2

Width	Right Hinged	Left Hinged	# of Doors
15"	5720060	5721060	1
18"	5120060	5121060	1
24"	5220060	5221060	1
30"	5320060	-	2
36"	5420060	-	2
42"	5620060	-	2
48"	5520060	-	2
513mm	5920060	5921060	1
1000mm	5820060	-	2

Sliding Framed Glass Door Units



Width	Item Number
30"	5320040
36"	5420040
42"	5620040
48"	5520040
1000mm	5820040





36" HIGH X 14" DEEP

Blind Corner Units





Hinged Door Blind Corner Unit

Hinged Framed

Glass Door Blind Corner Unit

Hinged Door Width Blind Corner Units		nits	Hinged Framed Glass Door Blind Corner Units	
	Door on Right	Door on Left	Door on Right	Door on Left
30"	5320360	5320460	5320330	5320430
36"	5420360	5420460	5420330	5420430
42"	5620360	5620460	5620330	5620430
48"	5520360	5520460	5520330	5520430
1000mm	5820360	5820460	5820330	5820430
Right door corner units shown left door reversed				

- Right door corner units shown, left door reversed
- Blind portion of cabinet is 14" wide

		_	
—			
•			
Top V	'iew		

WALL CABINETS

48" HIGH X 14" DEEP

Open Units



Width	Item Number
15"	5730010
18"	5130010
24"	5230010
30"	5330010
36"	5430010
42"	5630010
48"	5530010
513mm	5930010
1000mm	5830010

Hinged Framed Glass Door Units





up to 24" wide

30" and wider

Width	Right Hinged	Left Hinged	# of Doors
15"	5730030	5731030	1
18"	5130030	5131030	1
24"	5230030	5231030	1
30"	5330030	-	2
36"	5430030	-	2
42"	5630030	-	2
48"	5530030	-	2
513mm	5930030	5931030	1
1000mm	5830030	-	2





48" HIGH X 14" DEEP

Hinged Door Units





Width	Right Hinged	Left Hinged	# of Doors
15"	5730060	5731060	1
18"	5130060	5131060	1
24"	5230060	5231060	1
30"	5330060	-	2
36"	5430060	-	2
42"	5630060	-	2
48"	5530060	-	2
513mm	5930060	5931060	1
1000mm	5830060	-	2

Sliding Glass Door Units



Width	Item Number
30"	5330020
36"	5430020
42"	5630020
48"	5530020
1000mm	5830020

Sliding Framed Glass Door Units



Width	Item Number
30"	5330040
36"	5430040
42"	5630040
48"	5530040
1000mm	5830040

Sliding Door Units



Width	Item Number
30"	5330050
36"	5430050
42"	5630050
48"	5530050
1000mm	5830050

Blind Corner Units








Hinged Framed Glass Door Blind Corner Unit

Width	Hinged Door Blind Corner Units		Hinged Framed Glass Door Blind Corner Units	
	Door on Right	Door on Left	Door on Right	Door on Left
30"	5330360	5330460	5330330	5330430
36"	5430360	5430460	5430330	5430430
42"	5630360	5630460	5630330	5630430
48"	5530360	5530460	5530330	5530430
1000mm	5830360	5830460	5830330	5830430

- Right door corner units shown, left door reversed
- · Blind portion of cabinet is 14" wide

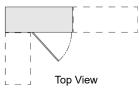





24" HIGH X 18" DEEP

Open Units



Hinged Framed Glass Door Units





up to 24" wide

Hinged Door Units





up to 24" wide

Sliding Glass Door Units



Width	Item Number
30"	5380020
36"	5480020
42"	5680020
48"	5580020
1000mm	5880020

Sliding Door Units



Width	Item Number
30"	5380050
36"	5480050
42"	5680050
48"	5580050
1000mm	5880050

Width **Item Number** 15" 5780010 18" 5180010 24" 5280010 30" 5380010 36" 5480010 42" 5680010 48" 5580010 513mm 5980010 1000mm 5880010

Width	Right Hinged	Left Hinged	# of Doors
15"	5780030	5781030	1
18"	5180030	5181030	1
24"	5280030	5281030	1
30"	5380030	-	2
36"	5480030	-	2
42"	5680030	-	2
48"	5580030	-	2
513mm	5980030	5981030	1
1000mm	5880030	-	2

Width	Right Hinged	Left Hinged	# of Doors
15"	5780060	5781060	1
18"	5180060	5181060	1
24"	5280060	5281060	1
30"	5380060	-	2
36"	5480060	-	2
42"	5680060	-	2
48"	5580060	-	2
513mm	5980060	5981060	1
1000mm	5880060	-	2

Sliding Framed Glass Door Units



Width	Item Number
30"	5380040
36"	5480040
42"	5680040
48"	5580040
1000mm	5880040





24" HIGH X 18" DEEP

Blind Corner Units






Hinged Door Blind Corner Unit

Hinged Framed Glass Door Blind Corner Unit

Width	Hinged Door Blind Corner U	nged Door nd Corner Units		Hinged Framed Glass Door Blind Corner Units	
	Door on Right	Door on Left	Door on Right	Door on Left	
30"	5380360	5380460	5380330	5380430	
36"	5480360	5480460	5480330	5480430	
42"	5680360	5680460	5680330	5680430	
48"	5580360	5580460	5580330	5580430	
1000mm	5880360	5880460	5880330	5880430	

- · Right door corner units shown, left door reversed
- Blind portion of cabinet is 18" wide

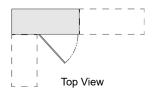

WALL CABINETS

30" HIGH X 18" DEEP

Open Units



Width	Item Number
15"	5740010
18"	5140010
24"	5240010
30"	5340010
36"	5440010
42"	5640010
48"	5540010
513mm	5940010
1000mm	5840010

Hinged Framed Glass Door Units





up to 24" wide 30" and wider

Width	Right Hinged	Left Hinged	# of Doors
15"	5740030	5741030	1
18"	5140030	5141030	1
24"	5240030	5241030	1
30"	5340030	-	2
36"	5440030	-	2
42"	5640030	-	2
48"	5540030	-	2
513mm	5940030	5941030	1
1000mm	5840030	-	2

Hinged Door Units





Width	Right Hinged	Left Hinged	# of Doors
15"	5740060	5741060	1
18"	5140060	5141060	1
24"	5240060	5241060	1
30"	5340060	-	2
36"	5440060	-	2
42"	5640060	-	2
48"	5540060	-	2
513mm	5940060	5941060	1
1000mm	5840060	-	2





30" HIGH X 18" DEEP

Sliding Glass Door Units



Width	Item Number
30"	5340020
36"	5440020
42"	5640020
48"	5540020
1000mm	5840020

Sliding Framed Glass Door Units



Width	Item Number
30"	5340040
36"	5440040
42"	5640040
48"	5540040
1000mm	5840040


Sliding Door Units



Width	Item Number
30"	5340050
36"	5440050
42"	5640050
48"	5540050
1000mm	5840050

Blind Corner Units



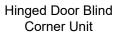


Hinged Framed Glass Door Blind Corner Unit



Hinged Door Blind Corner Units		Hinged Framed Glass Door Blind Corner Units	
Door on Right	Door on Left	Door on Right	Door on Left
5340360	5340460	5340330	5340430
5440360	5440460	5440330	5440430
5640360	5640460	5640330	5640430
5540360	5540460	5540330	5540430
5840360	5840460	5840330	5840430
	Blind Corner U Door on Right 5340360 5440360 5640360 5540360 5840360	Blind Corner UnitsDoor on RightDoor on Left5340360534046054403605440460564036056404605540360554046058403605840460	Blind Corner Units Blind Corner Units Door on Right Door on Left Door on Right 5340360 5340460 5340330 5440360 5440460 5440330 5640360 5640460 5640330 5540360 5540460 5540330

- · Right door corner units shown, left door reversed
- Blind portion of cabinet is 18" wide

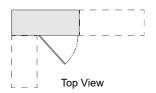





36" HIGH X 18" DEEP

Open Units



Width	Item Number
15"	5750010
18"	5150010
24"	5250010
30"	5350010
36"	5450010
42"	5650010
48"	5550010
513mm	5950010
1000mm	5850010

Hinged Framed Glass Door Units





up to 24" wide

Width	Right Hinged	Left Hinged	# of Doors
15"	5750030	5751030	1
18"	5150030	5151030	1
24"	5250030	5251030	1
30"	5350030	-	2
36"	5450030	-	2
42"	5650030	-	2
48"	5550030	-	2
513mm	5950030	5951030	1
1000mm	5850030	-	2

Left Hinged

5751060

5151060

5251060

5951060

of Doors

1

1

1

2

2

2

1

2

Hinged Door Units





Sliding Glass Door Units

1	

Width	Item Number
30"	5350020
36"	5450020
42"	5650020
48"	5550020
1000mm	5850020

Sliding Framed Glass Door Units

Width

15"

18"

24"

30"

36"

42"

48"

513mm

1000mm

Right Hinged

5750060

5150060

5250060

5350060

5450060

5650060

5550060

5950060

5850060



Width	Item Number
30"	5350040
36"	5450040
42"	5650040
48"	5550040
1000mm	5850040

Sliding Door Units



Width	Item Number
30"	5350050
36"	5450050
42"	5650050
48"	5550050
1000mm	5850050





QUALITY BY DESIGN

WALL CABINETS

36" HIGH X 18" DEEP

Blind Corner Units





Hinged Framed Glass Door Blind Corner Unit

Hinged Door Blind Corner Units		Hinged Framed Glass Door Blind Corner Units	
Door on Right	Door on Left	Door on Right	Door on Left
5350360	5350460	5350330	5350430
5450360	5450460	5450330	5450430
5650360	5650460	5650330	5650430
5550360	5550460	5550330	5550430
5850360	5850460	5850330	5850430
	Blind Corner U Door on Right 5350360 5450360 5650360 5550360 5850360	Blind Corner UnitsDoor on RightDoor on Left53503605350460545036054504605650360565046058503605850460	Blind Corner Units Blind Corner Units Door on Right Door on Left Door on Right 5350360 5350460 5350330 5450360 5450460 5450330 5650360 5650460 5650330 5550360 5550460 5550330

- Right door corner units shown, left door reversed
- Blind portion of cabinet is 18" wide

Top V	iew	

WALL CABINETS

Open Units

Corner Unit



Width	Item Number
15"	5760010
18"	5160010
24"	5260010
30"	5360010
36"	5460010
42"	5660010
48"	5560010
513mm	5960010
1000mm	5860010

Hinged Framed Glass Door Units





30" and wider

Width	Right Hinged	Left Hinged	# of Doors
15"	5760030	5761030	1
18"	5160030	5161030	1
24"	5260030	5261030	1
30"	5360030	-	2
36"	5460030	-	2
42"	5660030	-	2
48"	5560030	-	2
513mm	5960030	5961030	1
1000mm	5860030	_	2





48" HIGH X 18" DEEP

Hinged Door Units





Width	Right Hinged	Left Hinged	# of Doors
15"	5760060	5761060	1
18"	5160060	5161060	1
24"	5260060	5261060	1
30"	5360060	-	2
36"	5460060	-	2
42"	5660060	-	2
48"	5560060	-	2
513mm	5960060	5961060	1
1000mm	5860060	-	2

Sliding Glass Door Units



Width	Item Number
30"	5360020
36"	5460020
42"	5660020
48"	5560020
1000mm	5860020

Sliding Framed Glass Door Units



Width	Item Number
30"	5360040
36"	5460040
42"	5660040
48"	5560040
1000mm	5860040

Sliding Door Units



Width	Item Number
30"	5360050
36"	5460050
42"	5660050
48"	5560050
1000mm	5860050

Blind Corner Units








Hinged Framed Glass Door Blind Corner Unit

Width	Hinged Door Blind Corner Units		Hinged Framed Glass Door Blind Corner Units	
	Door on Right	Door on Left	Door on Right	Door on Left
30"	5360360	5360460	5360330	5360430
36"	5460360	5460460	5460330	5460430
42"	5660360	5660460	5660330	5660430
48"	5560360	5560460	5560330	5560430
1000mm	5860360	5860460	5860330	5860430

- Right door corner units shown, left door reversed
- Blind portion of cabinet is 18" wide

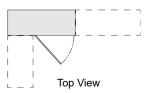





DIAGONAL CORNER WALL CABINETS

14" DEEP & 18" DEEP


Hinged Framed Glass Corner Door Units




Hinged Door Corner Units



Top View 14" Deep

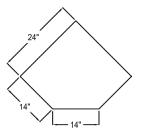

Top View 18" Deep

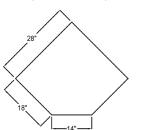

Height	Right Hinged	Left Hinged	Depth
24"	5070030	5071030	14"
30"	5010030	5011030	14"
36"	5020030	5021030	14"
48"	5030030	5031030	14"
24"	5080030	5081030	18"
30"	5040030	5041030	18"
36"	5050030	5051030	18"
48"	5060030	5061030	18"

Height	Right Hinged	Left Hinged	Depth
24"	5070060	5071060	14"
30"	5010060	5011060	14"
36"	5020060	5021060	14"
48"	5030060	5031060	14"
24"	5080060	5081060	18"
30"	5040060	5041060	18"
36"	5050060	5051060	18"
48"	5060060	5061060	18"

- · Left hinged corner door units shown, right door is reversed
- 24" high units receive one shelf
- 30" high units receive two shelves
- 36" high units receive two shelves
- 48" high units receive three shelves





Flipper Door Units



Width	Item Number	Depth
30"	53A0300	14"
36"	54A0300	14"
48"	55A0300	14"
30"	53B0300	18"
36"	54B0300	18"
48"	55B0300	18"

- Flipper door units are 16-1/2" high (no shelf is provided)
- · The flipper door on these units swings up and into the cabinet body
- · Order this cabinet to match width of support structure frame
- Add option code FR to suspend this wall cabinet from Sigma Cart[™] and Sigma Flex[™]
- Add option code HR to suspend this wall cabinet from Sigma Frame™


Wall Cabinet "J" Mounting Rail



Width	Item Number
24"	WHR0024
30"	WHR0030
36"	WHR0036
42"	WHR0042
48"	WHR0048
60"	WHR0060
72"	WHR0072

- A "J" Mounting Rail is required to hang metal above-counter storage cabinets on to Sigma Flex™ or Sigma Frame™ products
- Order "J" rails to match width of support structure frame and Wall Cabinet Hanger Rail
- "J" Rail ships loose


"Z" Wall Mounting Rail

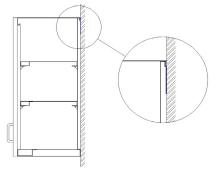



- Add option code "ZR" to any wall cabinet to have a "Z" Mounting Rail shipped with the wall cabinet
- The "Z" Mounting Rail is 16 gauge zinc coated steel and mounts to the wall to easily install and align wall cabinets
- · Cabinets must still be anchored to wall to support rated load
- "Z" Wall Mounting Rail ships loose

Wall Cabinet Hanging Rails

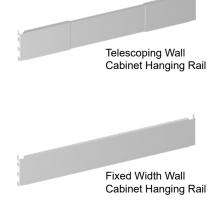

	Sigma Flex™		Sigma F	rame™
Width	Telescoping Rail	Fixed Width Rail	Telescoping Rail	Fixed Width Rail
24"	TWR0024	FWR0024	AWR0024	SWR0024
30"	TWR0030	FWR0030	AWR0030	SWR0030
36"	TWR0036	FWR0036	AWR0036	SWR0036
42"	TWR0042	FWR0042	AWR0042	SWR0042
48"	TWR0048	FWR0048	AWR0048	SWR0048
54"	TWR0054	FWR0054	AWR0054	SWR0054
60"	TWR0060	FWR0060	AWR0060	SWR0060
72"	TWR0072	FWR0072	AWR0072	SWR0072

- Wall Cabinet Hanger Rails are required to hang above-counter storage cabinets (must be customized to hang on cores with panels)
- Order rails to match width of support structure frame (Telescoping Rails are not suitable for cabinets smaller than 24" wide)
- If ordering rails larger than the wall cabinet, two rails are required for support.
- Telescoping Rail is adjustable up to 5"
- Wall cabinets must be ordered with option code "FR" to suspend from Sigma Flex™ and Sigma Cart™ and option code "HR" to suspend from Sigma Frame™





WALL CABINET OPTIONS

Option	Description
00	Modified door/drawer configuration
02	Label holders on all doors/drawers
03	Numeric or alpha identification disks on all doors (customer details required)
07	Cut-outs for electrical and/or plumbing services (customer drawing required)
15	Finished back
16	Ship handles not installed
18	Horizontally mounted handles on doors
40	One extra cabinet shelf
41	Two extra cabinet shelves
45	Duolite safety glass (replaces float glass)
46	Tempered glass (replaces float glass)
47	Customer specified glass (replaces float glass)
4A	Glass tint 44% (Grey Smoke) (replaces float glass)
51	Component made from 304 stainless steel
52	Component made from 316 stainless steel
53	Galvenneal steel cabinet construction (Sigma Barrier™)
56	Other customer specified material
59	Customer specified runners (customer details required)
DC2	Magnetic door catch
DS	Double sided wall cabinet
FM	Wall case modified to mount on Sigma Fusion™
FCW	Front corners of doors and drawers welded
FR	Special mounting provisions for wall cabinets to suspended from Sigma Flex™
HR	Special mounting provisions for wall cabinets to suspended from Sigma Frame™
OL	Overlay design
PDC	Positive door catch (standard on overlay cabinets)
SK	Skid item
SL	Seismic lip
WF	Wood fronts - see specifications
WO	Wood fronts provided and installed by others
ZR	"Z" mounting rail

See Steel Specification & Options section for information on the following options:

- Door and drawer handles
- · Door hinges
- Drawer body and suspension
- Lock and key options





FLOOR CABINETS

84" HIGH X 14" DEEP

Open Units



Width	Item Number
18"	6110010
24"	6210010
30"	6310010
36"	6410010
42"	6610010
48"	6510010
513mm	6910010
1000mm	6810010
	1 16

Left Hinged

6110230

6210230

6910230

Fixed center shelf

Width

18"

24"

30"

36"

42"

48"

513mm

Hinged Framed Glass Door Units





up to 24" wide

30" and wider

1000mm 6810030 • Add option code 80 to order full height doors

Right Hinged

6110030

6210030

6310030

6410030

6610030

6510030

6910030

Hinged Door Units





up to 24" wide

30" and wider

Width	Right Hinged	Left Hinged	# of Doors
18"	6110060	6110260	2
24"	6210060	6210260	2
30"	6310060	-	4
36"	6410060	-	4
42"	6610060	-	4
48"	6510060	-	4
513mm	6910060	6910260	2
1000mm	6810060	-	4
A .1 .1 4! .		A Could be a South A selection of	

- · Add option code 80 to order full height doors
- Fixed center shelf

Sliding Door Units



Width

30"

36"

42"

48"

1000mm

- Add option code 80 to order full height doors
- · Fixed center shelf

Sliding Glass Door Units



	٠	Fixed	center	shelf
--	---	-------	--------	-------

er shelf	

Fixed center shelf

of Doors

2

4

4

4

4

2

4

Item Number	Width	Item Number
6310050	30"	6310020
6410050	36"	6410020
6610050	42"	6610020
6510050	48"	6510020
6810050	1000mm	6810020

Width	Item Number	
30"	6310040	
36"	6410040	
42"	6610040	
48"	6510040	
1000mm	6810040	

Sliding Framed Glass Door Units

[·] Fixed center shelf





FLOOR CABINETS

84" HIGH X 18" DEEP

Open Units



Width	Item Number
18"	6120010
24"	6220010
30"	6320010
36"	6420010
42"	6620010
48"	6520010
513mm	6920010
1000mm	6820010
	55256.5

Left Hinged

6120230

6220230

6920230

of Doors

2

4

4

4

4

2

4

· Fixed center shelf

Width

24"

30"

36"

42"

48"

513mm

1000mm

Hinged Framed Glass Door Units





up to 24" wide

30" and wider

· Add option code 80 to order full height doors

Right Hinged

6120030

6220030

6320030

6420030

6620030

6520030

6920030

6820030

Hinged Door Units





up to 24" wide

30" and wider

Width	Right Hinged	Left Hinged	# of Doors
18"	6120060	6120260	2
24"	6220060	6220260	2
30"	6320060	-	4
36"	6420060	-	4
42"	6620060	-	4
48"	6520060	-	4
513mm	6920060	6920260	2
1000mm	6820060	-	4

- · Add option code 80 to order full height doors
- · Fixed center shelf

Sliding Door Units



- Add option code 80 to order full height doors
- · Fixed center shelf

Sliding Glass Door Units

· Fixed center shelf

Sliding Framed Glass Door Units

· Fixed center shelf

Width	Item Number	
30"	6320050	
36"	6420050	
42"	6620050	
48"	6520050	
1000mm	6820050	

Width	Item Number
30"	6320020
36"	6420020
42"	6620020
48"	6520020
1000mm	6820020

Width	Item Number	
30"	6320040	
36"	6420040	
42"	6620040	
48"	6520040	
1000mm	6820040	

[·] Fixed center shelf





FLOOR CABINETS

84" HIGH X 22" DEEP

Open Units



Item Number
6130010
6230010
6330010
6430010
6630010
6530010
6930010
6830010

Left Hinged

6130230

6230230

6930230

of Doors

2

2

4

4

4

4

2

4

Right Hinged

6130030

6230030

6330030

6430030

6630030

6530030

6930030

6830030

Width

18"

24"

30"

36"

42"

48"

513mm

1000mm

Hinged Framed Glass Door Units





up to 24" wide

30" and wider

· Add option code 80 to order full height doors · Fixed center shelf

Hinged Door Units





up to 24" wide

30" and wider

width	Right Hinged	Left Hinged	# of Doors
18"	6130060	6130260	2
24"	6230060	6230260	2
30"	6330060	-	4
36"	6430060	-	4
42"	6630060	-	4
48"	6530060	-	4
513mm	6930060	6930260	2
1000mm	6830060	-	4

- · Add option code 80 to order full height doors
- · Fixed center shelf

Sliding Door Units



- Add option code 80 to order full height doors
- · Fixed center shelf

Sliding Glass Door Units



· Fixed center shelf

Sliding Framed Glass Door Units



Width	Item Number	
30"	6330040	
36"	6430040	
42"	6630040	
48"	6530040	
1000mm	6830040	

Width	Item Number
30"	6330050
36"	6430050
42"	6630050
48"	6530050
1000mm	6830050

Width	Item Number	
30"	6330020	
36"	6430020	
42"	6630020	
48"	6530020	
1000mm	6830020	

Fixed center shelf





QUALITY BY DESIGN

WARDROBE CABINETS

84" HIGH X 22" DEEP


Wardrobe Units







30" and wider

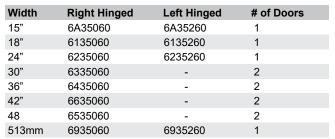

- · Wardrobe Units have solid hinged double doors and include one fixed shelf and a clothes rod
- · Units must have a minimum depth of 22" to hold hangers



	-5		
		1	
_			
		-	

		à.
	-	

Width	Item Number
30"	6335660
36"	6435660
42"	6635660
48"	6535660

- · Wardrobe Units have solid hinged double doors and includes one fixed shelf and clothes rod on the left side with one fixed top shelf and four adjustable shelves on the right side
- · Units must have a minimum depth of 22" to hold hangers

Width	Item Number
30"	6335760
36"	6435760
42"	6635760
48"	6535760


- · Wardrobe Units have solid hinged double doors and include one fixed shelf with clothes rod on the right side with one fixed top shelf and four adjustable shelves on the left side
- · Units must have a minimum depth of 22" to hold hangers





Width	Item Number	
30"	6338060	
36"	6438060	
42"	6638060	
48"	6538060	

- Tool Storage Units have solid hinged double doors and include one adjustable shelf at bottom, one pegboard on inside of each door, one pegboard on inside back and two rotating pegboards
- The pegboards on these units are metal and fit most standard tool pegs (pegs are not included)

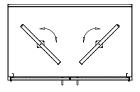

Top View - Tool Storage Units shown with rotating pegboards





INSTRUMENT STORAGE CABINETS

84" HIGH X 18" DEEP

Instrument Storage Units





30" and wider

Width	Right Hinged	Left Hinged	# of Doors
18"	61200A0	61202A0	2
24"	62200A0	62202A0	2
30"	63200A0	-	4
36"	64200A0	-	4
42"	66200A0	-	4
48"	65200A0	-	4
513mm	69200A0	69202A0	2
1000mm	68200A0	-	4

- The upper section on the Instrument Storage Units has hinged framed glass doors with three adjustable shelves
- The lower section is 34-3/4" high, with hinged solid doors, and one adjustable shelf

Instrument Storage With Pull-Out Boards



Width	Item Number
30"	63200B0
36"	64200B0
42"	66200B0
48"	65200B0
1000mm	68200B0

- The upper section on the Instrument Storage Units has sliding framed glass doors with three adjustable shelves
- The lower section is 34-3/4" high, with hinged solid doors, one adjustable shelf, and two painted steel pull-out boards with full extension runners



Width	Item Number
30"	63200C0
36"	64200C0
42"	66200C0
48"	65200C0
1000mm	68200C0

- The upper section on the Instrument Storage Units has sliding framed glass doors with three adjustable shelves
- The lower section is 34-3/4" high, with sliding solid doors, one adjustable shelf, and two painted steel pull-out boards with full extension runners

Instrument Storage Units With Compartments







Width	Right Hinged	Left Hinged	# of Doors
24"	6227030	6227230	1
36"	6427030	-	2
48"	6527030	-	2

- The 24" wide Instrument Storage Units have a hinged framed glass door with 15 fixed width compartments
- The 36" wide Instrument Storage Unit has hinged framed glass double doors with 20 fixed width compartments
- The 48" wide Instrument Storage Unit has hinged framed glass double doors with 30 fixed width compartments
- · All units include four height adjustable shelves





FLOOR CABINET OPTIONS

Option	Description
00	Modified door/drawer configuration
02	Label holders on all doors/drawers
03	Numeric or alpha identification disks on all doors/drawers (customer details required)
04	Standard louver configuration on doors (6 louvers at 3" wide x 5-1/2" long)
05	Non-standard louver configuration on doors (customer drawing required)
07	Cut-outs for electrical and/or plumbing services (customer drawing required)
80	Security proof panel
09	Stainless steel waste flap
13	Pull out writing board - cold rolled steel & high pressure laminate writing surface*
15	Finished back
16	Ship handles not installed
17	Nylon levelers
18	Horizontally mounted handles
40	One extra cabinet shelf
41	Two extra cabinet shelves
42	One interior cabinet shelf on full extension runners (replaces fixed shelf)
43	Two interior cabinet shelves on full extension runners (replaces fixed shelves)
44	Three interior cabinet shelves on full extension runners (replaces fixed shelves)
45	Duolite safety glass (replaces float glass)
46	Tempered glass (replaces float glass)
47	Customer supplied glass (replaces float glass)
4A	Glass tint 44%, Gray Smoke (replaces float glass)
51	Component made from 304 stainless steel
52	Component made from 316 stainless steel
53	Galvenneal steel cabinet construction (Sigma Barrier™)
56	Other customer specified material
59	Customer specified runners (customer details required)
70	Front toe kick
71	Toe kick - front & left side
72	Toe kick - front & right side
73	Toe kick - front & 2 sides
74	Toe kick - front & back
74A	Toe kick - 1-1/2" on back

Option	Description
75	Toe kick - front, back & two sides
76	Replace toe kick with four legs
79	Other customer specified base or toe kick
80	Full height doors
81	Coat rod for wardrobe floor cabinets
90	Sink supports - for under counter mounted sinks
BF	Replace toe kick with 4" bullet feet
C25	Replace toe kick with 250lb casters
DC2	Magnetic door catch
F26	File drawer on 200lb full extension runners
FCW	Front corners welded
101	100lb full extension runners with interlock
126	200lb full extension runners with interlock
NB	No base, cabinet to sit on floor, levelers included
NC	No base, cabinet to sit on 4" high curb, levelers not included
NF	No cabinet floor and toe kick is integral with door
OL	Overlay design
PDC	Positive door catch (standard on overlay cabinets)
SK	Skid item
S01	100lb full extension soft close drawer runners
S26	200lb full extension soft close drawer runners
T1	Replace one shelf with one pull out tray on full extension runners with 2" high lip
T2	Replace two shelf with two pull out tray on full extension runners with 2" high lip
Т3	Replace three shelf with three pull out tray on full extension runners with 2" high lip
WF	Wood fronts - see specifications
WO	Wood fronts provided and installed by others

See Steel Specification & Options section for information on the following options:

- Door and drawer handles
- · Door hinges
- · Drawer body and suspension
- · Lock and key options





ACID STORAGE CABINETS

34-3/4" HIGH X 18"/22" DEEP

See page E18 for vent hole locations.

Acid Storage Units (18" Deep)





Width	Right Hinged	Left Hinged	# of Doors	
15"	1931011	1931211	1	
18"	1131011	1131211	1	
24"	1231011	1231211	1	
30"	1331011	-	2	
36"	1431011	-	2	
42"	1831011	-	2	
48"	1531011	-	2	
513mm	1631011	1631211	1	
1000mm	1731011	-	2	

- · Polyethylene liner and phenolic resin shelf are chemically resistant
- · Add option code VK to order a Corrugated Polyolefin Vent kit for these Units





Width	Right Hinged	Left Hinged	# of Doors	
15"	1936022	1936222	1	
18"	1136022	1136222	1	
24"	1236022	1236222	1	
30"	1336022	-	2	
36"	1436022	-	2	
42"	1836022	-	2	
48"	1536022	-	2	
513mm	1636022	1636222	1	
1000mm	1736022	-	2	

- Polyethylene liner and phenolic resin shelf are chemically resistant
- · Add option code VK to order a Corrugated Polyolefin Vent kit for these Units

Acid Storage Units (22" Deep)





Width	Right Hinged	Left Hinged	# of Doors	
15"	1911011	1911211	1	
18"	1111011	1111211	1	
24"	1211011	1211211	1	
30"	1311011	-	2	
36"	1411011	-	2	
42"	1811011	-	2	
48"	1511011	-	2	
513mm	1611011	1611211	1	
1000mm	1711011	_	2	

- Polyethylene liner and phenolic resin shelf are chemically resistant
- Add option code VK to order a Corrugated Polyolefin Vent kit for these Units





Width	Right Hinged	Left Hinged	# of Doors	
15"	1916022	1916222	1	
18"	1116022	1116222	1	
24"	1216022	1216222	1	
30"	1316022	-	2	
36"	1416022	-	2	
42"	1816022	-	2	
48"	1516022	-	2	
513mm	1616022	1616222	1	
1000mm	1716022	-	2	

- Polyethylene liner and phenolic resin shelf are chemically resistant
- · Add option code VK to order a Corrugated Polyolefin Vent kit for these Units

See page E18 for vent hole locations.





ACID STORAGE CABINETS

28-3/4" HIGH X 18"/22" DEEP

Acid Storage Units (18" Deep)





Width	Right Hinged	Left Hinged	# of Doors	
15"	1941011	1941211	1	
18"	1141011	1141211	1	
24"	1241011	1241211	1	
30"	1341011	-	2	
36"	1441011	-	2	
42"	1841011	-	2	
48"	1541011	-	2	
513mm	1641011	1641211	1	
1000mm	1741011	-	2	

- · Polyethylene liner and phenolic resin shelf are chemically resistant
- Add option code VK to order a Corrugated Polyolefin Vent kit for these Units





Width	Right Hinged	Left Hinged	# of Doors	
15"	1946022	1946222	1	
18"	1146022	1146222	1	
24"	1246022	1246222	1	
30"	1346022	-	2	
36"	1446022	-	2	
42"	1846022	-	2	
48"	1546022	-	2	
513mm	1646022	1646222	1	
1000mm	1746022	-	2	

- · Polyethylene liner and phenolic resin shelf are chemically resistant
- · Add option code VK to order a Corrugated Polyolefin Vent kit for these Units

Acid Storage Units (22" Deep)





Width	Right Hinged	Left Hinged	# of Doors	
15"	1921011	1921211	1	
18"	1121011	1121211	1	
24"	1221011	1221211	1	
30"	1321011	-	2	
36"	1421011	-	2	
42"	1821011	-	2	
48"	1521011	-	2	
513mm	1621011	1621211	1	
1000mm	1721011	-	2	

- Polyethylene liner and phenolic resin shelf are chemically resistant
- · Add option code VK to order a Corrugated Polyolefin Vent kit for these Units





Width	Right Hinged	Left Hinged	# of Doors	
15"	1926022	1926222	1	
18"	1126022	1126222	1	
24"	1226022	1226222	1	
30"	1326022	-	2	
36"	1426022	-	2	
42"	1826022	-	2	
48"	1526022	-	2	
513mm	1626022	1626222	1	
1000mm	1726022	-	2	

- · Polyethylene liner and phenolic resin shelf are chemically resistant
- Add option code VK to order a Corrugated Polyolefin Vent kit for these Units







ACID STORAGE CABINETS

31-3/4" HIGH X 22" DEEP

of Doors

1

1

2

2

2

2

1

2

See page E18 for vent hole locations.

Acid Storage Units



up to 24" wide





1000mm

• Polyethylene liner and phenolic resin shelf are chemically resistant

1781011

· Add option code VK to order a Corrugated Polyolefin Vent kit for these Units

1681211

Left Hinged

1981211

1181211

1281211





Width **Right Hinged** Left Hinged # of Doors 15" 1986022 1986222 1 18" 1186022 1186222 24" 1286022 1286222 1 30" 1386022 2 1486022 2 36" 42" 2 1886022 2 48" 1586022 513mm 1686022 1686222 1 1000mm 1786022 2

- Polyethylene liner and phenolic resin shelf are chemically resistant
- · Add option code VK to order a Corrugated Polyolefin Vent kit for these Units

ACID STORAGE CABINETS


84" HIGH X 18"/22" DEEP

Acid Storage Units



Width	18" Deep	22" Deep
30"	6322480	6332480
36"	6422480	6432480
42"	6622480	6632480
48"	6522480	6532480

- The plastic liner and all shelves, including four adjustable shelves and one fixed center shelf, are chemically resistant
- 3" high x 4" diameter stub for exhaust duct connection located 4-3/8" from rear of cabinet
- · Add option code 04 to add a standard louver configuration to this unit

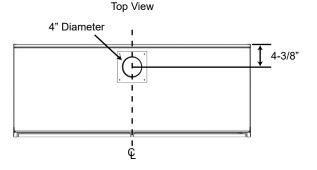

See page E19 for vent hole locations.





INSULATED SOLVENT CABINETS

34-3/4" HIGH X 18"/22" DEEP

Insulated Solvent Units (18" Deep)



Width	Right Hinged	Left Hinged	# of Doors
15"	1932011	1932111	1
18"	1132011	1132111	1
24"	1232011	1232111	1
30"	1332011	-	2
36"	1432011	-	2
42"	1832011	-	2
48"	1532011	-	2
513mm	1632011	1632111	1
1000mm	1732011	-	2
0 - 15 - 1 1 -			

- · Self closing door(s) provided
- · Add option code FL to order self closing door(s) with fusible link





AMMADDE E	FLANMADLE automates
" wide	30" and wider

Width	Right Hinged	Left Hinged	# of Doors	
15"	1932022	1932122	1	
18"	1132022	1132122	1	
24"	1232022	1232122	1	
30"	1332022	-	2	
36"	1432022	-	2	
42"	1832022	-	2	
48"	1532022	-	2	
513mm	1632022	1632122	1	
1000mm	1732022	-	2	
· Includes s	ittina heiaht UL lis	ted Insulated Solve	nt Cabinet	

- Includes sitting height UL listed Insulated Solvent Cabinet
- Self closing door(s) provided
- Add option code FL to order self closing door(s) with fusible link
- · Upper apron assembles are 6" high. Not UL approved with electrical installed
- · If Variac is required, upper apron increases in height, contact Mott for details
- · Dimensional modifications may invalidate the UL listing

Insulated Solvent Units (22" Deep)









Width	Right Hinged	Left Hinged	# of Doors	
15"	1912011	1912111	1	
18"	1112011	1112111	1	
24"	1212011	1212111	1	
30"	1312011	-	2	
36"	1412011	-	2	
42"	1812011	-	2	
48"	1512011	-	2	
513mm	1612011	1612111	1	
1000mm	1712011	-	2	

- Self closing door(s) provided
- · Add option code FL to order self closing door(s) with fusible link

Width	Right Hinged	Left Hinged	# of Doors	
15"	1912022	1912122	1	
18"	1112022	1112122	1	
24"	1212022	1212122	1	
30"	1312022	-	2	
36"	1412022	-	2	
42"	1812022	-	2	
48"	1512022	-	2	
513mm	1612022	1612122	1	
1000mm	1712022	-	2	

- · Includes sitting height UL listed Insulated Solvent Cabinet
- · Self closing door(s) provided
- · Add option code FL to order self closing door(s) with fusible link
- · Upper apron assembles are 6" high. Not UL approved with electrical installed
- · If Variac is required, upper apron increases in height, contact Mott for details
- · Dimensional modifications may invalidate the UL listing





INSULATED SOLVENT CABINETS

28-3/4" HIGH X 18"/22" DEEP

See page E19 for vent hole locations.

Insulated Solvent Units (18" Deep)



Width	Right Hinged	Left Hinged	# of Doors	
15"	1942011	1942111	1	
18"	1142011	1142111	1	
24"	1242011	1242111	1	
30"	1342011	-	2	
36"	1442011	-	2	
42"	1842011	-	2	
48"	1542011	-	2	
513mm	1642011	1642111	1	
1000mm	1742011	-	2	

- · Self closing door(s) provided
- · Add option code FL to order self closing door(s) with fusible link

Insulated Solvent Units (22" Deep)



Width	Right Hinged	Left Hinged	# of Doors	
15"	1922011	1922111	1	
18"	1122011	1122111	1	
24"	1222011	1222111	1	
30"	1322011	-	2	
36"	1422011	-	2	
42"	1822011	-	2	
48"	1522011	-	2	
513mm	1622011	1622111	1	
1000mm	1722011	-	2	

- · Self closing door(s) provided
- · Add option code FL to order self closing door(s) with fusible link

INSULATED SOLVENT CABINETS

31-3/4" HIGH X 22" DEEP

Insulated Solvent Units



Width	Right Hinged	Left Hinged	# of Doors	
15"	1982011	1982111	1	
18"	1182011	1182111	1	
24"	1282011	1282111	1	
30"	1382011	-	2	
36"	1482011	-	2	
42"	1882011	-	2	
48"	1582011	-	2	
513mm	1682011	1682111	1	
1000mm	1782011	-	2	

- · Self closing door(s) provided
- Add option code FL to order self closing door(s) with fusible link

INSULATED SOLVENT CABINETS

44" HIGH X 18" DEEP

Insulated Solvent Units



- Solvent Units are designed to OSHA and NFPA 30 standards
- Self closing doors provided
- Add option code FL to order self closing doors with fusible link
- These units include threaded vent connections in the back (standard location) with flame arresting screen and a three point locking handle
- For threaded vent connections on both sides with flame arresting screen use option code VS
- Includes two adjustable shelves

Width	Item Number
36"	6444160
42"	6644160
48"	6544160





QUALITY BY DESIGN

INSULATED SOLVENT CABINETS

65" HIGH X 18"/22" DEEP

See page E19 for vent hole locations.

Insulated Solvent Units (18" Deep)



Width	Item Number
36"	6454160
42"	6654160
48"	6554160

- · Solvent Units are designed to OSHA and NFPA 30 standards
- · Self closing doors provided
- · Add option code FL to order self closing doors with fusible link
- These units include threaded vent connections in the back (standard location) with flame arresting screen and a three point locking handle
- For threaded vent connections on both sides with flame arresting screen use option code VS
- · Includes three adjustable shelves

Insulated Solvent Units (22" Deep)



Width	Item Number
36"	6464160
42"	6664160
48"	6564160

- · Solvent Units are designed to OSHA and NFPA 30 standards
- · Self closing doors provided
- · Add option code FL to order self closing doors with fusible link
- These units include threaded vent connections in the back (standard location) with flame arresting screen and a three point locking handle
- · For threaded vent connections on both sides with flame arresting screen use option code VS
- · Includes three adjustable shelves








SOLVENT CABINETS

34-3/4" HIGH X 18"/22" DEEP

rage Units (18" Deen) See page E19 for vent hole locations.

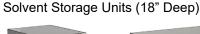





Width **Right Hinged** # of Doors Left Hinged 1932211 1932311 15" 1 18" 1132211 1132311 24" 1232211 1232311 1 30" 1332211 2 36" 2 1432211 42" 1832211 2 48" 1532211 2 1632311 513mm 1632211 1 1000mm 1732211 2

Solvent Storage Units (22" Deep)





Width Right Hinged Left Hinged # of Doors 15" 1912311 1912211 1 1112311 18" 1112211 24" 1212311 1212211 1 30" 1312211 2 36" 1412211 2 42" 1812211 2 48" 1512211 2 513mm 1612211 1612311 1 1000mm 1712211 2

- · Solvent Units are not insulated
- · These units are designed to meet OSHA and NFPA 30 standards
- Self closing door(s) provided
- Add option code FL to order self closing door(s) with fusible link

SOLVENT CABINETS

28-3/4" HIGH X 18"/22" DEEP

Solvent Storage Units (18" Deep)





width	Right Hinged	Left Hinged	# of Doors	
15"	1942211	1942311	1	
18"	1142211	1142311	1	
24"	1242211	1242311	1	
30"	1342211	-	2	
36"	1442211	-	2	
42"	1842211	-	2	
48"	1542211	-	2	
513mm	1642211	1642311	1	
1000mm	1742211	-	2	

Solvent Storage Units (22" Deep)





Width	Right Hinged	Left Hinged	# of Doors	
15"	1922211	1922311	1	
18"	1122211	1122311	1	
24"	1222211	1222311	1	
30"	1322211	-	2	
36"	1422211	-	2	
42"	1822211	-	2	
48"	1522211	-	2	
513mm	1622211	1622311	1	
1000mm	1722211	_	2	

- · Solvent Units are not insulated
- These units are designed to meet OSHA and NFPA 30 standards
- Self closing door(s) provided
- Add option code FL to order self closing door(s) with fusible link





QUALITY BY DESIGN

SOLVENT CABINETS

31-3/4" HIGH X 22" DEEP

See page E19 for vent hole locations.

Solvent Storage Units





Width	Right Hinged	Left Hinged	# of Doors
15"	1982211	1982311	1
18"	1182211	1182311	1
24"	1282211	1282311	1
30"	1382211	-	2
36"	1482211	-	2
42"	1882211	-	2
48"	1582211	-	2
513mm	1682211	1682311	1
1000mm	1782211	-	2

- · Solvent Units are not insulated
- · These units are designed to meet OSHA and NFPA 30 standards
- · Self closing door(s) provided
- · Add option code FL to order self closing door(s) with fusible link

SOLVENT CABINETS

44" HIGH X 18"DEEP

Solvent Storage Units



- · Solvent Units are not insulated
- · These units are designed to OSHA and NFPA 30 standards
- · Self closing doors provided
- · Add option code FL to order self closing doors with fusible link
- These units include threaded vent connections in back with flame arresting screen and a three point locking handle
- · Add option code VS for vent holes on side of cabinet
- Includes two adjustable shelves

Width	Item Number
36"	6443160
42"	6643160
48"	6543160

SOLVENT CABINETS

65" HIGH X 18"/22" DEEP

Solvent Storage Units (18" Deep)



Width	Item Number
36"	6453160
42"	6653160
48"	6553160

· Solvent Units are not insulated.

- These units are designed to OSHA and NFPA 30 standards
- Self closing doors provided
- · Add option code FL to order self closing doors with fusible link
- These units include threaded vent connections in back with flame arresting screen and a three point locking handle
- · Add option code VS for vent holes on side of cabinet
- Includes two adjustable shelves

Solvent Storage Units (22" Deep)



Width	Item Number
36"	6463160
42"	6663160
48"	6563160

- · Solvent Units are not insulated.
- These units are designed to OSHA and NFPA 30 standards
- Self closing doors provided
- · Add option code FL to order self closing doors with fusible link
- These units include threaded vent connections in back with flame arresting screen and a three point locking handle
- · Add option code VS for vent holes on side of cabinet
- · Includes two adjustable shelves







PUMP STORAGE CABINETS

34-3/4" HIGH X 18" DEEP

See page E19 for vent hole locations. **Pump Storage Units**















Width	Right Hinged	Left Hinged	# of Doors
15"	193P011	193P211	1
18"	113P011	113P211	1
24"	123P011	123P211	1
30"	133P011	-	2
36"	143P011	-	2
42"	183P011	-	2
48"	153P011	-	2
513mm	163P011	163P211	1
1000mm	173P011	-	2

- See "Note 1" on page E11
- · Add option code PF for optional pump storage cabinet fan

Width	Right Hinged	Left Hinged	# of Doors
15"	193P022	193P222	1
18"	113P022	113P222	1
24"	123P022	123P222	1
30"	133P022	-	2
36"	143P022	-	2
42"	183P022	-	2
48"	153P022	-	2
513mm	163P022	163P222	1
1000mm	173P022	-	2

- See "Note 1" on page E11
- Add option code PF for optional pump storage cabinet fan

PUMP STORAGE CABINETS

34-3/4" HIGH X 22" DEEP

Pump Storage Units







30" and wider





Width **Right Hinged** Left Hinged # of Doors 15" 191P011 191P211 1 18" 111P011 111P211 121P211 24" 121P011 1 30" 2 131P011 36" 141P011 2 2 42" 181P011 48" 151P011 2 513mm 161P011 161P211 1 1000mm 171P011 2

- See "Note 1" on page E11
- · Add option code PF for optional pump storage cabinet fan

Width	Right Hinged	Left Hinged	# of Doors
15"	191P022	191P222	1
18"	111P022	111P222	1
24"	121P022	121P222	1
30"	131P022	-	2
36"	141P022	-	2
42"	181P022	-	2
48"	151P022	-	2
513mm	161P022	161P222	1
1000mm	171P022	-	2

- See "Note 2" on page E11
- · Add option code PF for optional pump storage cabinet fan

Pump Storage





28-3/4" HIGH X 18" DEEP

PUMP STORAGE CABINETS

Pump Storage Units











up to 24" wide



Width	Right Hinged	Left Hinged	# of Doors
15"	194P011	194P211	1
18"	114P011	114P211	1
24"	124P011	124P211	1
30"	134P011	-	2
36"	144P011	-	2
42"	184P011	-	2
48"	154P011	-	2
513mm	164P011	164P211	1
1000mm	174P011	-	2

- See "Note 1" on page E11
- Add option code PF for optional pump storage cabinet fan

Width	Right Hinged	Left Hinged	# of Doors
15"	194P022	194P222	1
18"	114P022	114P222	1
24"	124P022	124P222	1
30"	134P022	-	2
36"	144P022	-	2
42"	184P022	-	2
48"	154P022	-	2
513mm	164P022	164P222	1
1000mm	174P022	-	2

- See "Note 2" on page E11
- Add option code PF for optional pump storage cabinet fan

PUMP STORAGE CABINETS

28-3/4" HIGH X 22" DEEP

Pump Storage Units







30" and wider





Width	Right Hinged	Left Hinged	# of Doors
15"	192P011	192P211	1
18"	112P011	112P211	1
24"	122P011	122P211	1
30"	132P011	-	2
36"	142P011	-	2
42"	182P011	-	2
48"	152P011	-	2
513mm	162P011	162P211	1
1000mm	172P011	-	2

- See "Note 1" on page E11
- · Add option code PF for optional pump storage cabinet fan

Width	Right Hinged	Left Hinged	# of Doors
15"	192P022	192P222	1
18"	112P022	112P222	1
24"	122P022	122P222	1
30"	132P022	-	2
36"	142P022	-	2
42"	182P022	-	2
48"	152P022	-	2
513mm	162P022	162P222	1
1000mm	172P022	-	2

- See "Note 2" on page E11
- · Add option code PF for optional pump storage cabinet fan







PUMP STORAGE CABINETS

31-3/4" HIGH X 22" DEEP

Pump Storage Units





30" and wider

Width Right Hinged Left Hinged # of Doors 15" 198P011 198P211 1 18" 118P011 118P211 24" 128P011 128P211 1 30" 138P011 2 36" 148P011 2 42" 2 188P011 48" 158P011 2 168P011 168P211 513mm 1 1000mm 178P011 2

Note 1:

- · Pump Storage Cabinets with full doors are insulated with sound absorption foam for noise reduction.
- Pump hose connection and 120V/20amp,(15amp in Canada) duplex receptacle is located on the rear of the cabinet and includes 8 feet of electrical whip. Duplex and cover plate are black.
- · These units include a pull-out tray and removable access panel at the back of the cabinet.
- Add option code PF for optional pump cabinet fan. Fan details: 115 VAC, 120 X 120 X 38mm, 81 CFM
- · Vibration isolators are available, contact Mott for details.





30" and wider

Width	Right Hinged	Left Hinged	# of Doors
15"	198P022	198P222	1
18"	118P022	118P222	1
24"	128P022	128P222	1
30"	138P022	-	2
36"	148P022	-	2
42"	188P022	-	2
48"	158P022	-	2
513mm	168P022	168P222	1
1000mm	178P022	-	2

Note 2:

- · Pump Storage Cabinets with a false panel are insulated with sound absorption foam for noise reduction.
- Pump hose connection is located on the rear of these cabinets.
- · Cabinets are supplied with switch and outlet pre-wired to junction box on back panel and includes 8 feet of electrical whip.
- Switch shown on the false panel controls interior 120V/20amp,(15amp in Canada) duplex receptacle for remote ON/OFF control.
- Duplex and cover plate are black.
- · These units include a pull-out tray and removable access panel at the back of the cabinet.
- Add option code PF for optional pump cabinet fan. Fan details: 115 VAC, 120 X 120 X 38mm, 81 CFM.
- · Vibration isolators are available, contact Mott for details.

PUMP STORAGE CABINET OPTIONS

Pump Storage Cabinet With No Floor



False Panel Unit With NFD Option Shown

Option Code

NFD

- Pump cabinet with no floor, toe kick integral to the door and dolly with 1" high lip with swivel casters with brake
- NFD option can be applied on all Pump Storage Cabinets with full height doors and false panel units
- Roll-out dolly to match material of cabinet

Option Code

NF

- Pump cabinet with no floor and toe kick integral to the door
- NF option can be applied on all Pump Storage Cabinets with full height doors and false panel units





QUALITY BY DESIGN

CONTROL PANEL CABINETS

34-3/4" HIGH X 22" DEEP

Control Panel Units





Right Hinged	Left Hinged	# of Doors
1015022	1015222	1
1915022	1915222	1
1115022	1115222	1
1215022	1215222	1
1315022	-	2
1415022	-	2
1815022	-	2
1515022	-	2
1615022	1615222	1
1715022	-	2
	1015022 1915022 1115022 1215022 1315022 1415022 1815022 1515022 1615022	1015022 1015222 1915022 1915222 1115022 1115222 1215022 1215222 1315022 - 1415022 - 1815022 - 1515022 - 1615022 1615222

[·] Includes an adjustable shelf

CONTROL PANEL CABINETS

28-3/4" HIGH X 22" DEEP

Control Panel Units





Width	Right Hinged	Left Hinged	# of Doors
12"	1025022	1025222	1
15"	1925022	1925222	1
18"	1125022	1125222	1
24"	1225022	1225222	1
30"	1325022	-	2
36"	1425022	-	2
42"	1825022	-	2
48"	1525022	-	2
513mm	1625022	1625222	1
1000mm	1725022	-	2

Does not include a shelf

CONTROL PANEL CABINETS

31-3/4" HIGH X 22" DEEP

Control Panel Units







30" and wider

Width	Right Hinged	Left Hinged	# of Doors
12"	1085022	1085222	1
15"	1985022	1985222	1
18"	1185022	1185222	1
24"	1285022	1285222	1
30"	1385022	-	2
36"	1485022	-	2
42"	1885022	-	2
48"	1585022	-	2
513mm	1685022	1685222	1
1000mm	1785022	-	2

[·] Includes an adjustable shelf





GLASSWARE DRYING CABINETS

34-3/4" HIGH X 18" DEEP

Glass Drying Units





Width	Item Number	
30"	1330822	
36"	1430822	
42"	1830822	
48"	1530822	
1000mm	1730822	

- Glassware Drying Units have double wall insulated construction
- These units include one perforated stainless steel shelf, stainless steel drip tray in cabinet bottom, a pilot light and 750 watt heating element with adjustable temperature up to 121° C / 250° F
- The controls are removable for easy element replacement from front of cabinet
- Operates on 110V, AC with a grounded cord. Includes 6' of electrical cord with NEMA 5-15P

GLASSWARE DRYING CABINETS

34-3/4" HIGH X 22" DEEP

Glass Drying Units





Width	Item Number
30"	1310822
36"	1410822
42"	1810822
48"	1510822
1000mm	1710822

- Glassware Drying Units have double wall insulated construction
- These units include one perforated stainless steel shelf, stainless steel drip tray in cabinet bottom, a pilot light and 750 watt heating element with adjustable temperature up to 121° C / 250° F
- The controls are removable for easy element replacement from front of cabinet
- Operates on 110V, AC with a grounded cord. Includes 6' of electrical cord with NEMA 5-15P

GLASSWARE DRYING CABINETS

31-3/4" HIGH X 22" DEEP

Glass Drying Units





Width	Item Number
30"	1380822
36"	1480822
42"	1880822
48"	1580822
1000mm	1780822

- Glassware Drying Units have double wall insulated construction
- These units include one perforated stainless steel shelf, stainless steel drip tray in cabinet bottom, a pilot light and 750 watt heating element with adjustable temperature up to 121° C / 250° F
- The controls are removable for easy element replacement from front of cabinet
- Operates on 110V, AC with a grounded cord. Includes 6' of electrical cord with NEMA 5-15P





SECURITY CABINETS

24"/30" HIGH X 8"/14" DEEP

Double Door Security Cabinets (8" Deep)



Width	Right Hinged	Left Hinged	Height
15"	5700260	5701260	24"
18"	5100260	5101260	24"
24"	5200260	5201260	24"
15"	5790260	5791260	30"
18"	5190260	5191260	30"
24"	5290260	5291260	30"

- · Double door cabinet complete with locks on both doors
- · Inside lock is keyed separately from outside lock
- · Shelves adjustable in 1" increments
- · Add option code RL to order the optional remote LED light that indicates the door is unlocked

Double Door Security Cabinets (14" Deep)



Width	Right Hinged	Left Hinged	Height
15"	5770260	5771260	24"
18"	5170260	5171260	24"
24"	5270260	5271260	24"
15"	5710260	5711260	30"
18"	5110260	5111260	30"
24"	5210260	5211260	30"

- · Double door cabinet complete with locks on both doors
- · Inside lock is keyed separately from outside lock
- · Shelves adjustable in 1" increments
- · Add option code RL to order the optional remote LED light that indicates the door is unlocked

GOGGLE CABINETS

30" HIGH X 10-1/2" DEEP

Goggle Cabinet





Width **Item Number**

- 24" 52C0260
- · Goggle cabinets expose safety goggles to germicidal ultraviolet
- · When the doors are closed and locked, the inside of the cabinet is flooded with germicidal ultra violet light
- · A safety interlock device automatically shuts off the ultraviolet light when the doors are opened
- Dimensions: 24" (60.96cm) Wide, 30" (76.2cm) High, 10-1/2" (26.67cm) Deep
- Cabinet can be mounted on a wall or table and includes two safety door switches and key locks
- · Holds 30 pairs of goggles (not included)
- · Operates on 110V, AC with a grounded cord
- Includes 6' of electrical cord with NEMA 5-15P
- · Cabinet is ONLY available in Bright White #693226







RECESSED OPERATING ROOM CABINETS

60" HIGH X 18" DEEP

Construction Features:

- Powder coated steel construction. Also available in type 304 stainless steel, or type 316 stainless steel
- All exposed welds are ground and polished to match a #4 finish
- Cabinets will fit into an opening of: cabinet width minus 1" to 2" x 57-3/4" high x 18"minimum depth

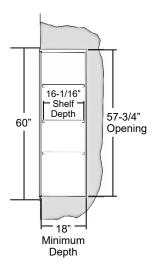

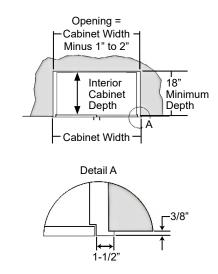

Hinged Framed Glass Door Unit



Width	Right Hinged	Left Hinged	# of Doors
18"	4170030	4170230	1
24"	4270030	4270230	1
30"	4370030	-	2
36"	4470030	-	2
42"	4670030	-	2
48"	4570030	-	2

- These units include three shelves, adjustable in ½"increments
- Glass doors are supplied with 6mm laminated safety glass

Drawer Desk Unit



Width	Item Number
18"	4176000
24"	4276000
30"	4376000
36"	4476000
42"	4676000
48"	4576000

These units include one metal pullout writing board

Hinged Door Unit



Width	Right Hinged	Left Hinged	# of Doors
18"	4170060	4170260	1
24"	4270060	4270260	1
30"	4370060	-	2
36"	4470060	-	2
42"	4670060	-	2
48"	4570060	-	2

• These units include three shelves, adjustable in ½" increments

Recessed Cabinet Joint



Item Number 4070000

- 57-3/4"high
- Used to join two recessed cabinets side by side





QUALITY BY DESIGN

of Doors

1

1

2

2 2

2

1

SIGMA SEAL™ CABINETS

34-3/4" HIGH X 28" DEEP

Left Hinged

11B0211

12B0211

16B0211

Sigma Seal[™] Units





Width Right Hinged 18" 11B0011 24" 12B0011 30" 13B0011 36" 14B0011 42" 18B0011 48" 15B0011 513mm 16B0011

- · All stainless steel construction
- · Doors include gasket
- Front of the units are water tight for wash down applications

SIGMA SEAL™ CABINETS

28-3/4" HIGH X 28" DEEP

Sigma Seal[™] Units





Width	Right Hinged	Left Hinged	# of Doors
18"	11C0011	11C0211	1
24"	12C0011	12C0211	1
30"	13C0011	-	2
36"	14C0011	-	2
42"	18C0011	-	2
48"	15C0011	-	2
513mm	16C0011	16C0211	1

- · All stainless steel construction
- Doors include gasket
- · Front of the units are water tight for wash down applications





MOBILE SINK CABINETS

MOBILE SINK CABINET

Mobile Sink Cabinets






Width	Right Hinged	Left Hinged	# of Doors
24"	MS24011	MS24111	1
30"	MS30011	-	2
36"	MS36011	-	2
42"	MS42011	-	2
48"	MS48011	-	2


This cabinet provides a flexible supply of water in any research, school or care facility and can be used all through the facility when access to water is not available.

- Available in either type 304 stainless steel or powder coated steel with protective Sigma Barrier™ coating.
- Durable stainless steel work top with raised perimeter marine edge with integrally welded sink allows for easy cleaning.
- Work top provides a 1" overhang on all four sides.
- Automatic faucet with infrared sensor provides hands-free
- Plate mounted, non-marking 5" total lock swivel castors with santoprene rubber wheels allows the user to quickly and easily move to another work station, as tasks change.
- · 5 US gallon fresh water container, 5 US gallon waste water container with drain spigot for easy release of waste water into the floor drain.
- Automatic pump provides up to 1.1 US gallons per minute with maximum current draw of 1.9 amperes at 12 volt.
- Onboard sealed 12volt battery and 1.5 amp automatic charger including 6' cord and plug. Manual power switch on cabinet side.
- Doors on front and back of cabinet provides convenient access for servicing.
- Shipped assembled and skidded.
- Use option code WS for supply hose to hook up with facility water supply, if fresh water tank and pump is not required.


Top View of Work Top

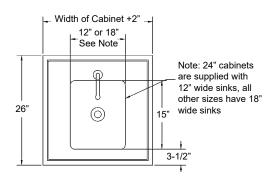

Side View

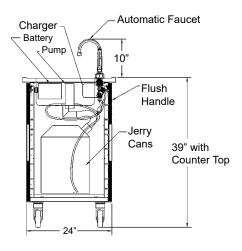

Front View

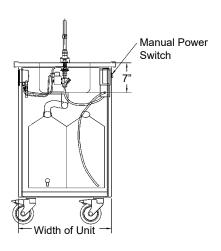

V931






ACID CABINETS


CONSTRUCTION DETAILS

Acid Cabinet Details

- · Polyethylene liner and phenolic resin shelf are chemically resistant in acid cabinets.
- · The removable liner is secured with capped stainless steel screws and includes a vent connection.
- Vent location applies to 28-3/4", 31-3/4", and 34-3/4" height full door (Fig. 1) cabinets. For false panel units (Fig. 2) dimension A varies with the height of the cabinet. Cabinets are marked "ACID" in red letters.
- Note: acid cabinets must be vented. Mott recommends venting acid cabinets with a Polyolefin vent kit (PVSS or PVSA for field retrofit
 see fume hood section).

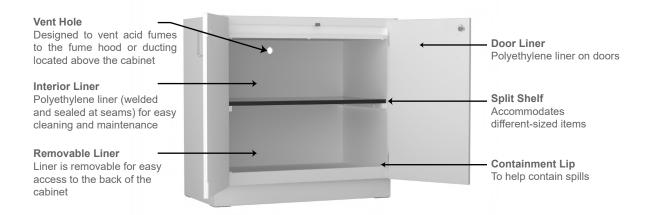

Cabinet Vent Hole Locations Acid Cabinet Rear View

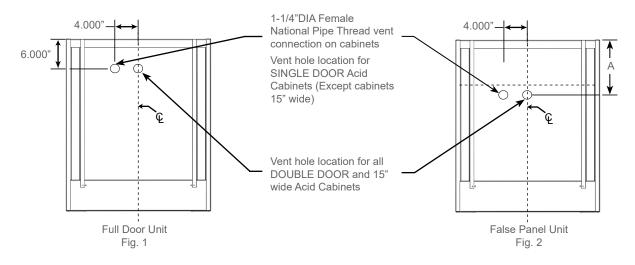






SOLVENT/PUMP CABINETS

CONSTRUCTION DETAILS

Insulated & Non-Insulated Solvent Cabinet Details

- Insulated solvent storage cabinets are fabricated with 2" insulated double walls, meet National Fire Protection Association standards and undercounter sizes are approved to UL1275.
- Uninsulated safety cabinets meet National Fire Protection Association standards.
- Insulated and non-insulated solvent cabinet floors have one 2" deep liquid tight spill containment pan.
- Doors are mounted on continuous hinges with a three point locking self latching handle.
- Each cabinet has two 2" diameter vents in the rear of the cabinet, complete with flame arrestor and vent covers. Vent location applies to 28-3/4", 31-3/4", 34-3/4", 44" and 65" height cabinets.
- Cabinets are provided with one full depth shelf, adjustable in 1/2" increments.
- Cabinets are marked "FLAMMABLE KEEP FIRE AWAY" in red letters.
- Self closing doors are standard, with option to add a fusible link.
- Note: NFPA30 Fire code does not require solvent cabinets to be vented. Do not vent solvent cabinets to fume hood. If cabinet must be vented use fireproof pipe. Contact Mott for details.


Cabinet Vent Hole Locations Solvent Cabinet Rear View

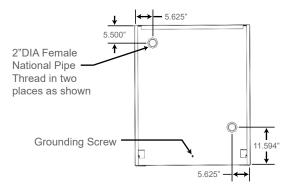

Pump Storage Cabinet Details

- Pump cabinets are insulated with 1" thick acoustic foam which provides sound absorption and noise reduction over a wide range of frequencies.
- Foam is incorporated to the interior of the cabinet including the door and it is contained with a mechanical inner steel frame. The foam is fire resistant and complies with MVSS 302 AND UL94 HF-1 requirements.
- Two vent holes are incorporated in the cabinet back for connections to centralized lab exhaust systems. Vent location applies to 28-3/4", 31-3/4", and 34-3/4" height full door (Fig. 1) cabinets. For false panel units (Fig. 2) dimension A varies with the height of the cabinet. Custom vent cutouts are available.
- Back access panel is easily removable to provide convenient access to venting or electrical connections after installation.
- One 120 Volt / 20 Amp (15 Amp in Canada) electrical duplex receptacle provides power to the vacuum pump (final power connection required).
- The switch that controls receptacles is located on the front panel on false panel units; full height door units are typically wired through the fume hood.

Cabinet Vent Hole Locations Pump Storage Cabinet Rear View

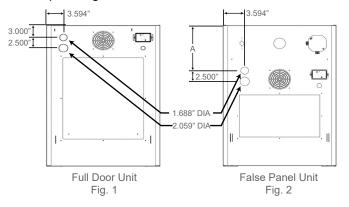





SPECIALTY CABINET OPTIONS

Option	Description
00	Modified door/drawer configuration
02	Label holders on all doors/drawers
03	Numeric or alpha identification disks on all doors/drawers (customer details required)
04	Standard louver configuration on doors (6 louvers at 3" wide x 5-1/2" long)
05	Non-standard louver configuration on doors (customer drawing required)
07	Cut-outs for electrical and/or plumbing services (customer drawing required)
80	Security proof panel
09	Stainless steel waste flap
10	One fixed drawer divider per drawer (centered in drawer)
11	Multiple drawer dividers per drawer (customer drawing required)
14	Ventilation grill in toe kick base
15	Finished back on cabinet
16	Ship handles not installed
17	Nylon levelers
18	Horizontally mounted handles
26	200lb full extension drawer runners
40	One extra cabinet shelf
41	Two extra cabinet shelves
42	One interior cabinet shelf on full extension runners (replaces fixed shelf)
43	Two interior cabinet shelf on full extension runners (replaces fixed shelf)
44	Three interior cabinet shelf on full extension runners (replaces fixed shelf)
50	304 stainless steel drawer bodies in power coated cold rolled steel cabinet
51	Entire 304 stainless steel component
52	Entire 316 stainless steel component
53	Entire galvenneal steel cabinet construction (Sigma Barrier™)
56	Other customer specified material
59	Customer specified runners (customer details required)
70	Front toe kick
71	Toe kick - front & left side
72	Toe kick - front & right side
73	Toe kick - front & 2 sides
74	Toe kick - front & back

Notes: *Designed for use on solvent cabinets with self-closing. A fusible link holds doors open during use and melts at 165° F to automatically close the doors in the event of a fire.

Option	Description
75	Toe kick - front, back & two sides
76	Replace toe kick with four legs
77	Replace toe kick on base cabinets with four casters
78	Replace toe kick with customer specified casters
79	Other customer specified base or toe kick
BL	BASE label (use this option code to replace standard label on acid cabinets)
CL	CORROSIVES label (use this option code to replace standard label on acid cabinets)
CSA	ESA inspection required with a CSA label
ESA	ESA inspection required
C25	Replace toe kick on base cabinets with 250lb casters
C35	Replace toe kick on base cabinets with 350lb casters
FL	Fusible link for self closing doors for insulated solvent unit*
F1	Welded finished top
F2	Welded finished top & back
LA	Bilingual label (French/English)
NB	No base, cabinet to sit on floor, levelers included
NC	No base, cabinet to sit on 4" high curb, levelers not included
NF	No cabinet floor and toe kick is integral with door
NFD	Pump cabinet unit with no cabinet floor and toe kick is integral with the door, dolly included
OL	Overlay design (not available on Solvent cabinets)
PF	Fan for pump cabinet
PP	Pre-plumbed (option code not to be used for fume hoods and customer details required)
PT	Single Polypropylene tray for acid cabinets
PW	Pre-wired (option code not to be used for fume hoods, customer details required)
P1	Plumbing - mount only
P2	Factory pre-plumbed up
RL	Remote LED light that indicates the door is unlocked, for Double Door Security Cabinets
SK	Skid item
T1	Replace one shelf with one pull out tray on full extension runners with 2" high lip
T2	Replace two shelf with two pull out tray on full extension runners with 2" high lip
Т3	Replace three shelf with three pull out tray on full extension runners with 2" high lip
VK	Corrugated Polyolefin vent kit for acid cabinets
VS	Vent holes on the side of the cabinet
WF	Wood fronts - see specifications
WO	Wood fronts provided and installed by others





SUSPENDED CABINETS

24" HIGH X 22" DEEP

Open Units



Width	Item Number
15"	1950000
18"	1150000
24"	1250000
30"	1350000
36"	1450000
42"	1850000
48"	1550000
513mm	1650000
1000mm	1750000

Note: Add option code 15 to order a finished back on these units.

Full Door Units





up to 24" wide

30" and wider

Door/Drawer Units







de

Width	Right Hinged	Left Hinged	# of Doors
15"	1950011	1950211	1
18"	1150011	1150211	1
24"	1250011	1250211	1
30"	1350011	-	2
36"	1450011	-	2
42"	1850011	-	2
48"	1550011	-	2
513mm	1650011	1650211	1
1000mm	1750011	-	2

[·] Cabinets receive one shelf

Width	Right Hinged	Left Hinged	# of Doors
15"	1950022	1950222	1
18"	1150022	1150222	1
24"	1250022	1250222	1
30"	1350022	-	2
36"	1450022	-	2
42"	1850022	-	2
48"	1550022	-	2
513mm	1650022	1650222	1
1000mm	1750022	-	2

• Cabinets up to 24" wide receive one shelf, while cabinets 25" and wider do not have a shelf

Width	Right Hinged	Left Hinged
30"	1350041	1350014
36"	1450041	1450014
42"	1850041	1850014
48"	1550041	1550014
1000mm	1750041	1750014

- · Left hand Door/Drawer unit shown, right hand reversed
- · Cabinet receives one shelf

Sliding Door Units



Width	Item Number	
30"	1350711	
36"	1450711	
42"	1850711	
48"	1550711	
1000mm	1750711	
 Cabinets 	receive one shelf	





Width

15"

18"

24"

30"

36"

42"

48"

513mm

Width

15"

18" 24"

30"

36"

42"

48"

513mm

Width

15"

18"

24"

30"

36"

42"

48"

513mm 1000mm

1000mm

1000mm

Item Number

1950055

1150055

1250055

1350055

1450055

1850055

1550055

1650055

1750055

1950066 1150066

1250066

1350066

1450066

1850066

1550066

1650066

1750066

1950044

1150044

1250044

1350044

1450044

1850044

1550044 1650044

1750044

Item Number

Item Number

SUSPENDED CABINETS

24" HIGH X 22" DEEP

Note: Add option code 15 to order a finished back on these units.

Drawer Units







30" and wider







30" and wider







30" and wider



Width	Item Number
30"	1350097
36"	1450097
42"	1850097
48"	1550097
1000mm	1750097



Width	Item Number
30"	1350098
36"	1450098
42"	1850098
48"	1550098
1000mm	1750098

File Drawer Units







30" and wider

Width	Item Number
15"	1950122
18"	1150122
24"	1250122
30"	1350122
36"	1450122
42"	1850122
48"	1550122
513mm	1650122
1000mm	1750122





SUSPENDED CABINETS

File Drawer Units





up to 24" wide

30" and wider





up to 24" wide

30" and wider

Sink Units





Acid Storage Units





24" HIGH X 22" DEEP

Note: Add option code 15 to order a finished back on these units.

Width	Item Number
15"	1950177
18"	1150177
24"	1250177
30"	1350177
36"	1450177
42"	1850177
48"	1550177
513mm	1650177
1000mm	1750177

Width	Item Number
15"	1950188
18"	1150188
24"	1250188
30"	1350188
36"	1450188
42"	1850188
48"	1550188
513mm	1650188
1000mm	1750188

Width	Right Hinged	Left Hinged	# of Doors
15"	1953011	1953211	1
18"	1153011	1153211	1
24"	1253011	1253211	1
30"	1353011	-	2
36"	1453011	-	2
42"	1853011	-	2
48"	1553011	-	2
513mm	1653011	1653211	1
1000mm	1753011	-	2

Width	Right Hinged	Left Hinged	# of Doors
15"	1951011	1951211	1
18"	1151011	1151211	1
24"	1251011	1251211	1
30"	1351011	-	2
36"	1451011	-	2
42"	1851011	-	2
48"	1551011	-	2
513mm	1651011	1651211	1
1000mm	1751011	-	2

- · Polyethylene liner and phenolic resin shelf are chemically resistant in acid cabinets
- Add option code VK to order a Corrugated Polyolefin Vent Kit for these units
- · Note: Acid Cabinets must be vented. Mott recommends venting acid cabinets with a Polyolefin Vent Kit
- · Contact Mott for more details







SUSPENDED CABINETS

24" HIGH X 22" DEEP

2

Solvent Storage Units





30" and wider

Width **Right Hinged** Left Hinged # of Doors 15" 1952211 1952311 1 18" 1152211 1152311 24" 1252211 1252311 1 30" 1352211 2 36" 1452211 2 42" 2 1852211 48" 1552211 2 1652211 1652311 513mm 1

Note: Add option code 15 to order a finished back on these units.

1752211 · Solvent units are not insulated

1000mm

- Self closing doors are standard, use option FL to add a fusible link
- Solvent Units are designed to meet OSHA and NFPA 30 standards
- Contact Mott for more details

SUSPENDED CABINETS

18-3/4" HIGH X 22" DEEP

Open Units



Width	Item Number
15"	1960000
18"	1160000
24"	1260000
30"	1360000
36"	1460000
42"	1860000
48"	1560000
513mm	1660000
1000mm	1760000

Full Door Units





Width	Right Hinged	Left Hinged	# of Doors
15"	1960011	1960211	1
18"	1160011	1160211	1
24"	1260011	1260211	1
30"	1360011	-	2
36"	1460011	-	2
42"	1860011	-	2
48"	1560011	-	2
513mm	1660011	1660211	1
1000mm	1760011	-	2

Door/Drawer Units





Width	Right Hinged	Left Hinged	# of Doors
15"	1960022	1960222	1
18"	1160022	1160222	1
24"	1260022	1260222	1
30"	1360022	-	2
36"	1460022	-	2
42"	1860022	-	2
48"	1560022	-	2
513mm	1660022	1660222	1
1000mm	1760022	-	2





Width

15"

18"

24"

30"

36"

42"

48"

513mm

Width

1000mm

SUSPENDED CABINETS

18-3/4" HIGH X 22" DEEP

Note: Add option code 15 to order a finished back on these units.

Door/Drawer Units



Width	Right Hinged	Left Hinged
30"	1360041	1360014
36"	1460041	1460014
42"	1860041	1860014
48"	1560041	1560014
1000mm	1760041	1760014

[·] Right hand Door/Drawer Unit shown, left hand reversed

Item Number

1960044

1160044

1260044

1360044

1460044

1860044

1560044

1660044

1760044

Right Hinged

Drawer Units





up to 24" wide

30" and wider

Drawer Units





up to 24" wide

30" and wider

Width **Item Number** 15" 1960177 18" 1160177 24" 1260177 30" 1360177 36" 1460177 42" 1860177 1560177 48" 513mm 1660177 1000mm 1760177

Sink Units





15" 1963011 1963211 1 18" 1163011 1163211 1 24" 1263011 1263211 1 30" 1363011 2 36" 1463011 2 -42" 1863011 2 48" 1563011 2 1663011 1663211 513mm 1 1000mm 1763011 2

Left Hinged

of Doors





SUSPENDED COMPONENTS

Pencil Drawers and Pull Out Writing Tablets



Hanging CPU Holder



Width	Pencil Drawer	Pull Out Writing Tablet
15"	PEN0015	PUL0015
18"	PEN0018	PUL0018
20"	PEN0020	PUL0020
24"	PEN0024	PUL0024
30"	PEN0030	PUL0030

- Pencil Drawers and Pull Out Writing Tablets suspends from tables with hanging rails and are supplied with full extension runners
- When ordering a table with aprons, add option code HU to suspend Pencil Drawer
- Pencil Drawers and Pull Out Writing Tablets do not fit on a 22" deep table

Item Number

CPU1010

- Suspends from tables with hanging rails and is supplied with full extension runners to slide out CPU unit
- Inside clear dimension of CPU holder is 9-3/16"
- · Hanging CPU Holder does not fit on a 22" deep table
- When ordering a table with aprons, add option code HU to suspend CPU holder





Option

84

Description

Stainless steel bin insert

SUSPENDED CABINET OPTIONS

Option	Description
00	Modified door/drawer configuration
02	Label holders on all doors/drawers
03	Numeric or alpha identification disks on all doors/drawers (customer details required)
04	Standard louver configuration on doors (6 louvers at 3" wide x 5-1/2" high)
05	Non-standard louver configuration on doors (customer drawing required)
07	Cut-outs for electrical and/or plumbing services (customer drawing required)
08	Security proof panel(s)
09	Stainless steel waste flap
10	One fixed drawer divider per drawer (centered in drawer)
11	Multiple drawer dividers per drawer (customer drawing required)
12	Pull out writing board * - cold rolled steel
13	Pull out writing board - cold rolled steel & high pressure laminate writing surface
15	Finished back on cabinet
16	Ship handles not installed
18	Horizontally mounted handles on doors
26	200lb full extension drawer runners
40	One extra cabinet shelf
41	Two extra cabinet shelves
42	One interior cabinet shelf on full extension runners (replaces fixed shelf)
48	Keyboard tray (replaces drawer)
50	304 stainless steel drawer bodies in power coated cold rolled steel cabinet
51	Component made from 304 stainless steel
52	Component made from 316 stainless steel
53	Galvenneal steel cabinet construction (Sigma Barrier™)
54	304 stainless steel drawer bodies in galvenneal steel cabinet (Sigma Barrier™)
56	Other customer specified material
59	Customer specified runners (customer details required)

85	PVC bin insert
89	Lazy Susan storage unit for corner cabinets
90	Sink supports - for under counter mounted sinks
BL	BASE label (use this option code to replace standard label on acid cabinet)
CL	CORROSIVES label (use this option code to replace standard label on acid cabinet)
DC2	Magnetic door catch
F1	Welded finished top
F2	Welded finished top & back
F26	200lb full extension runners on file drawers
FCW	Front corners of doors and drawers welded
FL	Fusible link with self closer for insulated solvent unit**
101	100lb full extension runners with interlock
126	200lb full extension runners with interlock
LA	Bilingual label (French/English)
LF	Longitudinal file hangers on file drawers
OL	Overlay design (not available on Solvent Cabinets)
	Pre-plumbed (option code not to be used for fume
PP	hoods and customer details required)
PP PW	
	hoods and customer details required)
PW	hoods and customer details required) Pre-wired (customer details required)
PW PDC	hoods and customer details required) Pre-wired (customer details required) Positive door catch (standard on overlay cabinets)
PW PDC PFE	hoods and customer details required) Pre-wired (customer details required) Positive door catch (standard on overlay cabinets) Pull out writing board * on full extension runners
PW PDC PFE RB	hoods and customer details required) Pre-wired (customer details required) Positive door catch (standard on overlay cabinets) Pull out writing board * on full extension runners Removable back on drawer cabinet
PW PDC PFE RB SK	hoods and customer details required) Pre-wired (customer details required) Positive door catch (standard on overlay cabinets) Pull out writing board * on full extension runners Removable back on drawer cabinet Skid item Silk screen printing (to replace standard labels on acid
PW PDC PFE RB SK SP	hoods and customer details required) Pre-wired (customer details required) Positive door catch (standard on overlay cabinets) Pull out writing board * on full extension runners Removable back on drawer cabinet Skid item Silk screen printing (to replace standard labels on acid cabinet and/or flammable cabinet
PW PDC PFE RB SK SP S01	hoods and customer details required) Pre-wired (customer details required) Positive door catch (standard on overlay cabinets) Pull out writing board * on full extension runners Removable back on drawer cabinet Skid item Silk screen printing (to replace standard labels on acid cabinet and/or flammable cabinet 100lb full extension soft close drawer runners
PW PDC PFE RB SK SP S01 S26	hoods and customer details required) Pre-wired (customer details required) Positive door catch (standard on overlay cabinets) Pull out writing board * on full extension runners Removable back on drawer cabinet Skid item Silk screen printing (to replace standard labels on acid cabinet and/or flammable cabinet 100lb full extension soft close drawer runners 200lb full extension soft close drawer runners Replace one shelf with one pull out tray on full exten-
PW PDC PFE RB SK SP S01 S26 T1	hoods and customer details required) Pre-wired (customer details required) Positive door catch (standard on overlay cabinets) Pull out writing board * on full extension runners Removable back on drawer cabinet Skid item Silk screen printing (to replace standard labels on acid cabinet and/or flammable cabinet 100lb full extension soft close drawer runners 200lb full extension soft close drawer runners Replace one shelf with one pull out tray on full extension runners with 2" high lip Replace two shelf with two pull out tray on full extension
PW PDC PFE RB SK SP S01 S26 T1 T2	hoods and customer details required) Pre-wired (customer details required) Positive door catch (standard on overlay cabinets) Pull out writing board * on full extension runners Removable back on drawer cabinet Skid item Silk screen printing (to replace standard labels on acid cabinet and/or flammable cabinet 100lb full extension soft close drawer runners 200lb full extension soft close drawer runners Replace one shelf with one pull out tray on full extension runners with 2" high lip Replace two shelf with two pull out tray on full extension runners with 2" high lip Replace three shelf with three pull out tray on full
PW PDC PFE RB SK SP S01 S26 T1 T2 T3	hoods and customer details required) Pre-wired (customer details required) Positive door catch (standard on overlay cabinets) Pull out writing board * on full extension runners Removable back on drawer cabinet Skid item Silk screen printing (to replace standard labels on acid cabinet and/or flammable cabinet 100lb full extension soft close drawer runners 200lb full extension soft close drawer runners Replace one shelf with one pull out tray on full extension runners with 2" high lip Replace two shelf with two pull out tray on full extension runners with 2" high lip Replace three shelf with three pull out tray on full extension runners with 2" high lip

Wood fronts - see specifications

Wood fronts provided and installed by others

See Steel Specification & Options section for information on the following options:

- · Door and drawer handles
- Door hinges
- · Drawer body and suspension
- · Lock and key options

Note: *30" wide cabinets and larger receive two pull out writing boards, one pull out board is provided on full width drawer cabinets. Note: **Designed for use on solvent cabinets with self-closing. A fusible link holds doors open during use and melts at 165° F to automatically close the doors in the event of a fire.

WF

WO





MISCELLANEOUS ITEMS

Front Filler Scribes



	Standir	ng Ht.	Sittin	g Ht.	ADA	Ht.	84"	High
Width	Right Hand	Left Hand						
1"	FFR1001	FFL1001	FFR2001	FFL2001	FFR8001	FFL8001	FFR3001	FFL3001
2"	FFR1002	FFL1002	FFR2002	FFL2002	FFR8002	FFL8002	FFR3002	FFL3002
4"	FFR1004	FFL1004	FFR2004	FFL2004	FFR8004	FFL8004	FFR3004	FFL3004
6"	FFR1006	FFL1006	FFR2006	FFL2006	FFR8006	FFL8006	FFR3006	FFL3006
8"	FFR1008	FFL1008	FFR2008	FFL2008	FFR8008	FFL8008	FFR3008	FFL3008
10"	FFR1010	FFL1010	FFR2010	FFL2010	FFR8010	FFL8010	FFR3010	FFL3010
12"	FFR1012	FFL1012	FFR2012	FFL2012	FFR8012	FFL8012	FFR3012	FFL3012

- The Front Filler Scribe is used to close a space between a base cabinet and a wall.
- · A double angle is provided to enclose the raw edge of the panel.
- Filler panels can be ordered in 1" increments, up to 24" wide. Change the last two digits to the dimension required.
- · Right hand shown.
- Use option code "OL" when ordering Architectural Overlay design.

Front Filler Panels



Width	Standing Ht.	Sitting Ht.	ADA Ht.	84" High
1"	FFC1001	FFC2001	FFC8001	FFC3001
2"	FFC1002	FFC2002	FFC8002	FFC3002
4"	FFC1004	FFC2004	FFC8004	FFC3004
6"	FFC1006	FFC2006	FFC8006	FFC3006
8"	FFC1008	FFC2008	FFC8008	FFC3008
10"	FFC1010	FFC2010	FFC8010	FFC3010
12"	FFC1012	FFC2012	FFC8012	FFC3012
24"	FFC1024	FFC2024	FFC8024	FFC3024

- The Front Filler Panel is used to close a space between two base cabinets.
- Filler panels can be ordered in 1" increments, up to 24" wide. Change the last two digits to the dimension required.
- Use option code "OL" when ordering Architectural Overlay design.

Corner Filler Panels

Corner Filler Panel Scribes



Width & Depth	Standing Ht.	Sitting Ht.	ADA Ht.
1" x 1"	CFC1001	CFC2001	CFC8001
2" x 2"	CFC1002	CFC2002	CFC8002
4" x 4"	CFC1004	CFC2004	CFC8004
6" x 6"	CFC1006	CFC2006	CFC8006

- The Corner Filler Panel is used to close a corner space between two base
- Filler panels can be ordered in 1" increments, up to 24" wide. Change the last two digits to the dimension required.
- Use option code "OL" when ordering Architectural Overlay design.

Width & Depth	Standing Ht.	Sitting Ht.	ADA Ht.	84" High
1" x 1"	CF\$1001	CF\$2001	CFS8001	CF\$3001
2" x 2"	CF\$1002	CF\$2002	CF\$8002	CF\$3002
4" x 4"	CF\$1004	CFS2004	CF\$8004	CF\$3004
6" x 6"	CF\$1006	CF\$2006	CF\$8006	CF\$3006


- The Corner Filler Panel Scribes are used to close a corner space between two base cabinets.
- Filler Panel Scribes can be ordered in 1" increments, up to 24" wide. Change the last two digits to the dimension required.
- Two double angles are provided to enclose the raw edges of the panel.
- Use option code "OL" when ordering Architectural Overlay design.





MISCELLANEOUS ITEMS


Adjustable Corner Fillers

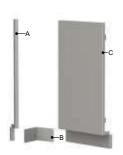

Item	Standing Ht.	Sitting Ht.	ADA Ht.	
A	FFC1001	FFC2001	FFC8001	
В	9900001	9900001	9900001	
С	FFC1018	FFC2018	FFC8018	

- The items A, B, & C combine to make the Adjustable Corner Filler Panel.
- This panel is used to close a corner space between two base cabinets. It offers the added advantage of being adjustable to enclose spaces up to 14" wide.
- Use option code "OL" when ordering Architectural Overlay design.

End Cover Panel Scribes

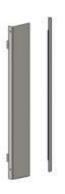

Width	Standing Ht.	Sitting Ht.	ADA Ht.	84"High
1"	EC\$1001	ECS2001	EC\$5001	EC\$3001
4"	EC\$1004	ECS2004	EC\$5004	EC\$3004
7"	EC\$1007	ECS2007	EC\$5007	EC\$3007
8"	EC\$1008	EC\$2008	EC\$5008	EC\$3008
10"	EC\$1010	EC\$2010	EC\$5010	EC\$3010
11"	EC\$1011	EC\$2011	EC\$5011	EC\$3011
120mm	ECS100A	ECS200A	ECS500A	ECS300A
266mm	ECS100B	ECS200B	ECS500B	ECS300B

- The End Cover Panel Scribes are used to close spaces between the back of a base cabinet and a wall. The Scribe feature allows for easy width adjustment
- Filler panels can be ordered in 1" increments, up to 24" wide. Change the last two digits to the dimension required.

End Cover Channels



Width	Standing Ht.	Sitting Ht.	ADA Ht.	84"High
1"	ECC1001	ECC2001	ECC5001	ECC3001
2"	ECC1002	ECC2002	ECC5002	ECC3002
4"	ECC1004	ECC2004	ECC5004	ECC3004
7"	ECC1007	ECC2007	ECC5007	ECC3007
10"	ECC1010	ECC2010	ECC5010	ECC3010
13"	ECC1013	ECC2013	ECC5013	ECC3013
14"	ECC1014	ECC2014	ECC5014	ECC3014
120mm	ECC100A	ECC200A	ECC500A	ECC300A

- The End Cover Channel performs the same task as the End Cover Panel Scribe but is of a fixed width.
- Filler panels can be ordered in 1" increments, up to 24" wide. Change the last two digits to the dimension required.





MISCELLANEOUS ITEMS

Finished Back Panels



Width	Standing Ht.	Sitting Ht.	ADA Ht.
15"	FBP1015	FBP2015	FBP6015
18"	FBP1018	FBP2018	FBP6018
24"	FBP1024	FBP2024	FBP6024
30"	FBP1030	FBP2030	FBP6030
36"	FBP1036	FBP2036	FBP6036
42"	FBP1042	FBP2042	FBP6042
48"	FBP1048	FBP2048	FBP6048

- Finished Back Panels are used to close off the back of a base cabinet when the back is in plain view.
- Add option code 15 for a welded back on a cabinet.

Extended Finished Back Panels



	Standing Ht.		Sittii	Sitting Ht.		A Ht.
Width	4"Depth	18"Depth	4"Depth	18"Depth	4"Depth	18"Depth
15"	FEA1504	FEA1518	FEC1504	FEC1518	FED1504	FED1518
18"	FEA1804	FEA1818	FEC1804	FEC1818	FED1804	FED1818
24"	FEA2404	FEA2418	FEC2404	FEC2418	FED2404	FED2418
30"	FEA3004	FEA3018	FEC3004	FEC3018	FED3004	FED3018
36"	FEA3604	FEA3618	FEC3604	FEC3618	FED3604	FED3618
42"	FEA4204	FEA4218	FEC4204	FEC4218	FED4204	FED4218
48"	FEA4804	FEA4818	FEC4804	FEC4818	FED4804	FED4818
60"	FEA6004	FEA6018	FEC6004	FEC6018	FED6004	FED6018
72"	FEA7204	FEA7218	FEC7204	FEC7218	FED7204	FED7218

- Extended Finished Back Panels are used to close off the back of a base cabinet when the back is in plain view.
- Panels are available in 1" increments from 4" 18" Deep. Change the last two digits to the dimension required.
- Corner gussets are welded into top and bottom corners.

Sing	le	Ang	les



Double Angles



1" x 1"	2" x 2"	
SAG110C	SAG220C	
SAG110B	SAG220B	
SAG110A	SAG220A	
SAG1118	SAG2218	
SAG1124	SAG2224	
SAG1136	SAG2236	
SAG1148	SAG2248	
SAG1160	SAG2260	
SAG1172	SAG2272	
SAG1184	SAG2284	
	SAG110C SAG110B SAG110A SAG1118 SAG1124 SAG1136 SAG1148 SAG1160 SAG1172	SAG110C SAG220C SAG110B SAG220B SAG110A SAG220A SAG1118 SAG2218 SAG1124 SAG2224 SAG1136 SAG2236 SAG1148 SAG2248 SAG1160 SAG2260 SAG1172 SAG2272

Length	1" x 1"	
28-1/2"	DAG110C	
31-1/2"	DAG110B	
34-1/2"	DAG110A	
18"	DAG1118	
24"	DAG1124	
36"	DAG1136	
48"	DAG1148	
60"	DAG1160	
72"	DAG1172	
84"	DAG1184	





MISCELLANEOUS ITEMS


Front Wall Filler Scribes - 14" Deep



- The Wall Filler Scribe is used to close a space between a wall cabinet and a wall. It can be cut to fit gaps of various sizes.
- · A double angle is provided to enclose the raw edge of the panel.
- Filler panels can be ordered in 1" increments, up to 24" wide. Change the last two digits to the dimension required.
- Use option code "OL" when ordering Architectural Overlay design.

Width	24" High	30" High	36" High	48" High
1"	WFS5301	WFS6301	WFS7301	WF\$8301
2"	WFS5302	WFS6302	WF\$7302	WF\$8302
3"	WFS5303	WFS6303	WF\$7303	WFS8303
4"	WFS5304	WFS6304	WF\$7304	WFS8304
6"	WFS5306	WFS6306	WF\$7306	WF\$8306
8"	WFS5308	WFS6308	WF\$7308	WFS8308
10"	WFS5310	WFS6310	WF\$7310	WFS8310
12"	WFS5312	WFS6312	WF\$7312	WFS8312
15"	WFS5315	WFS6315	WF\$7315	WFS8315
18"	WFS5318	WFS6318	WFS7318	WFS8318
24"	WF\$5324	WF\$6324	WFS7324	WF\$8324

Front Wall Fillers - 14" Deep

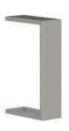

- The Wall Filler Scribe is used to close a space between two wall cabinets.
- Filler panels can be ordered in 1" increments, up to 24" wide. Change the last two digits to the dimension required.
- . Use option code "OL" when ordering Architectural Overlay design.

Width	24" High	30" High	36" High	48" High
1"	WFP5301	WFP6301	WFP7301	WFP8301
2"	WFP5302	WFP6302	WFP7302	WFP8302
3"	WFP5303	WFP6303	WFP7303	WFP8303
4"	WFP5304	WFP6304	WFP7304	WFP8304
6"	WFP5306	WFP6306	WFP7306	WFP8306
8"	WFP5308	WFP6308	WFP7308	WFP8308
10"	WFP5310	WFP6310	WFP7310	WFP8310
12"	WFP5312	WFP6312	WFP7312	WFP8312
15"	WFP5315	WFP6315	WFP7315	WFP8315
18"	WFP5318	WFP6318	WFP7318	WFP8318
24"	WFP5324	WFP6324	WFP7324	WFP8324


Front Wall Filler Scribes - 18" Deep



- The Wall Filler Scribe is used to close a space between a wall cabinet and a wall. It can be cut to fit gaps of various
- A double angle is provided to enclose the raw edge of the panel.
- Filler panels can be ordered in 1" increments, up to 24" wide. Change the last two digits to the dimension required.
- Use option code "OL" when ordering Architectural Overlay design.

Width	24" High	30" High	36" High	48" High
1"	WFS5201	WFS6201	WFS7201	WFS8201
2"	WFS5202	WFS6202	WFS7202	WFS8202
3"	WFS5203	WFS6203	WFS7203	WFS8203
4"	WFS5204	WFS6204	WFS7204	WFS8204
6"	WFS5206	WFS6206	WFS7206	WFS8206
8"	WFS5208	WFS6208	WFS7208	WFS8208
10"	WFS5210	WFS6210	WFS7210	WFS8210
12"	WFS5212	WFS6212	WFS7212	WFS8212
15"	WFS5215	WFS6215	WFS7215	WFS8215
18"	WFS5218	WFS6218	WFS7218	WFS8218
24"	WF\$5224	WFS6224	WFS7224	WFS8224

Front Wall Fillers - 18" Deep

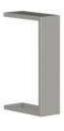

- The Wall Filler Scribe is used to close a space between two wall cabinets.
- Filler panels can be ordered in 1" increments, up to 24" wide. Change the last two digits to the dimension required.
- Use option code "OL" when ordering Architectural Overlay design.


Width	24" High	30" High	36" High	48" High
1"	WFP5201	WFP6201	WFP7201	WFP8201
2"	WFP5202	WFP6202	WFP7202	WFP8202
3"	WFP5203	WFP6203	WFP7203	WFP8203
4"	WFP5204	WFP6204	WFP7204	WFP8204
6"	WFP5206	WFP6206	WFP7206	WFP8206
8"	WFP5208	WFP6208	WFP7208	WFP8208
10"	WFP5210	WFP6210	WFP7210	WFP8210
12"	WFP5212	WFP6212	WFP7212	WFP8212
15"	WFP5215	WFP6215	WFP7215	WFP8215
18"	WFP5218	WFP6218	WFP7218	WFP8218
24"	WFP5224	WFP6224	WFP7224	WFP8224






MISCELLANEOUS ITEMS

Sloping Dust Caps

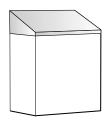


- The Sloping Dust Caps are used as a housekeeping aid on top of wall and floor cabinets.
- . 8-3/8"high.


Overlay Corner Dust Caps

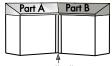

Front Wall Fillers

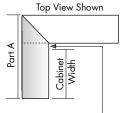

1" Front Wall Filler (WFS or WFP)*

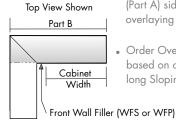

- The Overlay Corner Dust Caps are designed to close the corner space left by two wall cabinets. They are used in conjunction with a Sloping Dust Cap.
- 8-3/8" high.
- Order Overlay Corner Dust Cap based on cabinet width. Extra long sloping dust cap required see below to order.
- Cabinets must be equal depth.
- Order Front Wall Filler (WFP or WFS) in conjunction with an Overlap Corner Dust Cap. When ordering Overlay Corner
 - Dust Caps, always use left hand (Part A) side to be the corner overlaying dust cap.
- Order Overlay Corner Dust Cap based on cabinet width. Extra long Sloping Dust Cap required.

Width 14" Deep 18" Deep 22" Deep 15" SDC0315 SDC0215 SDC0115 18" SDC0318 SDC0218 SDC0118 24" SDC0324 SDC0224 SDC0124 30" SDC0330 SDC0230 SDC0130 36" SDC0336 SDC0236 SDC0136 42" SDC0342 SDC0242 SDC0142 48" SDC0348 SDC0248 SDC0148 60" SDC0360 SDC0260 SDC0160 72" SDC0372 SDC0272 SDC0172 513mm SDC030A SDC020A SDC010A 1000mm SDC030B SDC020B SDC010B

	Part A			
	Cabinet Width	14" Deep	18" Deep	22" Deep
	15"	CDL0315	CDL0215	CDL0115
	18"	CDL0318	CDL0218	CDL0118
	24"	CDL0324	CDL0224	CDL0124
_	30"	CDL0330	CDL0230	CDL0130
)	36"	CDL0336	CDL0236	CDL0136
	42"	CDL0342	CDL0242	CDL0142
	48"	CDL0348	CDL0248	CDL0148
	60"	CDL0360	CDL0260	CDL0160
	72"	CDL0372	CDL0272	CDL0172
	513mm	CDL030A	CDL020A	CDL010A
	1000mm	CDL030B	CDL020B	CDL010B

Part B			
Cabinet Width	14" Deep	18" Deep	22" Deep
15"	CDR0315	CDR0215	CDR0115
18"	CDR0318	CDR0218	CDR0118
24"	CDR0324	CDR0224	CDR0124
30"	CDR0330	CDR0230	CDR0130
36"	CDR0336	CDR0236	CDR0136
42"	CDR0342	CDR0242	CDR0142
48"	CDR0348	CDR0248	CDR0148
60"	CDR0360	CDR0260	CDR0160
72"	CDR0372	CDR0272	CDR0172
513mm	CDR030A	CDR020A	CDR010A
1000mm	CDR030B	CDR020B	CDR010B

Depth	Item Number				
14"	CDC0300				
18"	CDC0200				

- The Corner Wall Cabinet Dust Caps are used as a housekeeping aid on top of Corner Wall Cabinets.
- 8-3/8" high.

		Right Hand			Left Hand	
Width	14"Depth	18"Depth	22"Depth	14"Depth	18"Depth	22"Depth
1"	SDR0301	SDR0201	SDR0101	SDL0301	SDL0201	SDL0101
2"	SDR0302	SDR0202	SDR0102	SDL0302	SDL0202	SDL0102
3"	SDR0303	SDR0203	SDR0103	SDL0303	SDL0203	SDL0103
4"	SDR0304	SDR0204	SDR0104	SDL0304	SDL0204	SDL0104
5"	SDR0305	SDR0205	SDR0105	SDL0305	SDL0205	SDL0105

Corner Wall Cabinet Dust Caps



Sloping Dust Cap Scribes



- The scribe feature allows for easy width adjustment in the field.
- 8-3/8" high.
- Sloping Dust Caps comes with two ECS clips.
- · Right hand shown.





MISCELLANEOUS ITEMS

Light Valances



- · Left side return light valance shown. Right side return light valance is opposite of what is shown.
- 3" high.

	Right Side	e Return	Left Hand	Return
Width	12"Depth	16"Depth	12"Depth	16"Depth
18"	VPR3018	VPR4018	VPL3018	VPL4018
24"	VPR3024	VPR4024	VPL3024	VPL4024
30"	VPR3030	VPR4030	VPL3030	VPL4030
36"	VPR3036	VPR4036	VPL3036	VPL4036
42"	VPR3042	VPR4042	VPL3042	VPL4042
48"	VPR3048	VPR4048	VPL3048	VPL4048
60"	VPR3060	VPR4060	VPL3060	VPL4060
72"	VPR3072	VPR4072	VPL3072	VPL4072
96"	VPR3096	VPR4096	VPL3096	VPL4096



- . Light valance to be fastened to the bottom of a wall cabinet. Available with left side, right side, or both returns.
- Use 12" deep light valance for 14" deep cabinets.
- Use 16" deep light valance for 18" deep cabinets.
- . Light valance returns are 2" shorter in depth than the standard depth of the cabinet, to allow for clearance of opening and closing the doors above.
- . 3" high.

	With Side	e Returns
Width	12"Depth	16"Depth
18"	VP\$3018	VPS4018
24"	VPS3024	VPS4024
30"	VPS3030	VPS4030
36"	VPS3036	VPS4036
42"	VPS3042	VPS4042
48"	VPS3048	VPS4048
60"	VPS3060	VPS4060
72"	VPS3072	VPS4072
96"	VPS3096	VPS4096



- Front only light valance to be fastened to the bottom of a wall cabinet.
- . 3" high.

Width	Front Only	
18"	VPF1018	
24"	VPF1024	
30"	VPF1030	
36"	VPF1036	
42"	VPF1042	
48"	VPF1048	
60"	VPF1060	
72"	VPF1072	
96"	VPF1096	

Wall Cabinet Furring Panels



- Wall Cabinet Furring panels are used to close off the area between the top of a wall cabinet and the ceiling.
- Furring panels have a double scribe angle at top and a 2" wide flange at bottom.
- . 24" high.

Width	14" Deep	18" Deep	22" Deep
18"	FPS1018	FPS2018	FPS3018
24"	FPS1024	FPS2024	FPS3024
30"	FPS1030	FPS2030	FPS3030
36"	FPS1036	FPS2036	FPS3036
42"	FPS1042	FPS2042	FPS3042
48"	FPS1048	FPS2048	FPS3048
60"	FPS1060	FPS2060	FPS3060
72"	FPS1072	FPS2072	FPS3072
96"	FPS1096	FPS2096	FPS3096
1000mm	FPS100B	FPS200B	FPS300B





MISCELLANEOUS ITEMS

Apron Units - 22" Deep



- Aprons provide a work area from a seated position.
- Aprons are available with or without drawers.
- 15" wide drawer.
- 3-3/4" high, three piece adjustable construction.
- Use option code "OL" when ordering Architectural Overlay design.

Width	Plain Front	With Drawers	# of Drawers
24"	AUC0124	AUA0124	1
30"	AUC0130	AUA0130	1
36"	AUC0136	AUA0136	1
42"	AUC0142	AUA0142	1
48"	AUC0148	AUA0148	2
54"	AUC0154	AUA0154	2
60"	AUC0160	AUA0160	2
66"	AUC0166	AUA0166	2
72"	AUC0172	AUA0172	2
513mm	AUC010A	AUA010A	1
1000mm	AUC010B	AUA010B	1

Apron Units 28" Deep



- Aprons provide a work area from a seated position.
- Aprons are available with or without drawers.
- 15" wide drawer.
- 3-3/4" high, three piece adjustable construction.
- Use option code "OL" when ordering Architectural Overlay design.

Width	Plain Front	With Drawers	# of Drawers
24"	AUC0424	AUA0424	1
30"	AUC0430	AUA0430	1
36"	AUC0436	AUA0436	1
42"	AUC0442	AUA0442	1
48"	AUC0448	AUA0448	2
54"	AUC0454	AUA0454	2
60"	AUC0460	AUA0460	2
66"	AUC0466	AUA0466	2
72"	AUC0472	AUA0472	2
513mm	AUC040A	AUA040A	1
1000mm	AUC040B	AUA040B	1

Welded Box Aprons - 22" Deep



- Aprons provide a work area from a seated position.
- · Aprons are available with or without drawers.
- 15" wide drawer.
- 3-3/4" high, one piece welded construction.
- . Use option code "OL" when ordering Architectural Overlay design.

Width	Plain Front	With Drawers	# of Drawers
24"	BUC0124	BUA0124	1
30"	BUC0130	BUA0130	1
36"	BUC0136	BUA0136	1
42"	BUC0142	BUA0142	1
48"	BUC0148	BUA0148	2
54"	BUC0154	BUA0154	2
60"	BUC0160	BUA0160	2
66"	BUC0166	BUA0166	2
72"	BUC0172	BUA0172	2
513mm	BUC010A	BUA010A	1
1000mm	BUC010B	BUA010B	1

Welded Box Aprons - 28" Deep



- Aprons provide a work area from a seated position.
- · Aprons are available with or without drawers.
- 15" wide drawer.
- 3-3/4" high, one piece welded construction.
- Use option code "OL" when ordering Architectural Overlay design.

Width	Plain Front	With Drawers	# of Drawers
24"	BUC0424	BUA0424	1
30"	BUC0430	BUA0430	1
36"	BUC0436	BUA0436	1
42"	BUC0442	BUA0442	1
48"	BUC0448	BUA0448	2
54"	BUC0454	BUA0454	2
60"	BUC0460	BUA0460	2
66"	BUC0466	BUA0466	2
72"	BUC0472	BUA0472	2
513mm	BUC040A	BUA040A	1
1000mm	BUC040B	BUA040B	1





MISCELLANEOUS ITEMS

Apron Service Covers



- Apron Service Cover panels are used to conceal open spaces and services below aprons.
- · Attaches to adjacent cabinets and meets the bottom of the apron.
- · Assembly includes single angles for attachment to adjacent cabinets. Upper panel is removable to allow access to the services. The lower panel stays fixed in order to apply the coved base and overlapping the lower and upper panel allows to accommodate for uneven floors.

Width	Standing Ht.	Sitting Ht.	ADA Ht.
15"	SCP3015	SCP4015	SCP6015
18"	SCP3018	SCP4018	SCP6018
24"	SCP3024	SCP4024	SCP6024
30"	SCP3030	SCP4030	SCP6030
36"	SCP3036	SCP4036	SCP6036
42"	SCP3042	SCP4042	SCP6042
48"	SCP3048	SCP4048	SCP6048
54"	SCP3054	SCP4054	SCP6054
60"	SCP3060	SCP4060	SCP6060
66"	SCP3066	SCP4066	SCP6066
72"	SCP3072	SCP4072	SCP6072
513mm	SCP300A	SCP400A	SCP600A
1000mm	SCP300B	SCP400B	SCP600B

Front Rails





- Front Rails are commonly used to finish a space at a sitting area. They are used when a full frame apron is not required.
- In a "cut to fit" condition, use part number DAG000F for slip filler.

3-3/4" High	2″ High	
FRA0018	FRB0018	
FRA0024	FRB0024	
FRA0030	FRB0030	
FRA0036	FRB0036	
FRA0042	FRB0042	
FRA0048	FRB0048	
FRA0054	FRB0054	
FRA0060	FRB0060	
FRA0066	FRB0066	
FRA0072	FRB0072	
FRA000A	FRB000A	
FRA000B	FRB000B	
	FRA0018 FRA0024 FRA0030 FRA0036 FRA0042 FRA0048 FRA0054 FRA0060 FRA0066 FRA0066 FRA0072 FRA000A	FRA0018 FRB0018 FRA0024 FRB0024 FRA0030 FRB0030 FRA0036 FRB0036 FRA0042 FRB0042 FRA0048 FRB0048 FRA0054 FRB0054 FRA0060 FRB0060 FRA0066 FRB0066 FRA0072 FRB000A

Double Angle Sleeves



Length	Item Number	
3-3/4"	DAG000F	
3"	DAG0003	
2"	DAG0002	

• Double Angle Sleeves are used to cover the cut edge of a "cut to fit" front rail.

Front Rail Service Covers



- Front Rail Service Cover panels are used to conceal open spaces and services behind front rails.
- · Attaches to adjacent cabinets and meets the bottom of the work surface.
- · Assembly includes single angles for attachment to adjacent cabinets. Upper panel is removable to allow access to the services. The lower panel stays fixed in order to apply the coved base and overlapping the lower and upper panel allows to accommodate for uneven floors.

Width	Standing Ht.	Sitting Ht.	ADA Ht.
15"	SCP1015	SCP2015	SCP5015
18"	SCP1018	SCP2018	SCP5018
24"	SCP1024	SCP2024	SCP5024
30"	SCP1030	SCP2030	SCP5030
36"	SCP1036	SCP2036	SCP5036
42"	SCP1042	SCP2042	SCP5042
48"	SCP1048	SCP2048	SCP5048
54"	SCP1054	SCP2054	SCP5054
60"	SCP1060	SCP2060	SCP5060
66"	SCP1066	SCP2066	SCP5066
72"	SCP1072	SCP2072	SCP5072
513mm	SCP100A	SCP200A	SCP500A
1000mm	SCP100B	SCP200B	SCP500B





MISCELLANEOUS ITEMS

Box Turret Supports



- Box Turret Supports are used to support reagent shelves. They may be used with front rails to form a rigid reagent support structure.
- The side panel is removable in order to access
- · Counter top mounting clips are included.
- 11-1/2" deep.

Height	2" Wide	4" Wide	6" Wide
11"	9900062	9900010	9900012
17"	9900063	9900011	9900013
- A 1 - 1 - 1	C 0 0 / 4// L L ED 4	-1	

Notched for a 3-3/4"high FRA rail.

Height	2" Wide	4" Wide	6" Wide
11"	9990062	9990010	9990012
17"	9990063	9990011	9990013

. Notched for a 2"high FRB rail.

Fixed Box Turret Support



- Fixed Box Turret Supports are used to support reagent shelves. It may be used in conjunction with front rails to form a rigid reagent support structure.
- 2" wide.

Height	10" Deep	12" Deep	14" Deep	18" Deep
11"	BRS1010	BRS1012	BRS1014	BRS1018
14"	BRS2010	BRS2012	BRS2014	BRS2018
17"	BRS3010	BRS3012	BRS3014	BRS3018

Reagent Turret Supports



- Reagent Turret Supports are used to support reagent shelves.
- Stainless steel countertop mounting clips are included.
- 1-5/8" wide by 11-1/2" deep.


11" 9900006 17" 9900007 18" 9900008 19" 9900009	Height	Item Number	
18" 9900008	11"	9900006	
	17"	9900007	
19" 9900009	18"	9900008	
	19"	9900009	

Adjustable Box Turret Supports



- Adjustable Box Turret End Supports are used to support reagent shelves. It may be used in conjunction with front rails to form a rigid reagent support structure.
- Adjustable slots are provided on one side.
- 2" wide.

Adjustable Box Turret End Supports				
Height	10" Deep	12" Deep	14" Deep	18" Deep
24"	BRE1010	BRE1012	BRE1014	BRE1018
34"	BRE2010	BRE2012	BRE2014	BRE2018

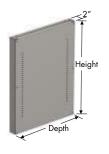

- Adjustable Box Turret Middle Supports are used to support reagent shelves. It may be used in conjunction with front rails to form a rigid reagent support structure.
- · Adjustable slots are provided on both sides.
- 2" wide.

Adjustable Box Turret Middle Supports				
Height	10" Deep	12" Deep	14" Deep	18" Deep
24"	BRM1010	BRM1012	BRM1014	BRM1018
34"	BRM2010	BRM2012	BRM2014	BRM2018





MISCELLANEOUS ITEMS

Support Legs



- Support Legs are used to support island counter tops or aprons where base cabinets are not required.
- 1-1/2" wide.
- Use option code "OL" when ordering Architectural Overlay design.

Depth	Standing Ht.	Sitting Ht.	ADA Ht.
28"	GLD1028	GLD2028	GLD3028
29"	GLD1029	GLD2029	GLD3029
34"	GLD1034	GLD2034	GLD3034
35"	GLD1035	GLD2035	GLD3035
46"	GLD1046	GLD2046	GLD3046
52"	GLD1052	GLD2052	GLD3052
58"	GLD1058	GLD2058	GLD3058



- Support Legs are used to support wall counter tops or aprons where base cabinets are not required.
- Lower portion depth is half of overall depth.
- Support Legs come complete with wall mounting channel.
- 1-1/2" wide.
- Use option code "OL" when ordering Architectural Overlay design.

Depth	Standing Ht.	Sitting Ht.	ADA Ht.
22"	GLA1022	GLA2022	GLA3022
23"	GLA1023	GLA2023	GLA3023
29"	GLA1029	GLA2029	GLA3029
35"	GLA1035	GLA2035	GLA3035

Gable Legs



- Gable Legs are used to support wall counter tops or aprons where base cabinets are not
- Gable Legs come complete with wall mounting channel.
- 1-1/2" wide.
- Use option code "OL" when ordering Architectural Overlay design.

Depth	Standing Ht.	Sitting Ht.	ADA Ht.
22"	GL\$1022	GL\$2022	GL\$3022
23"	GL\$1023	GL\$2023	GL\$3023
29"	GL\$1029	GL\$2029	GL\$3029
35"	GL\$1035	GLS2035	GL\$3035

Support Legs



- Support Legs are used to support counter tops where base cabinets are not required.
- · Leveling bolts and PVC boots are included.
- . 2" X 2" tube.
- Suitable for use with TAC, TAA, TCC, TAB, TCB & TCE Tables.

Table Leg Assembly					
Depth	Standing Ht.	Sitting Ht.	ADA Ht.		
22"	TLA1022	TLA2022	TLA5022		
23"	TLA1023	TLA2023	TLA5023		
28"	TLA1028	TLA2028	TLA5028		
29"	TLA1029	TLA2029	TLA5029		
35"	TLA1035	TLA2035	TLA5035		



Double Table Leg Assembly				
Depth	Standing Ht.	Sitting Ht.	ADA Ht.	
46"	TLA1046	TLA2046	TLA5046	
52"	TLA1052	TLA2052	TLA5052	
58"	TLA1058	TLA2058	TLA5058	





MISCELLANEOUS ITEMS

Table Legs



Height	Item Number	
28-3/4"	3200000	
34-3/4"	3100000	
31-3/4"	3500000	

- Table Legs are used to support counter tops where base cabinets are not required.
- · Leveling bolts and PVC boots are included.
- . 2" x 2' tube.

_		1
.>e	rvice	Leas
00		-093



Depth	Standing Ht.	Sitting Ht.	ADA Ht.	
7"	SLG1007	SLG2007	SLG3007	
11"	SLG1011	SLG2011	SLG3011	
14"	SLG1014	SLG2014	SLG3014	

- Service Legs are used to go behind cabinets when extra counter support is required.
- Service Legs can also be used by plumbers to attach piping within the service chase.
- Each assembly comes with 3/8" spring nuts & bolts.

Under Counter Support Brackets



Item Number

CBS0001

- Under Counter Support Brackets are used to support counter
- Suitable for 30" counter tops.
- Supports are 22" high by 28" deep welded from 1" x 2" 16 gauge steel tube.

Stainless Steel Electrical Pedestal Boxes





Single Gang	Double Gang	
PBS	PBD	


- . Double-faced deck-mounted stainless steel outlet box only.
- Complete with threaded nipple to attach to mounting surface.
- . Cover plate and receptacles are not included.
- . Not CSA or UL approved.






MISCELLANEOUS ITEMS


Full Access Island Service Drops

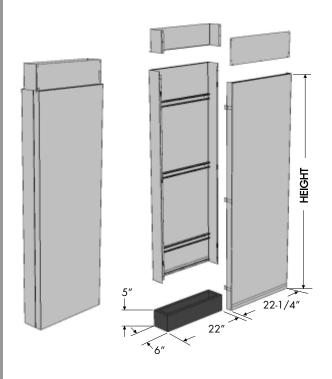


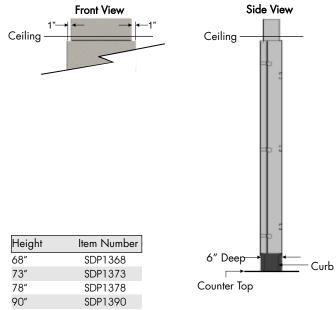

- Use to carry utilities from the counter top up through the ceiling.
- Service Drop offers full height (from counter top to ceiling) access to services.
- Front panel attaches with 'Z' clip.
- Designed to fit on a 22" wide \times 5" high \times 6" deep curb.

Full Access Wall Service Drops

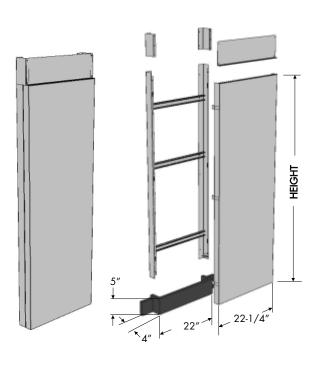

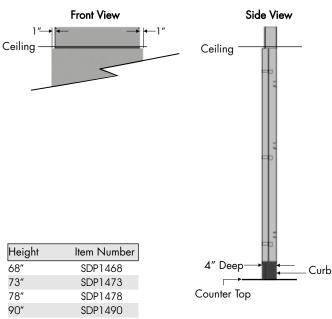


- Use to carry utilities from the counter top up through the ceiling.
- Service Drop offers full height (from counter top to ceiling) access to
- Front panel attaches with 'Z' clip.
- Designed to fit on a 22" wide x 5" high x 4" deep curb.






MISCELLANEOUS ITEMS


Service Drops

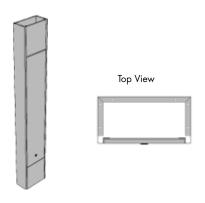

Radius Island Service Drops

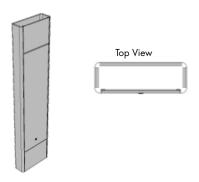

Radius Wall Service Drops

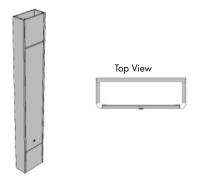

Service Drop Ceiling Trims



	6" Deep		8″ D	еер	
Height	12"Width	24"Width	12″Width	24"Width	
78″	SDP3078	SDP4078	SDP1078	SDP2078	
90"	SDP3090	SDP4090	SDP1090	SDP2090	

- Use to carry utilities from the counter top up through the ceiling, for wall and island applications.
- Each assembly consists of two 12" high fixed upper and lower panels and a lift off center panel.
- 78" high service drop has 54" center lift off panel.
- 90" high service drop has 66" center lift off panel.

Height	12"Width	24"Width
78″	SDP1178	SDP2178
90"	SDP1190	SDP2190

- Use to carry utilities from the counter top up through the ceiling, for island applications.
- Each assembly consists of two 12" high fixed upper and lower panels and a lift off center panel.
- 78" high service drop has 54" center lift off panel.
- 90" high service drop has 66" center lift off panel.
- 6" deep (maximum depth available).

Height	12″Width	24"Width
78″	SDP1278	SDP2278
90"	SDP1290	SDP2290

- Use to carry utilities from the counter top up through the ceiling, for wall applications.
- Each assembly consists of two 12" high fixed upper and lower panels and a lift off center panel.
- 78" high service drop has 54" center lift off panel.
- 90" high service drop has 66" center lift off panel.
- 6" deep (maximum depth available).

Depth	12"Width	24"Width	
6"	EDC1006	EDC2006	
8"	EDC1000	EDC2000	
J	LDCTOOO	LDOZOOO	

- Not for use with Radius and Full access Service Drops.
- Required when Service Drop extends past the ceiling tile.
- Suitable for a 6" and 8" deep Service Drop.





MISCELLANEOUS ITEMS

Cylinder Storage Racks



Width	Item Number	Capable of Holding	
16"	CSR1116	1 Cylinder	
30"	CSR1230	2 Cylinders	
44"	CSR1344	3 Cylinders	
58"	CSR1458	4 Cylinders	
72"	CSR1572	5 Cylinders	

- One cylinder deep storage rack.
- Holds up to 10" diameter cylinders.
- Welded steel rack 30"high x 16"deep.
- Secure to both walls and floor with $\frac{1}{2}$ " steel anchor bolts.
- Non-kink zinc plated chain.



Width	Item Number	Capable of Holding	
16"	CSR2116	2 Cylinder	
30"	CSR2230	4 Cylinders	
44"	CSR2344	6 Cylinders	
58"	CSR2458	8 Cylinders	
72"	CSR2572	10 Cylinders	

- Two cylinder deep storage rack.
- Holds up to 10" diameter cylinders.
- Welded steel rack 30"high x 26"deep.
- Secure to both walls and floor with $\frac{1}{2}$ " steel anchor bolts.
- Non-kink zinc plated chain.

1	
	6 cylinder rack showr

Width	Item Number	Capable of Holding	
16"	CSR3116	3 Cylinder	
30"	CSR3230	6 Cylinders	
44"	CSR3344	9 Cylinders	
58"	CSR3458	12 Cylinders	
72"	CSR3572	15 Cylinders	

- Three cylinder deep storage rack.
- · Holds up to 10" diameter cylinders.
- Welded steel rack 30"high x 38"deep.
- Secure to both walls and floor with $\frac{1}{2}$ " steel anchor bolts.
- Non-kink zinc plated chain.






MISCELLANEOUS ITEMS

Reagent Island Shelving Unit

Reagent Island Shelving Units are designed to provide support for shelving and are 1" wide by 2" deep. Units are available in custom lengths to suit your exact requirements.

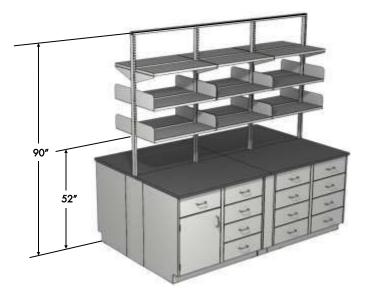


Reagent Uprights



Reagent Upper Rail



Reagent Upper Rail Jointer

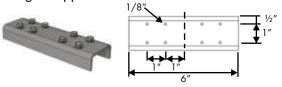

Item Number

RIV0090

- Used to create island assemblies where upright shelving is needed and relies on adjacent fixed casework for support.
- Reagent Uprights are 90" high by 1" wide by 2" deep tubing and are supplied with base support and hardware.
- Uprights are double slotted on both sides allowing shelving to be adjustable in 1" increments.
- Reagent Uprights are only suitable with "RIB" Shelf Brackets.
- Reagent Uprights are not suitable for suspending cabinets or work surfaces.

Item Number

RIH0120

- · Used for island assemblies where upright shelving is needed and relies on adjacent fixed casework for support.
- Reagent Upper Rails are 1" high by 2" deep by 120" long
- Maximum length available without a seam in upper rail is 120"
- · Available in custom lengths to suit your requirements.
- Supplied with two black plastic end caps.

Item Number

- Used for joining two Reagent Upper Rails together.
- Fits inside of 1" high by 2" deep tubing.
- Reagent Upper Rail Jointer is 6" long.
- Supplied with hardware.





MISCELLANEOUS ITEMS

Reagent Island Steel & Wood Shelves



- Shelves suitable for Reagent Island shelving units.
- Shelves are actual widths.
- Add 1" angles to shelves by using option codes:
 - AR (Rear Shelf Angle welded to shelf)
 - AF (Front Shelf Angle welded to shelf)
 - AB (Rear and Front Shelf Angle welded to shelf)

Steel Shelves	;				
Width	10"Deep	12"Deep	14"Deep	16"Deep	18" Deep
24"	RIS1024	RIS1224	RIS1424	RIS1624	RIS1824
30"	RIS1030	RIS1230	RIS1430	RIS1630	RIS1830
36"	RIS1036	RIS1236	RIS1436	RIS1636	RIS1836
42"	RIS1042	RIS1242	RIS1442	RIS1642	RIS1842
48"	RIS1048	RIS1248	RIS1448	RIS1648	RIS1848
60"	RIS1060	RIS1260	RIS1460	RIS1660	RIS1860



- Shelves suitable for Reagent Island shelving units.
- Shelves are actual widths.
- For wood shelving options see page N8.

Wood Shelves						
Width	10"Deep	12"Deep	14"Deep	16"Deep	18" Deep	
24"	RIS1024W	RIS1224W	RIS1424W	RIS1624W	RIS1824W	
30"	RIS1030W	RIS1230W	RIS1430W	RIS1630W	RIS1830W	
36"	RIS1036W	RIS1236W	RIS1436W	RIS1636W	RIS1836W	
42"	RIS1042W	RIS1242W	RIS1442W	RIS1642W	RIS1842W	
48"	RIS1048W	RIS1248W	RIS1448W	RIS1648W	RIS1848W	
60"	RIS1060W	RIS1260W	RIS1460W	RIS1660W	RIS1860W	

Reagent Island Shelf Brackets








Standard profile Radius profile Tapered profile

Depth	Standard Profile	Radius Profile	Tapered Profile	
10"	RIB0010	RIB3010	RIB4010	
12"	RIB0012	RIB3012	RIB4012	
14"	RIB0014	RIB3014	RIB4014	
16"	RIB0016	RIB3016	RIB4016	
18"	RIB0018	RIB3018	RIB4018	

- Brackets only suitable for Reagent components.
- Pair of shelf brackets are 6" high.

Reagent Island Inverted Shelf Brackets

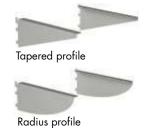

10" RIB1010 RIB1110 12" RIB1012 RIB1112 14" RIB1014 RIB1114 16" RIB1016 RIB1116			
12" RIB1012 RIB1112 14" RIB1014 RIB1114 16" RIB1016 RIB1116	Depth	Tapered Profile	Radius Profile
14" RIB1014 RIB1114 16" RIB1016 RIB1116	10"	RIB1010	RIB1110
16" RIB1016 RIB1116	12"	RIB1012	RIB1112
	14"	RIB1014	RIB1114
	16"	RIB1016	RIB1116
18" RIB1018 RIB1118	18"	RIB1018	RIB1118

- Brackets only suitable for Reagent components.
- Pair of shelf brackets are 4-1/4" high.


Reagent Island Shelf Angles



- · Shelf angles are only suitable for steel RIS shelves.
- 2" high with a 1" return and notched for shelf brackets.

Length	Item Number
20"	RSA2120
24"	RSA2124
30"	RSA2130
36"	RSA2136
42"	RSA2142
48"	RSA2148
60"	RSA2160

Reagent Island Shelf Retainer Rods

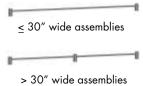

Letainer Rod Parts
ciamor itoa i ario
1 rod & 2 turrets
1 rod & 2 turrets
1 rod & 3 turrets


- 1" high retainer rod (5/16" dia.) Prevents articles from falling off the shelf.
- · Available in stainless steel or powder coated stainless steel and can be added to existing shelves. Ships loose.
- Threaded turrets includes hardware for mounting to shelf (shelf must have holes drilled to accept turrets).
- MI and MIW at the end of the part numbers represent option codes.






MISCELLANEOUS ITEMS

Sigma Flex™ Wall Pilasters

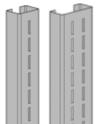

- Single Slotted Pilasters are 13/16" deep by 1-1/4"
- Double Slotted Pilasters are 13/16" deep by 2-1/2"
- Pilasters are adjustable in 1" increments. Accepts Sigma Flex TM shelves.
- Fastens to wall to give support structure where no services are required.
- Ensure wall is of sufficient strength to support intended load (significant blocking required).

Height	Single Slotted	Double Slotted	
18"	FPA1018	FPA0018	
24"	FPA1024	FPA0024	
30"	FPA1030	FPA0030	
36"	FPA1036	FPA0036	
42"	FPA1042	FPA0042	
48"	FPA1048	FPA0048	
60"	FPA1060	FPA0060	
72"	FPA1072	FPA0072	
96"	FPA1096	FPA0096	

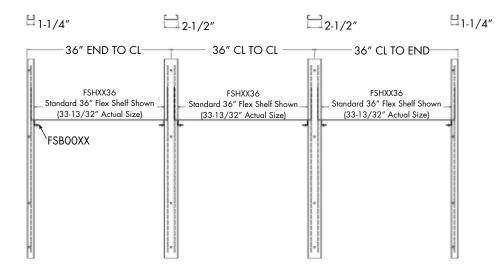

Sigma Flex™ Shelves



- Shelves suitable for Sigma Flex[™] support structure or Sigma Cart $^{\text{TM}}$.
- · Attach to shelf brackets with screws (screws not included).
- . Nominal width.
- For wood shelving options see page N8.
- · Add option code EW for exact width.



Steel Shelves						
Width	6"Deep	8"Deep	10"Deep	12"Deep	14" Deep	18" Deep
24"	FSH0624	FSH0824	FSH1024	FSH1224	FSH1424	FSH1824
30"	FSH0630	FSH0830	FSH1030	FSH1230	FSH1430	FSH1830
36"	FSH0636	FSH0836	FSH1036	FSH1236	FSH1436	FSH1836
42"	FSH0642	FSH0842	FSH1042	FSH1242	FSH1442	FSH1842
48"	FSH0648	FSH0848	FSH1048	FSH1248	FSH1448	FSH1848
54"	FSH0654	FSH0854	FSH1054	FSH1254	FSH1454	FSH1854
60"	FSH0660	FSH0860	FSH1060	FSH1260	FSH1460	FSH1860
66"	FSH0666	FSH0866	FSH1066	FSH1266	FSH1466	FSH1866
72"	FSH0672	FSH0872	FSH1072	FSH1272	FSH1472	FSH1872


Wood	Shelves					
Width	6"Deep	8"Deep	10"Deep	12"Deep	14" Deep	18" Deep
24"	FSH0624W	FSH0824W	FSH1024W	FSH1224W	FSH1424W	FSH1824W
30"	FSH0630W	FSH0830W	FSH1030W	FSH1230W	FSH1430W	FSH1830W
36"	FSH0636W	FSH0836W	FSH1036W	FSH1236W	FSH1436W	FSH1836W
42"	FSH0642W	FSH0842W	FSH1042W	FSH1242W	FSH1442W	FSH1842W
48"	FSH0648W	FSH0848W	FSH1048W	FSH1248W	FSH1448W	FSH1848W
54"	FSH0654W	FSH0854W	FSH1054W	FSH1254W	FSH1454W	FSH1854W
60"	FSH0660W	FSH0860W	FSH1060W	FSH1260W	FSH1460W	FSH1860W





MISCELLANEOUS ITEMS

Sigma Flex™ Shelf Retainer Rods

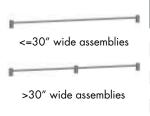

	Steel S	Steel Shelves		Wood Shelves		
Width	Powder Coated	Stainless Steel	Powder Coated	Stainless Steel	Retainer Rod Parts	
24"	PRRO024-FX	SRR0024-FX	PRRO024-FXW	SRR0024-FXW	1 rod & 2 turrets	
30"	PRRO030-FX	SRR0030-FX	PRR0030-FXW	SRR0030-FXW	1 rod & 2 turrets	
36"	PRRO036-FX	SRR0036-FX	PRRO036-FXW	SRR0036-FXW	1 rod & 3 turrets	
42"	PRRO042-FX	SRR0042-FX	PRRO042-FXW	SRR0042-FXW	1 rod & 3 turrets	
48"	PRRO048-FX	SRR0048-FX	PRRO048-FXW	SRR0048-FXW	1 rod & 3 turrets	
54"	PRRO054-FX	SRR0054-FX	PRR0054-FXW	SRR0054-FXW	1 rod & 3 turrets	
60"	PRRO060-FX	SRR0060-FX	PRR0060-FXW	SRR0060-FXW	1 rod & 3 turrets	
66"	PRROO66-FX	SRR0066-FX	PRR0066-FXW	SRR0066-FXW	1 rod & 4 turrets	
72"	PRRO072-FX	SRR0072-FX	-	-	1 rod & 4 turrets	

- 1" high retainer rod (5/16" dia.) Prevents articles from falling off the shelf.
- Available in stainless steel or powder coated stainless steel and can be added to existing shelves.
- Threaded turrets includes hardware for mounting to shelf (shelf must have holes drilled to accept turrets).
- . Ships loose.

Depth

8"

10"

12"

14"

18"

• FX and FXW at the end of the part numbers represent option codes.

Item Number FIB0006

FIB0008

FIB0010

FIB0012

FIB0014

FIB0018

Sigma Flex™ Top Shelf Brackets



- Pair of top shelf brackets can be used for shelves of different materials.
- Adjustable in 1" increments.
- Only compatible with Sigma Flex™ or Sigma Cart $^{\text{TM}}$.

Sigma Flex™ Lower Shelf Brackets



- Pair of lower shelf brackets can be used for shelves of different materials.
- . Adjustable in 1" increments.
- . Only compatible with Sigma $\mathsf{Flex}^\mathsf{TM}$ or Sigma CartTM.

Depth	Item Number
6"	FSB0006
8"	FSB0008
9"	FSB0009
10"	FSB0010
12"	FSB0012
14"	FSB0014
18"	FSB0018

Sigma Flex™ Radius Top Shelf Brackets



- · Pair of radius shelf brackets can be used for shelves of different materials.
- · Adjustable in 1" increments.
- · Only compatible with Sigma Flex™ or Sigma Cart™.

Depth	Item Number
6"	FIB3006
8"	FIB3008
10"	FIB3010
12"	FIB3012
14"	FIB3014
18"	FIB3018

Sigma Flex™ Radius Outer Shelf Brackets



- · Pair of radius outer shelf brackets can be used for shelves of different
- · Adjustable in 1" increments.
- Only compatible with Sigma Flex™ or Sigma CartTM.

Item Number
FSB3006
FSB3008
FSB3009
FSB3010
FSB3012
FSB3014
FSB3018

Sigma Flex™ Tapered Top Shelf Brackets



- Pair of tapered shelf brackets can be used for shelves of different materials.
- Adjustable in 1" increments.
- Only compatible with Sigma $Flex^{TM}$ or Sigma Cart TM .

Depth	Item Number
6"	FIB4006
8"	FIB4008
10"	FIB4010
12"	FIB4012
14"	FIB4014
18"	FIB4018

Sigma Flex[™] Tapered Outer Shelf Brackets



- · Pair of tapered outer shelf brackets can be used for shelves of different materials.
- · Adjustable in 1" increments.
- . Only compatible with Sigma Flex $^{\text{TM}}$ or Sigma $Cart^{TM}$.

Item Number
FSB4006
FSB4008
FSB4009
FSB4010
FSB4012
FSB4014
FSB4018






MISCELLANEOUS ITEMS

Sigma Frame™ Wall Pilasters

- Pilasters only suitable for Sigma Frame™, Inverted Shelf Brackets, "SBP" and "WPA" components.
- 1-5/8" wide by 13/16" deep.
- Pilasters are adjustable in 1" increments.
- Fastens to wall to give support structure where no services are
- Ensure wall is of sufficient strength to support intended load (significant blocking required).

Height	Single Slotted	Double Slotted	
18"	WPA1018	WPA0018	
24"	WPA1024	WPA0024	
30"	WPA1030	WPA0030	
36"	WPA1036	WPA0036	
42"	WPA1042	WPA0042	
48"	WPA1048	WPA0048	
60"	WPA1060	WPA0060	
72"	WPA1072	WPA0072	
96"	WPA1096	WPA0096	

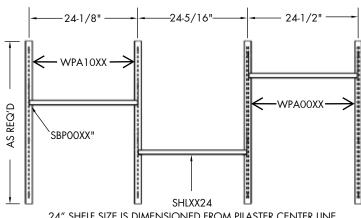

24" SHELF SIZE IS DIMENSIONED FROM PILASTER CENTER LINE * 24" LONG SPACING EXAMPLE SHOWN

Sigma Frame™ Steel & Wood Shelves For Wall Pilasters



- Shelves are actual widths.
- Shelves suitable for Sigma Frame™.
- Add 1" angles to shelves by using option codes:
 - AR (Rear Shelf Angle welded to shelf)
 - AF (Front Shelf Angle welded to shelf)
 - AB (Rear and Front Shelf Angle welded to shelf)

1	

- Shelves are actual widths.
- Shelves suitable for Sigma Frame $^{\text{TM}}$.
- Shelves are 1" thick.
- For wood shelving options see page N8.

Steel Shelves					
Width	10"Deep	12"Deep	14"Deep	16"Deep	18" Deep
24"	SHL1024	SHL1224	SHL1424	SHL1624	SHL1824
30"	SHL1030	SHL1230	SHL1430	SHL1630	SHL1830
36"	SHL1036	SHL1236	SHL1436	SHL1636	SHL1836
42"	SHL1042	SHL1242	SHL1442	SHL1642	SHL1842
48"	SHL1048	SHL1248	SHL1448	SHL1648	SHL1848
60"	SHL1060	SHL1260	SHL1460	SHL1660	SHL1860
513mm	SHL100A	SHL120A	SHL140A	SHL160A	SHL180A
1000mm	SHL100B	SHL120B	SHL140B	SHL160B	SHL180B
1512mm	SHL100C	SHL120C	SHL140C	SHL160C	SHL180C
2000mm	SHL100D	SHL120D	SHL140D	SHL160D	SHL180D


Wood Shelv	es				
Width	10"Deep	12"Deep	14"Deep	16"Deep	18" Deep
24"	SHL1024W	SHL1224W	SHL1424W	SHL1624W	SHL1824W
30"	SHL1030W	SHL1230W	SHL1430W	SHL1630W	SHL1830W
36"	SHL1036W	SHL1236W	SHL1436W	SHL1636W	SHL1836W
42"	SHL1042W	SHL1242W	SHL1442W	SHL1642W	SHL1842W
48"	SHL1048W	SHL1248W	SHL1448W	SHL1648W	SHL1848W
60"	SHL1060W	SHL1260W	SHL1460W	SHL1660W	SHL1860W
513mm	SHL100AW	SHL120AW	SHL140AW	SHL160AW	SHL180AW
1000mm	SHL100BW	SHL120BW	SHL140BW	SHL160BW	SHL180BW
1512mm	SHL100CW	SHL120CW	SHL140CW	SHL160CW	SHL180CW





MISCELLANEOUS ITEMS

Sigma Frame™ Shelf Retainer Rods

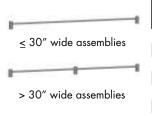

	Steel S	helves		Wood S	helves	
Width	Powder Coated	Stainless Steel		Powder Coated	Stainless Steel	Retainer Rod Parts
24"	PRRO024-MI	SRR0024-MI		PRR0024-MIW	SRR0024-MIW	1 rod & 2 turrets
30"	PRRO030-MI	SRR0030-MI		PRR0030-MIW	SRR0030-MIW	1 rod & 2 turrets
36"	PRRO036-MI	SRR0036-MI		PRR0036-MIW	SRR0036-MIW	1 rod & 3 turrets
42"	PRRO042-MI	SRR0042-MI		PRR0042-MIW	SRR0042-MIW	1 rod & 3 turrets
48"	PRRO048-MI	SRR0048-MI		PRRO048-MIW	SRR0048-MIW	1 rod & 3 turrets
60"	PRRO060-MI	SRR0060-MI		PRR0060-MIW	SRR0060-MIW	1 rod & 3 turrets
513mm	PRROOOA-MI	SRR000A-MI		PRROOOA-MIW	SRR000A-MIW	1 rod & 2 turrets
1000mm	PRROOOB-MI	SRROOOB-MI		PRROOOB-MIW	SRROOOB-MIW	1 rod & 3 turrets
1512mm	PRROOOC-MI	SRR000C-MI		PRROOOC-MIW	SRR000C-MIW	1 rod & 3 turrets
2000mm	PRROOOD-MI	SRR000D-MI		PRROOOD-MIW	SRR000D-MIW	1 rod & 4 turrets
			-		1 14	

- 1" high retainer rod (5/16" dia.) Prevents articles from falling off the shelf.
- Available in stainless steel or powder coated stainless steel and can be added to existing shelves.
- Threaded turrets includes hardware for mounting to shelf (shelf must have holes drilled to accept turrets).
- Ships loose.
- MI and MIW at the end of the part numbers represent option codes.

Sigma Frame™ Shelf Brackets For Wall Pilasters



Standard profile




Radius profile



Depth Standard Profile Radius Profile Tapered Profile 10" SBP0010 SBP3010 SBP4010 12" SBP0012 SBP3012 SBP4012 14" SBP0014 SBP3014 SBP4014 16" SBP0016 SBP3016 SBP4016 18" SBP0018 SBP3018 SBP4018

- Brackets only suitable for Sigma Frame $^{\text{TM}}$ and "WPA" components.
- Pair of shelf brackets are 6" high.
- Use option code NN to add alignment tabs on SBP brackets.

Sigma Frame™ Inverted Shelf Brackets For Wall Pilasters

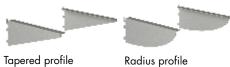




Brackets only suitable for Reagent components.

Pair of shelf brackets are 4-1/4" high.

Depth	Tapered Profile	Radius Profile
10"	SBP1010	SBP1110
12"	SBP1012	SBP1112
14"	SBP1014	SBP1114
16"	SBP1016	SBP1116
18"	SBP1018	SBP1118



888.635.2080 | newenglandlab.com

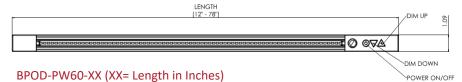

Brightedge Portable - Premium Quality Magnetic LED Light Fixture

- Delivers a minimum average of 50 fc to work surface located 20" below mounting position with uniform illumination
- Maximum wattage: 4.5 watts per linear foot
- Light output from fixture: 435 lumens per linear foot
- Incremental dimming
- Minimum Color Rendering Index (CRI): 95
- Correlated Color temperature to be 3500K, color shift over life of fixture to be maximum +/- 100K, other light temperatures available
- Minimum estimated life to 70% (TM-21) of light output: > 50,000 hours
- When motion is not detected for 30 minutes, the light shall go off but turn on again with motion detection
- If motion is not sensed within 60 minutes, the motion sensor shuts off, requiring the physical switch to turn the light back on
- White magnetic wire management clips, white cords and white power supply
- Standard Lengths: 23", 29", 35", 38", 50", 62" with additional lengths available
- · Light fixtures can be interconnected up to 60 Watts and have either independent or connected controls








MISCELLANEOUS ITEMS

PEGBOARDS

Epoxy Pegboards

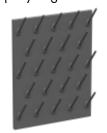

- Black epoxy, 1" thick construction.
- Accepts 2" and 4" deep drip troughs. Must be ordered separately.
- Add option code "15" for finished back.
- 5" long, 3/8" diameter, white polypropylene pegs included. For black pegs use option code "PB" For gray pegs use option code "PG".

Phenolic Pegboards



- Black phenolic, 1" thick construction.
- Accepts 2" and 4" deep drip troughs. Must be ordered separately.
- 5" long, 3/8" diameter, white polypropylene pegs included. For black pegs use option code "PB" For gray pegs use option code "PG".

Stainless Steel Pegboards



- 20 gauge, type 304 stainless steel with #4 finish.
- Add option "ZR" for "Z" wall mounting rail.
- Provided with a 2" or 4" deep drip trough integral to the pegboard.
- Add option code "15" for stainless steel finished back.
- 5" long, 3/8" diameter, white polypropylene pegs included. For black pegs use option code "PB" For gray pegs use option code "PG".

Width & Height	Ероху	Phenolic	Stainless Steel With 2" Deep Trough	Stainless Steel With 4" Deep Trough	Number Of Peas
18" X 24"	EPB1824	PPB1824	SPC1824	SPB1824	23
22" X 30"	EPB2230	PPB2230	SPC2230	SPB2230	33
22" X 36"	EPB2236	PPB2236	SPC2236	SPB2236	39
24" X 30"	EPB2430	PPB2430	SPC2430	SPB2430	33
24" X 36"	EPB2436	PPB2436	SPC2436	SPB2436	39
30" X 36"	EPB3036	PPB3036	SPC3036	SPB3036	53
36" X 30"	EPB3630	PPB3630	SPC3630	SPB3630	51



Drip Troughs



Drain Screen Insert For Drip Troughs



Peg Color	5" Long	
White	PEGW105	
Black	PEGB105	
Gray	PEGG105	

• 5" long, 3/8" diameter.

Width	2" Deep	4" Deep	
18"	DTC0218	DTC0418	
22"	DTC0222	DTC0422	
24"	DTC0224	DTC0424	
30"	DTC0230	DTC0430	
36"	DTC0236	DTC0436	

- 2" and 4" deep drip troughs with 3/8" drain hole. 3 feet of clear PVC drain hose included.
- 20 gauge, type 304 stainless steel with a #4 finish.
- Mounts to the underside of pegboard.

Width	2" Deep	4" Deep	
18"	DSI0218	DSI0418	
22"	DSI0222	DSI0422	
24"	DSI0224	DSI0424	
30"	DSI0230	DSI0430	
36"	DSI0236	DSI0436	

- 2" and 4" deep drain screen inserts for drip troughs.
- 20 gauge, type 304 stainless steel with a #4 finish.





MISCELLANEOUS ITEMS

PEGBOARD ACCESSORIES

Pegboard Hanging Brackets For Sigma Flex™ Wall Pilasters



Item Number

BPB0000

- \bullet Pair of hanging brackets for pegboards, only compatible with Sigma Flex $^{\text{TM}}$ Wall Pilasters.
- Attaches with fasteners, mounting hardware included.
- Ships loose.

Sigma Flex™ Wall Pilasters





Height	Single Slotted	Double Slotted	
18"	FPA1018	FPA0018	
24"	FPA1024	FPA0024	
30"	FPA1030	FPA0030	
36"	FPA1036	FPA0036	
42"	FPA1042	FPA0042	
48"	FPA1048	FPA0048	
60"	FPA1060	FPA0060	
72"	FPA1072	FPA0072	
96"	FPA1096	FPA0096	

- Single Slotted Pilasters are 13/16" deep by 1-1/4" wide.
- Double Slotted Pilasters are 13/16" deep by 2-1/2" wide.
- Pilasters are adjustable in 1" increments.
- Fastens to wall to give support structure where no services are required.
- Ensure wall is of sufficient strength to support intended load (significant blocking required).





MISCELLANEOUS ITEMS

METAL FRAMING SYSTEM

Continuous Slot Metal Framing System


- Installs quickly
- Completely adjustable
- Re-usable
- Eliminates all welding, drilling and use of special tools



1. Place the spring nut anywhere along the continuous slotted channel.




2. Turn clockwise 90° to line up the grooves in the nut with the inturned edges of the channel.

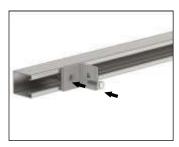

3. Inset the bolt through the fitting and into the spring nut.



4. Added channels are able to be bolted to the existing fitting by following procedure described in steps 1-3.

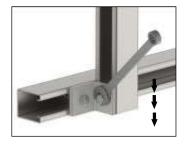

5. Using a wrench to tighten the bolt will lock the serrated teeth of the channel nut onto the inturned edges of the channel.






MISCELLANEOUS ITEMS

METAL FRAMING SYSTEM

1-5/8" x 1-5/8" Channel




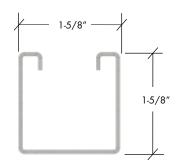

Channels are accurately and carefully cold rolled formed to size from low-carbon strip steel. Cold rolled and stainless steel are offered in 12 gauge.


1-5/8" x 13/16" Channel



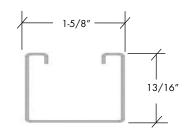


Channels are accurately and carefully cold rolled formed to size from low-carbon strip steel. Cold rolled and stainless steel are offered only 16 gauge.

Channel Details

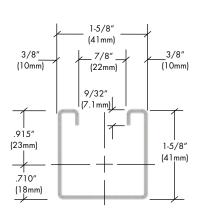

- Wt/100 Ft: 378 Lbs (281 kg/100m).
- Allowable Moment 5,070 In-Lbs (570 N*m).
- 12 gauge Nominal Thickness .105" (2.7mm).


Channels 1-5/8" x 1-5/8"

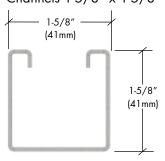

Length	12 Gauge
24"	SMF1224
48"	SMF1248
72"	SMF1272
96"	SMF1296
120"	SMF1220

- Width: 1-5/8"
- Height: 1-5/8"
- Available in powder coated steel or stainless steel.

Channels 1-5/8" x 13/16"

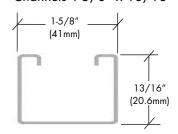

Length	16 Gauge
24"	SMF1624
48"	SMF1648
72"	SMF1672
96"	SMF1696
120"	SMF1620

- . Width: 1-5/8"
- . Height: 13/16"
- · Available in powder coated steel or stainless steel.





MISCELLANEOUS ITEMS

METAL FRAMING SYSTEM

Channel Nuts

Channel nuts are machined from steel bars and are case hardened. Nuts are designed with shaped ends that allow them to be turned clockwise within the channel after insertion through the opening. The channel nut is held in place when the two grooves with the serrated teeth engage the inturns of the channel after the bolt is tightened. This stops any movement of the nut and bolt within the framing member. Channel nuts conform to ASTM Specification A1011 SS GR 33 (material only). Screws conform to SAE J429 GR.





Miscellaneous Hardware









Item Number

SMFN100

- •Thread: 3/8"-16 • Thickness: 1/4"
- For use with 1-5/8" X 13/16" Channel
- Box of 100 pcs.

Item Number

SMFN200

- Thread: 3/8"-16
- Thickness: 3/8"
- For use with 1-5/8" X 1-5/8" Channel
- Box of 100 pcs.

Hex Head Bolt 1"

SMFB100

Hex Head Bolt 1-1/4"

SMFB200

- •Thread 3/8"-16
- ·Box of 100 pcs.

Lock Washer

SMFLW00

- Size: 3/8"
- Box of 100 pcs.

Flat Washer

SMFFW00

- Size: 3/8"
- · Box of 100 pcs.

Hex Nut

SMFHN00

- . Size: 3/8"
- . Box of 100 pcs.

Threaded Rod

SMFTR00

- •Thread:3/8"-16
- Length: 6'
- Electro-galvanized






MISCELLANEOUS ITEMS

METAL FRAMING SYSTEM

Note: No modifications are possible on flat plate and 90° angle fittings.

Flat Plate Fittings

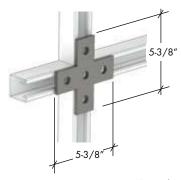

Item Number

SMFFP00


- .Thickness: 1/4"
- . Hole Diameter: 9/16" DIA.



Item Number

SMFFP03

- .Thickness: 1/4"
- . Hole Diameter: 9/16" DIA.

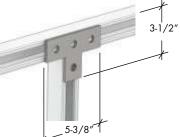

Item Number

SMFFP01

- .Thickness: 1/4"
- . Hole Diameter: 9/16" DIA.



Item Number

SMFFP04

- .Thickness: 1/4"
- . Hole Diameter: 9/16" DIA.



Item Number

SMFFP02

- .Thickness: 1/4"
- . Hole Diameter: 9/16" DIA.




Item Number

SMFFP05


- .Thickness: 1/4"
- . Hole Diameter: 9/16" DIA.

90° Angle Fittings

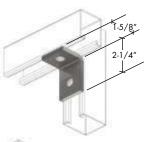

Two Hole Corner 1-5/8" x 2-1/4"


- SMFAF00
- .Thickness: 1/4"
- . Width: 1-5/8"
- . Hole Diameter: 9/16"
- . Hole Spacing: 1-7/8"
- .Edge Distance: 13/16"

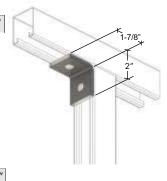

Two Hole Corner 1-7/8" x 2"

SMFAF10

- .Thickness: 1/4"
- . Width: 1-5/8"
- . Hole Diameter: 9/16"
- . Hole Spacing: 1-7/8"
- .Edge Distance: 13/16"

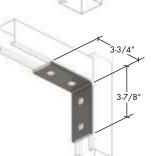

Four Hole Corner 3-3/4" x 3-7/8"

SMFAF20

- .Thickness: 1/4"
- .Width: 1-5/8"
- .Hole Diameter: 9/16"
- . Hole Spacing: 1-7/8"
- .Edge Distance: 13/16"







MISCELLANEOUS ITEMS

METAL FRAMING SYSTEM

Shelf Brackets

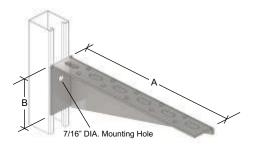


Right hand shown

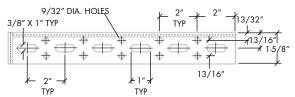

Length A Height B Right Hand Bracket Left Hand Bracket 12" 3-7/16" SMFRH12 SMFLH12 14" 3-15/16" SMFRH14 SMFLH14 16" 4-7/16" SMFRH16 SMFLH16 18" 4-15/16" SMFRH18 SMFLH18 20" 5-7/16" SMFRH20 SMFLH20 22" 5-15/16" SMFLH22 SMFRH22

- Thickness: 12 gauge
- Number of Slots = A/2
- Number of Rows of Holes = (A/2) 1
- No modification are possible.


Hat Channel Shelf Brackets For Wall Struts

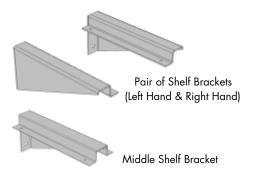

Depth	Pair of Brackets	Middle Bracket	
10"	TBP0010	TBD0010	
12"	TBP0012	TBD0012	
14"	TBP0014	TBD0014	
16"	TBP0016	TBD0016	
18"	TBP0018	TBD0018	

- Steel shelf brackets suitable for Wall Struts.
- 5" high.

Shelf Brackets For Wall Struts

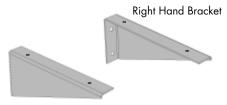

Left	Hand	Bracket
LCII	i iuiiu	DIUCKEI

Depth	Left Hand Bracket	Right Hand Bracket	
10"	TSL0010	TSR0010	
12"	TSL0012	TSR0012	
14"	TSL0014	TSR0014	
16"	TSL0016	TSR0016	
18"	TSL0018	TSR0018	

- Steel shelf brackets suitable for Wall Struts.
- 5" high.






MISCELLANEOUS ITEMS

METAL FRAMING SYSTEM

Note: No modifications are possible on post bases.

Post Bases

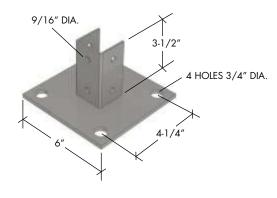

Item Number

SMFPB10


- .Thickness: 1/4"
- . Hole Diameter: 3/4" DIA.

Item Number

SMFPB20



. Hole Diameter: 3/4" DIA.

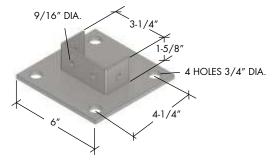


9/16" DIA 4 HOLES 3/4" DIA.

Item Number

SMFPB30

.Thickness: 1/4"

. Hole Diameter: 3/4" DIA.

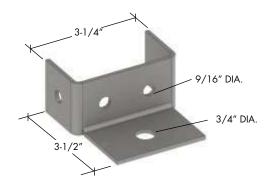

Item Number

SMFPB40

.Thickness: 1/4"

. Hole Diameter: 9/16" DIA.





MISCELLANEOUS ITEMS OPTIONS

Option Code	Description
00	Modified door/drawer configuration
01	Full extension drawer runners
02	Label holders on all doors/drawers
03	Numeric or alpha identification disks on all doors/drawers (customer details required)
07	Cut-outs for electrical and/or plumbing services (customer drawing required)
09	Stainless steel waste flap
10	One fixed drawer divider per drawer (centered in drawer)
11	Multiple drawer dividers per drawer (customer drawing required)
12	Pull out writing board - all cold rolled steel*
13	Pull out writing board - cold rolled steel & high pressure laminate writing surface*
14	Ventilation grill in toe kick base
16	Ship handles not installed
17	Nylon levelers
19	Replaces standard drawer(s) with full width drawer(s) on 5" high apron
26	200 lbs full extension drawer runners
48	Keyboard tray (replaces drawer)
79	Other customer specified base or toe kick
FCW	Front corners welded
OL	Overlay design
PFE	Pull out writing board on full extension runners*
PW	Pre-wired (option code not to be used for fume hoods, customer details required)
SK	Skid item
S01	Full extension soft close drawer runners
S26	2001bs full extension soft close drawers runners
WF	Wood fronts - see specifications
WO	Wood fronts provided and installed by others