


## OVERHEAD CEILING PANEL

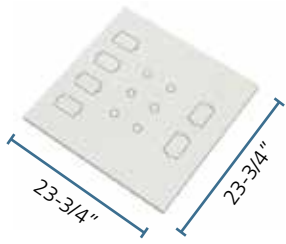

### service panel

12 services - OSP0002

- takes the place of a standard ceiling tile
- ability to mount/disconnect quick connect service fixtures, electrical and data outlets
- 6 service access holes and 6 electrical cutouts
- panel is compatible with most T grid acoustical suspended ceiling structures
- fixtures sold separately
- color: bright white only

BW

## UNDER COUNTER SUPPORT BRACKET



### steel bracket

29" deep - CBS0001

- reversible to support 30" d or 24" d work surfaces
- legs are 22" by 28" d welded from 1" x 2" 16-gauge steel tube
- recommended at 36" maximum intervals for most work surfaces
- color: warm gray and bright white

BW

## WALL SHELVING UPRIGHTS & BOOKEND BRACKETS



### double slotted uprights

36" high - WPA0036

- 1-5/8" w x 13/16" d
- slot pattern suitable for bookend brackets below
- pilasters should be mounted 1/2" on center
- color: warm gray and bright white

BW



### pair of brackets

12" deep - SBP0012

### inverted pair of brackets

12" deep - SBP1012

- pair of shelf brackets are 6" h
- color: warm gray or bright white

BW

## GABLE LEGS



### gable leg 34.75" high

22" deep - GLS1022 <sup>BW</sup>

29" deep - GLS1029 <sup>BW</sup>

### gable leg 28.75" high

22" deep - GLS2022

29" deep - GLS2029

- supports work surfaces or aprons where base cabinets are not required
- complete with wall mounting channel
- 1.5" w
- color: warm gray, limited pieces in bright white

## WALL AND REAGENT SHELVING



### plastic laminate

27" long - PLM1227FS

30" long - PLM1230FS

33" long - PLM1233FS

36" long - PLM1236FS

### phenolic resin

27" long - PHN1227FS

30" long - PHN1230FS

33" long - PHN1233FS

36" long - PHN1236FS

### wood veneer\*

27" long - WD1227FS

30" long - WD1230FS

33" long - WD1233FS

36" long - WD1236FS

\*specify veneer - wood shelves are made to order and may require additional lead time

- 12" deep x 1" thick shelving to span a maximum of 36"
- for wall and reagent shelving only
- finished edge on all sides
- color: black (other colors can be readily fabricated; please consult your representative for lead time)

plastic laminate and phenolic are now available in gray!

please note: gray phenolic will have a black edge (see below). Gray phenolic with gray core can be purchased for an additional fee.



## REAGENT UPRIGHTS



### slotted post assembly

90" high - RIV0090

- slot pattern suitable for bookend brackets
- adjustable in height on 1" increments
- should be mounted 1/2" on center
- 1" w x 2" d; double-sided, double-slotted post with threaded bolt and nut at top
- adjustable shelving slots begin above the work surface at approximately 16" above bench height
- each post sold with base support

BW

## REAGENT CROSSBARS



### crossbar tube

120" length - RIH0120

- 1" x 2" tube with no penetrations for posts
- field drill holes to locate post bolts
- crossbar comes with 2 plastic end caps
- tube insert for seaming lengths over 120" sold separately

### seam connector insert

insert - RIJ0000

- tube insert for seaming lengths over 120"

BW



### reagent shelving

construction: painted steel double-sided, double-slotted post with threaded top-bolt, floor mounting gusset and crossbar

#### installation:

- 1) mount gusset to floor
- 2) adjust post to be plumbed and level
- 3) cut crossbar to length
- 4) drill holes in crossbar for attachment to threaded top-bolts
- 5) secure assembly to base cabinets below
- 6) if possible, secure ends to walls or pipe chases for extra support

color: warm gray

spacing: 27-1/2", 30-1/2", 33-1/2" or 36-1/2" on center for 27", 30", 33" and 36" shelves respectively

shelving: p-lam or phenolic shelving and brackets (see pg. 14)

## SUPER ERECTA



### 18" deep 4-shelf unit

36" wide - M-EZ1836NC-4  
48" wide - M-EZ1848NC-4  
60" wide - M-EZ1860NC-4

### 24" deep 4-shelf unit

36" wide - M-EZ2436NC-4  
48" wide - M-EZ2448NC-4  
60" wide - M-EZ2460NC-4

- includes 4 shelves and 4 stationary posts with cap and leveling bolt
- chrome frame, posts and shelves
- open wire shelf design minimizes dust accumulation and allows free circulation of air

### metro shelving

800 lb. capacity per shelf up to 48" wide

600 lb. capacity per shelf over 48" wide

2,000 lb. maximum capacity per unit

all units 74" tall

easy assembly

nsf listed



## METROMAX i POLYMER



### 18" deep 5-shelf unit

24" wide - M-5X317GX3  
30" wide - M-5X327GX3  
36" wide - M-5X337GX3  
42" wide - M-5X347GX3  
48" wide - M-5X357GX3  
60" wide - M-5X367GX3  
72" wide - M-5X377GX3

### 24" deep 5-shelf unit

36" wide - M-5X537GX3  
42" wide - M-5X547GX3  
48" wide - M-5X557GX3  
60" wide - M-5X567GX3  
72" wide - M-5X577GX3

- includes 5 shelves and 4 stationary posts with cap and leveling bolt
- polymer frame, posts and shelves
- easy-to-clean removable polymer shelf mats
- Microban protection built into shelf mats, frames and posts



\*MICROBAN® and the MICROBAN® symbol are registered trademarks of the Microban Products Company, Huntensville, NC

## SUPER ADJUSTABLE 2



### 18" deep 5-shelf unit

36" wide shelving unit  
48" wide shelving unit  
60" wide shelving unit

### 24" deep 5-shelf unit

48" wide shelving unit  
60" wide shelving unit  
72" wide shelving unit

- includes 5 shelves and 4 stationary posts with cap and leveling bolt
- open-wire design minimizes dust accumulation and allows free circulation of air
- unique corner release system that allows shelves to be adjusted without tools
- type 304 stainless steel wire shelves and posts

