
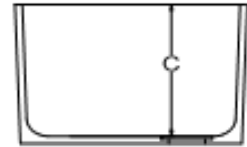


# Epoxy Resin Sinks:

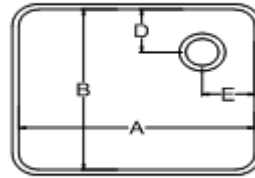

Lipped/Drop in sinks:

SINK	Wt. lbs.	"A" DIM	"B" DIM	"C" DIM	"D" DIM	"E" DIM	DRAIN
AESL-01C	12	9.00	6.00	6.00	3.00	4.50	CENTER
AESL-5	22	14.00	10.00	6.00	5.00	4.00	END
AESL-5C	22	14.00	10.00	6.00	5.00	7.00	CENTER
AESL-6C	32	12.00	12.00	12.00	6.00	6.00	CENTER
AESL-10C	22	16.00	8.00	7.00	4.00	8.00	CENTER
AESL-10E	22	16.00	8.00	7.00	4.00	4.00	END
AESL-13C	15	12.00	8.00	6.00	4.00	6.00	CENTER
AESL-15	30	16.00	12.00	8.00	6.00	4.00	END
AESL-15C	30	16.00	12.00	8.00	6.00	8.00	CENTER
AESL-15RC	30	16.00	12.00	8.00	4.00	4.00	CORNER
AESL-20	39	16.00	16.00	7.50	4.00	4.00	CORNER
AESL-22C	19	18.00	6.50	6.00	3.25	9.00	CENTER
AESL-23C	60	16.00	16.00	11.625	8.00	8.00	CENTER
AESL-24C	35	18.00	14.00	10.50	7.00	9.00	CENTER
AESL-25	39	18.00	15.00	8.00	4.00	4.00	CORNER
AESL-30	53	18.00	15.00	11.00	4.00	4.00	CORNER
AESL-30C	53	18.00	15.00	11.00	7.50	9.00	CENTER
AESL-30D *	130	37.625	15.00	11.00	4.00	4.00	CORNER
AESL-33E	64	21.00	17.00	10.00	8.50	4.00	END
AESL-45	70	21.375	15.375	10.875	4.00	4.00	CORNER
AESL-50C	53	24.00	16.00	8.00	8.00	12.00	CENTER
AESL-52	70	24.00	18.00	11.00	4.00	4.00	CORNER
AESL-54	53	25.00	15.00	8.00	4.00	4.00	CORNER
AESL-55	61	25.00	15.00	10.00	4.00	4.00	CORNER
AESL-57	80	25.00	15.00	13.63	4.00	4.00	CORNER
AESL-58	85	25.00	15.00	17.63	4.00	4.00	CORNER
AESL-59	68	28.00	15.00	12.00	4.00	4.00	CORNER
AESL-61	105	30.00	16.00	17.75	4.50	4.50	CORNER
AESL-65C	110	35.50	19.50	9.75	9.75	17.75	CENTER
AESL-68E	94	30.00	16.00	9.875	8.00	4.50	END
AESL-74C	140	30.00	17.875	17.75	8.937	15.00	CENTER
AESL-ADA-5	18	14.00	10.00	4.875	4.00	4.00	CORNER
AESL-ADA-7	20	14.00	14.00	4.875	4.00	4.00	CORNER
AESL-ADA-15	24	16.00	12.00	4.875	4.00	4.00	CORNER
AESL-ADA-25	35	18.00	15.00	4.875	4.00	4.00	CORNER
AESL-ADA-55	48	25.00	15.00	4.875	4.00	4.00	CORNER
AESL-ADA-26	43	18.00	15.00	11/4.875	4.50	4.50	CORNER
AESL-ADA-56	53	25.00	15.00	11/4.875	4.50	4.50	CORNER
*	DOUBLE BOWL SINK, 2-18" X 15" I.D. BOWLS						
	<b>SPECIALITY SINKS</b>						
AESL-R12	27	12.00 DIA.		7.875			CENTER
AESL-T48	42	45.562	6.375	3.500	3.188	27.125	OFF-CENTER

*Sink list as of 2/5/2020 – additional sinks being added weekly*



## Epoxy Resin Sinks:

**Material:** All resin sinks are manufactured from the same resin material formulation as the rest of American Epoxy Scientific's resin products and meet the same chemical and physical guidelines.

**Thickness:** Sinks are molded in one piece with a 1/2" minimum thickness.



**Lip Thickness:** Drop-in style sinks will have a lip thickness of 1/4" to allow for a stronger bond to the sink cutout

**Outlet Size:** All sinks shall have a standard outlet size of 3 1/2"

**Detail:** Sinks shall have a smooth bottom that slopes 1 degree to the outlet with a 1 1/2" covered radius in the bottom corners. Sinks can be modified to accept garbage disposals, fixture holes, or troughs when required.

**Supports:** All under mount sinks should be supported as per the furniture manufacturer's guidelines.

### Sink Outlet:

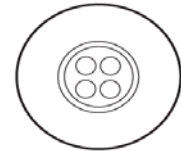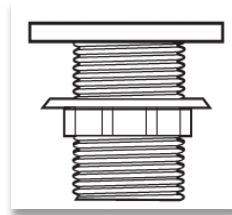

<b>Material</b>	Polypropylene
<b>Colors</b>	Black; White; Rock White; Charcoal Grey; Grey; Sand; Dark Khaki
<b>Lock Nut</b>	Black
<b>Diameter</b>	3-3/8

## Cup Sink Specification

**Material:** Polypropylene

**Color:** Black; Charcoal Grey; Rock White; Sand; Dark Khaki; Brite White

**Size:** 3"x9" and 3"x6"

Low profile cup sinks are also available in black.

**Accessories:** Standard cup sink covers available in polypropylene in black

Fume hood cup sink covers available in epoxy resin in matching colors to the fume hood top.

**Cup sinks are also available in the same epoxy resin formulation that our work-surfaces are made from.**

