

**American Epoxy
Scientific LLC.**

new england lab
laboratory furniture systems

888.635.2080 | newenglandlab.com

Under-mount sinks:

SINK	Wt. lbs.	"A" DIM	"B" DIM	"C" DIM	"D" DIM	"E" DIM	DRAIN
AESU-01C	11	9.00	6.00	6.00	3.00	4.50	CENTER
AESU-5	22	14.00	10.00	6.00	5.00	4.00	END
AESU-5C	22	14.00	10.00	6.00	5.00	7.00	CENTER
AESU-6C	32	12.00	12.00	12.00	6.00	6.00	CENTER
AESU-10C	21	16.00	8.00	7.00	4.00	8.00	CENTER
AESU-10E	21	16.00	8.00	7.00	4.00	4.00	END
AESU-13C	15	12.00	8.00	6.00	4.00	6.00	CENTER
AESU-15	30	16.00	12.00	8.00	6.00	4.00	END
AESU-15C	30	16.00	12.00	8.00	6.00	8.00	CENTER
AESU-15RC	30	16.00	12.00	8.00	4.00	4.00	CORNER
AESU-20	38	16.00	16.00	7.50	4.00	4.00	CORNER
AESU-22C	18	18.00	6.50	6.00	3.25	9.00	CENTER
AESU-23C	58	16.00	16.00	11.625	8.00	8.00	CENTER
AESU-24C	35	18.00	14.00	10.50	7.00	9.00	CENTER
AESU-25	39	18.00	15.00	8.00	4.00	4.00	CORNER
AESU-30	53	18.00	15.00	11.00	4.00	4.00	CORNER
AESU-30C	53	18.00	15.00	11.00	7.50	9.00	CENTER
AESU-33E	64	21.00	17.00	10.00	8.50	4.00	END
AESU-35E	70	21.00	18.00	10.00	10.50	4.00	END
AESU-40	50	21.375	15.50	7.00	4.00	4.00	CORNER
AESU-45	70	21.375	15.375	10.875	4.00	4.00	CORNER
AESU-49	61	22.000	14.500	10.625	4.00	4.00	CORNER
AESU-50C	53	24.00	16.00	8.000	8.00	12.00	CENTER
AESU-51C	75	24.00	16.00	11.75	8.00	12.00	CENTER
AESU-52	70	24.00	18.00	11.00	4.00	4.00	CORNER
AESU-54	51	25.00	15.00	8.00	4.00	4.00	CORNER
AESU-55	61	25.00	15.00	10.00	4.00	4.00	CORNER
AESU-57	80	25.00	15.00	13.625	4.00	4.00	CORNER
AESU-58	82	25.00	15.00	17.625	4.00	4.00	CORNER
AESU-59	65	28.00	15.00	12.00	4.00	4.00	CORNER
AESU-59E	65	28.00	15.00	12.00	7.50	11.50	CENTER
AESU-61	105	30.00	16.00	17.75	4.50	4.50	CORNER
AESU-65C	105	35.50	19.50	9.75	9.75	17.75	CENTER
AESU-67C		36.00	19.75	16.00	9.875	18.00	CENTER
AESU-68E	94	30.00	16.00	9.875	8.00	4.50	END
AESU-71C	130	43.00	15.00	11.00	7.50	21.50	CENTER
AESU-73	170	38.00	20.75	19.00	5.25	5.25	CORNER
AESU-74C	140	30.00	17.88	17.75	8.937	15.00	CENTER
AESU-75C	140	42.00	20.00	10.00	10.00	21.00	CENTER
AESE-47		22.00	13.625	12.50	9.50	11.00	FRONT-CENTER
END SINK W/SPLASH							

Sink list as of 2/5/2020- additional sinks being added weekly

Quality & Reliability



Epoxy Resin Sinks:

Material: All resin sinks are manufactured from the same resin material formulation as the rest of American Epoxy Scientific's resin products and meet the same chemical and physical guidelines.

Thickness: Sinks are molded in one piece with a ½" minimum thickness.



Lip Thickness: Drop-in style sinks will have a lip thickness of ¼" to allow for a stronger bond to the sink cutout

Outlet Size: All sinks shall have a standard outlet size of 3 ½"

Detail: Sinks shall have a smooth bottom that slopes 1 degree to the outlet with a 1 ½" covered radius in the bottom corners. Sinks can be modified to accept garbage disposals, fixture holes, or troughs when required.

Supports: All under mount sinks should be supported as per the furniture manufacturer's guidelines.

Sink Outlet:

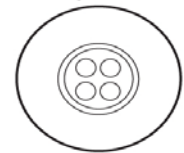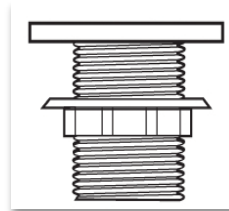

Material	Polypropylene
Colors	Black; White; Rock White; Charcoal Grey; Grey; Sand; Dark Khaki
Lock Nut	Black
Diameter	3-3/8

Cup Sink Specification

Material: Polypropylene

Color: Black; Charcoal Grey; Rock White; Sand; Dark Khaki; Brite White

Size: 3"x9" and 3"x6"

Low profile cup sinks are also available in black.

Accessories: Standard cup sink covers available in polypropylene in black

Fume hood cup sink covers available in epoxy resin in matching colors to the fume hood top.

Cup sinks are also available in the same epoxy resin formulation that our work-surfaces are made from.

